
MAX Knee

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MAX-F03YT623

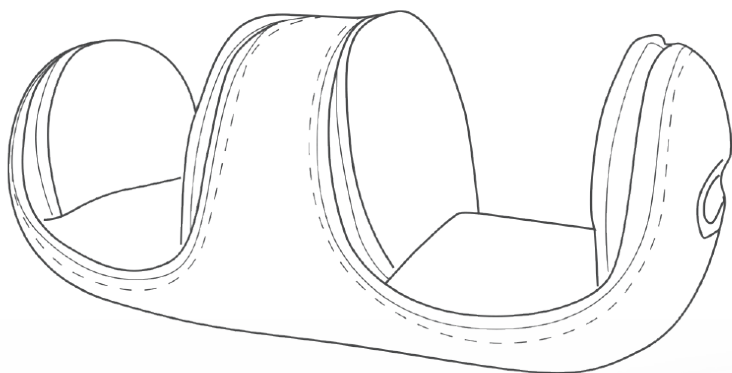maxcare®

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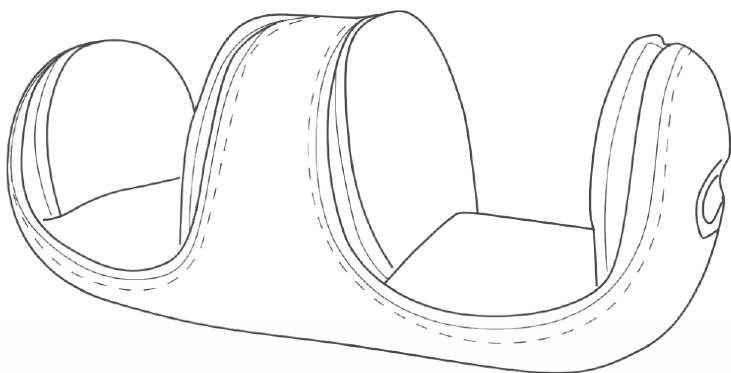
MAX-F03YT623

COMPANY PROFILE

After two decades of experience in the health industry, **maxcare** founder utilized his knowledge and expertise in massage therapy from an Asian-Pacific perspective to establish our brand - **maxcare**. Our massage therapy business was expanded to Hong Kong in 2012 and is conferred 2017 Hong Kong Top Brand Award. Based in Hong Kong, we are specializing in design & innovation, manufacturing, and retailing of massage and healthcare products.

Believing in the value of 'Technology improves our quality of life' and 'Caring for your family', **maxcare** has been innovating highly qualified products with international patents. Besides, **maxcare** products fulfill various international safety standards, namely EU CE Mark, US UL/ETL Mark, Germany GS Mark, TUV certification, and China CCC Mark.

Motivated by our value, **maxcare** founder has been leading us to be one of the most competitive brands in Hong Kong health industry. **maxcare** applies the contemporary fashion design – neat and modern – on the mechanical engineering of massage chairs. To deliver our value and offer incredible massage experience and healthcare solution to the public, **maxcare** devotes to the innovation, design, manufacturing the top-notch family massage products. By integrating the company philosophy – Spirit Of Healthy Life - into our products, **maxcare** brings health to everyone.

MAX Knee

Thank you for purchasing the Maxcare products. Before using the product, please read this manual thoroughly, particularly the Safety Precautions section, for the correct use of the product.

This manual should be well kept for easy reference. It enables you to perform a safe and efficient operation.

CONTENT

Safety Precautions	P.01 - P.03
Name and Function of Components	P.04
How to Use	P.04 - P.05
Cleaning and Maintenance	P.06
FAQs and Solutions	P.07
Technical Specifications	P.08

Safety Precautions

1. Important safety warning

- Please use well-grounding power supply suited to this product.
- Please pull out the plug after use or before cleaning to avoid injury or damage to the product.
- Do not use unauthorized accessories.
- Do not use this product outdoors.
- Do not use this product when the body is wet.
- Please read this manual instruction carefully before operation.
- Any other usage that is not listed in this manual instruction is forbidden.
- Recommended usage time of this product is 15 minutes/ session.
- Please check if any part of this products is damaged before use. If so, please do not use it and contact us for repair.
- Please keep this product away from aerosol products or where oxygen is administered.
- Please do not remove the cover.
- Please do not use it if the cover is broken or damaged.
- Please do not use it if the shell or surface is damaged.
- Please do not allow foreign objects into this product.
- Please do not allow water or liquid into this product.
- Please do not use this product under a blanket or pillow. Overheating may occur and cause fire, electric shock or injury.
- Please do not fall asleep while using this product.
- Please do not use it if drunk or feeling unwell.
- Please do not use it after taking pain killers or sedatives.
- Please do not adjust to the maximum intensity at the beginning of the massage to avoid the risk of being injured.
- Keep the plug clean from needles, rubbish or water.
- Never cover the air opening when in use. Keep it clean from dirt, lint, hair and alike.
- If the power cord is damaged, do not attempt to change on your own. Please contact us for repair for safety reasons.

2. Environment for usage

- Please do not use it in a moist environment such as bathrooms.
- Please do not use it when there is a drastic change in temperature.
- Please do not use it under heavy-dusty or caustic environment.
- Please do not use it where there is not enough space or not good ventilation.
- Please switch off the power and unplug the machine after use.
- Please do not use it when surrounding temperature is up to 50°C.
- Please do not use it when surrounding temperature is below -20°C. The lubricant can thicken quickly and increase the friction between parts.
- Turn off and pull out the plug during thunderstorms.
- Do not use outdoors.
- Do not place it near fire or directly under the sun.
- When moved from a low-temperature storage to a warm environment, the moisture inside the machine may affect, or damage, its functionality. Please rest the product within room temperature for 1 hour before use.

Safety Precautions

3. People who have the below symptoms, please consult a doctor before use

- People who have rashes or infectious skin diseases.
- People who have injuries or open wounds.
- People who have swollen or inflamed feet.
- People who are suffering from phlebitis, thrombosis, joint dysfunction, inflammation, weak bones or osteoporosis, etc.
- People who are suffering from heart disease or carrying electronic medical equipment such as heart pacemaker.
- People who have pins or screws inside the body, or artificial joints.
- People who have diabetes, sensory impairment, acute illnesses, heart diseases, or severe high blood pressure, etc.
- People who have a fever.
- People who have an abnormal or curved spine.
- Pregnant women or who has recently given birth.
- Children under 14-year old or people mentally unbalanced without being supervised are not allowed to use this machine.
- People who are told by doctor to have rest or people who feel unwell.
- The appliance has heating function. People who are insensitive to heat must be careful when using the appliance.
- People who are receiving medical treatment, especially those who are physically unfit.
- People who are at an increased risk of developing blood clots.
- People who aim for recovery.

4. Safety

- Please do not place heavy objects on this product.
- Please do not stand while using this product. Place and remove feet in and from this product while seated.
- Do not use other therapy devices while using this product.
- This product is not to be used for medical treatment purpose.
- Only one person at a time.
- Recommended usage time of this product is 15 minutes/ session.
- Check if the voltage is suitable to the specification of this product.
- Do not use wet hands to pull the plug.
- Do not allow liquid into this product, for the avoidance of electric shock or damage to this product.
- Do not pull or drag the power cord when unplugging.
- Do not damage the wires or change the circuit of this product.
- Do not alter the structure or parts of this product.
- Do not use wet cloths to clean the electric parts such as the switch and plug.
- Do not use any chemical solvent to clean this product. Use only a damp cloth to clean the surface.
- Stay away from this product during power cut to prevent injury in case power resumes suddenly.
- Stop using this product if there appears any irregularity, and please contact us immediately.

Safety Precautions

- Stop using this product if it causes discomfort, and please consult a doctor.
- This appliance is not intended for people (including children) with reduced physical, sensory or mental capabilities, or a lack of experience and knowledge, unless they are under supervision or instruction concerning the use of the appliance by a person responsible for their safety.
- Children should be supervised to ensure that they do not play with the appliance.
- If the power cord is damaged, it must be replaced by us in order to avoid a hazard.
- Do not carry this product by its power cord or use it as a handle.
- Keep the power cord away from hot or moving objects.
- Please cut off power and pull out the plug before moving.
- If the product is found damaged, broken, suffering from power leakage, or with exposed parts, please stop using it.

5. Maintenance


- This product should be maintained by us, users are prohibited to disassemble or maintain by themselves.
- Please always turn off the main power after use.
- Do not use this product if the socket is loose.
- If this product is to be left idle for long, please tie up the power cords and keep it in dry and dustless environment.
- Do not keep this product at a high temperature or near fire. Do not keep under direct sunlight.
- Please clean this product with dry cloths. Do not use thinner, benzene or alcohol.
- The mechanical components of this product are specially designed and made, no special maintenance required.
- Do not point sharp objects at this product.
- Please use it intermittently; do not keep it continuously running for long.

6. Solution of ordinary malfunctions

- It is normal if the motor makes noises during operation.
- If the rated working time is over, the power will automatically go off. If this product works continuously for too long, the temperature control will automatically turn off the machine, if so, please allow a 30-minute break before resuming.


 Product Parts and Function

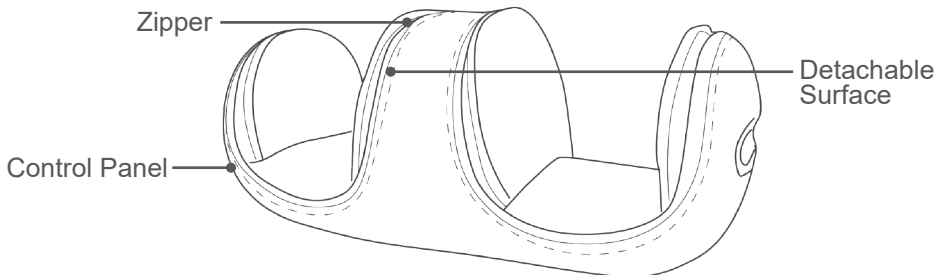 HOW TO USE

Connect to Power

Connect the plug to a power source, and the green light will be lit, which means it is ready for charging. When the charger plug and the massager charging port are connected, the signal light will turn red, implying it is charging. After 2-3 hours (actual charging time depends on the remaining battery level), the signal light will turn green, implying battery is full. Unplug the charger and use the massager.

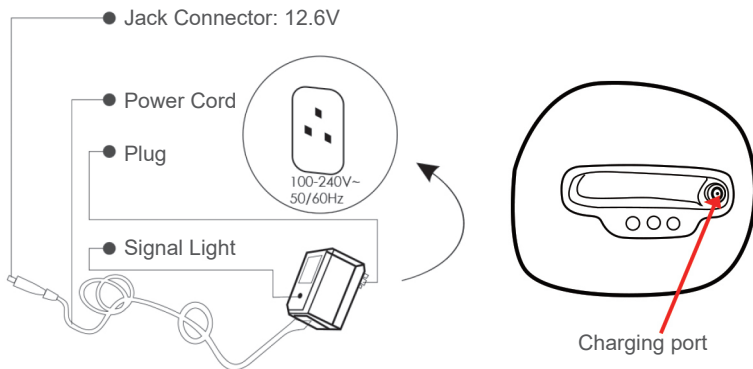 HOW TO USE**Steps**

- 1: When the charger is connected to the power source, check if the massager is charged. If the signal light is green, then it can be used; if the signal light is red, then battery level is low and it is charging.
- 2: When the battery level is low, the three buttons on the panel will start flashing, please charge.
- 3: Press the power button. If the massager is idle for more than 15 minutes, it will be automatically turned off.
- 4: Choose different programs to start the massage.

Remarks: It is not recommended to use this product while it is on charge as it may affect battery life.

Steps

-  **Power Button:** Press the power button once and enter Auto Mode 1. The signal light will be blue. Heating, vibration and air pressure will be on, with vibration being the lead. Press the power button again to enter Auto Mode 2. Heating, vibration and air pressure will be on, with air pressure being the lead. Press the power button again to turn off this product.
-  **Vibration Button:** Press the vibration button once to turn on vibration function. The signal light will be blue. Press the button again to turn off the function.
-  **Heating Button:** Press the heating button to turn on heating function. The signal light will be blue, and the temperature will be set to 40-50°C. Press the button again. The signal light will be red, and the temperature will be set to 45-55°C.

Remarks: When battery level is low, all three buttons will flash and please put the massager on charger.

 **Cleaning and Maintenance****Cleaning**

- Before cleaning, please make sure the switch has been turned off, and the power plug has been pulled out from socket.
- Please use a piece of clean and soft cotton cloth to clean this product. The cloth can be soaked in soap water first, but please dry it to damp before cleaning.

**CAUTION**

- Do not pull the plug from socket with wet hands.
- Rinsing this product with water is not allowed.
- Do not use corrosive cleaning agents like benzene or thinner to clean the product, otherwise, the color of product will fade, and the internal parts may be corrupted or broken.
- Do not disassemble the product when it is wet.

Maintenance

- When not in use, turn off the main power, and remove the plug from socket.
- Keep the product clean and away from high-temperature/ humid areas, or any environment with volatile gases.
- Keep the product in a carton, and put it in a safe, dry, well-ventilated and shaded place.
- If this product is to be left idle for long, properly cover it with a dust cover (self-purchase).
- Please keep the product away from direct sunlight, or high temperature, to avoid fading.

**CAUTION**

- Please keep the product away from direct sunlight, or high temperature, to avoid fading or product damage.

Kindly reminding

- If the product is used through the year, please check conditions as below:
 - * Burnt smell.
 - * Unstable power supply after contact with the power cord.
 - * Power cord heating up.
 - * Power cords damage.
 - * Other abnormalities.

**CAUTION**

- If above mentioned happens, please stop using the product.
- Switch off the power supply, pull out the plug, and contact us.
- Do not disassembly the product by yourself.

? FAQs and Solutions

Issues	Solutions
Product not functioning	<ol style="list-style-type: none">1. Check if the power cord is well2. Check if the power plug is connected to power source3. The product has automatically stopped working after long time of operation to protect the engine from overheating. Please stop using for 30 to 50 minutes
Control panel is not reacting and the product stops working	<ol style="list-style-type: none">1. Auto timer is up. Press the power button again to start massage.2. Turn off power and restart again.

If the abnormal situations persist after taking the above solutions, please contact us.

Technical Specifications

Product Name	: MAX Knee
Product Model	: MAX-F03YT623
Product Net Weight	: 3.5 kg
Input Voltage	: 100-240V
Rated Frequency	: 50/60Hz
Rated Power	: 12.6V 1.5A, 20W
Rated Time	: 15 min
Product Size	: 515x235x240mm

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Trademark holder: Maxcare Industrial Limited