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
Massage Mattress
with heat

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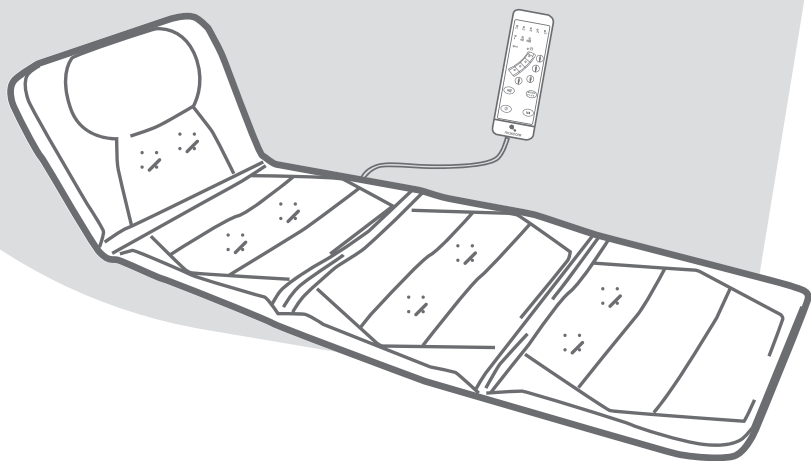maxcare®

SAFETY PRECAUTIONS 安全使用事項

On circumstances stated below, turn off the switch immediately and unplug the cord from the socket. 如有下述情況,請立即關閉開關,拔掉電源線插頭

Do not spill water and/or other liquids on the machine to prevent electric shock.
避免弄濕產品,以免觸電。

If you get unusual physical discomfort while using the machine, immediately discontinue the use of the machine and seek medical advice.
使用過程中,身體出現或感到不適時,請立即停止使用,並諮詢醫生。

Stop using and unplug the power cord if the machine produces a burnt smell and abnormal sound.
使用時發現任何異常的狀況(如冒煙、異味等),請立即停止使用,並切斷電源。

Unplug the power cord during power outage to prevent failure or an accident.
如果停電請關閉開關,切斷電源,以避免來電時產生危險。

If there is a threat of lightning, stop operation immediately and unplug the power cord.
請勿在行雷閃電時使用,並將電源插頭拔離電源插座。

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Do not wipe the product with solvents such as benzene, thinner, and diluents.
切勿使用苯或稀釋劑等具有腐蝕性的清潔劑來清潔產品。

Clean the product with a dry cloth.
請用柔軟的乾布清潔產品。

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Cases of which serious damage or death can occur when instructions are ignored.
錯誤使用,使用者可能會蒙受傷害

Do not use the product directly on exposed skin to prevent burn injuries, abrasion, etc.
請勿在裸露的肌膚上使用本產品,以免傷害皮膚。

If you get acute pain or other unusual physical discomfort while using the machine, immediately discontinue the use of the Machine.
使用過程中,身體出現或不適或痛楚時,請立即停止使用。

Do not use the product on wounded area, and for those with abnormal conditions such as acute diseases, malignant tumors, high fever, and skin diseases, consult with doctor before using the product.
皮膚有創傷的人、需安靜休養的人、患有惡性腫瘤、急性疾病、發燒、皮膚病患者等患者,使用前請先諮詢醫生,在醫生指導下方可使用。

This appliance is not intended for use by children unless they have been given supervision or instruction concerning use of the appliance by a person responsible for their safety.
未成年人士必須在正常成人監護下才使用。

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Read carefully before using the product 使用該產品前請仔細閱讀以下內容:

Strapping System 固定帶設計:

This product incorporates a unique strapping system which allows you to fasten it to the bed or other places. This stabilizes the product on the same position.

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The plug is located in the controller on the side of the product. Fasten the strap on an appropriately position and enjoy the massage after plugging the power cord and switch on the product.

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INTENSITY

Intensity button: adjusts the strength of your massage.
Choose from 3 levels of intensity.
強度按鈕: 調整振動強度,有3檔振動強度供您選擇。

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