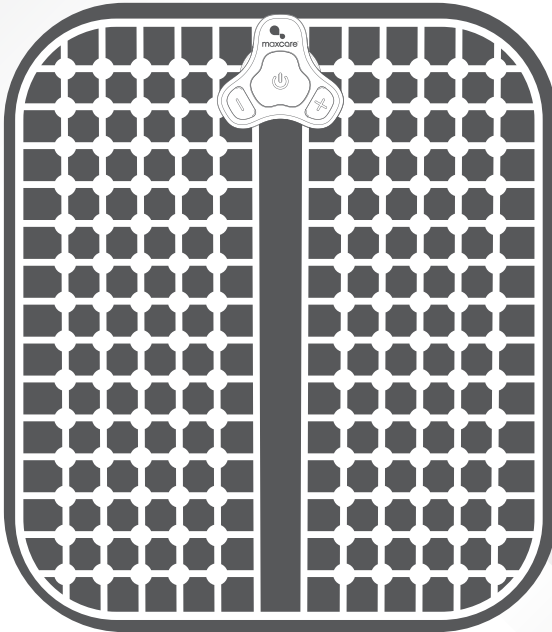


Mini Foot Pad

User Manual

MAX-M01MC0512





COMPANY PROFILE

After two decades of experience in the health industry, **maxcare** founder utilized his knowledge and expertise in massage therapy from an Asian-Pacific perspective to establish our brand - **maxcare**. Our massage therapy business was expanded to Hong Kong in 2012 and is conferred 2017 Hong Kong Top Brand Award. Based in Hong Kong, we are specializing in design & innovation, manufacturing, and retailing of massage and healthcare products.

Believing in the value of 'Technology improves our quality of life' and 'Caring for your family', **maxcare** has been innovating highly qualified products with international patents. Besides, **maxcare** products fulfill various international safety standards, namely EU CE Mark, US UL/ETL Mark, Germany GS Mark, TUV certification, and China CCC Mark.

Motivated by our value, **maxcare** founder has been leading us to be one of the most competitive brands in Hong Kong health industry. **maxcare** applies the contemporary fashion design – neat and modern – on the mechanical engineering of massage chairs. To deliver our value and offer incredible massage experience and healthcare solution to the public, **maxcare** devotes to the innovation, design, manufacturing the top-notch family massage products. By integrating the company philosophy – Spirit Of Healthy Life - into our products, **maxcare** brings health to everyone.

Mini Foot Pad

Thank you for purchasing the Maxcare products. Before using the product, please read this manual thoroughly, particularly the Safety Precautions section, for the correct use of the product.

This manual should be well kept for easy reference. It enables you to perform a safe and efficient operation.

CONTENT

Safety Precautions	P.01 - P.03
Name and Function of Components	P.04
How to Use	P.05
Cleaning and Maintenance	P.06
FAQs and Solutions	P.06
Technical Specifications	P.07



Safety Precautions



1. Important safety warning

- Keep children and pets away from the movable parts of this product.
- When in use, sliding or transporting, please ensure that there are no people or pets nearby to avoid accidents.
- Please use well-grounding power supply suited to this product.
- Please pull out the plug after use or before cleaning to avoid injury or damage to the product.
- Do not use unauthorized accessories.
- Do not use this product outdoors.
- Please read this manual instruction carefully before operation.
- Any other usage that is not listed in this manual instruction is forbidden.
- Recommended usage time of this product is 20 minutes/ session.
- Please check if any part of this products is damaged before use. If so, please do not use it and contact us for repair.
- Please do not use it if the cover is broken or damaged.
- Please do not use it if the shell or surface is damaged.
- Please do not drop anything into this product.
- Please do not fall asleep while using this product.
- Please do not use it while driving.
- Please do not use it if drunk or feeling unwell.
- Please do not use it within one hour after having meal.
- Please do not adjust to the maximum intensity at the beginning of the massage to avoid possible injury.
- Keep the plug clean from needles, rubbish or water.
- If the power cord is damaged, do not attempt to change on your own. Please contact us for repair for safety reasons.



2. Environment for usage

- Please don't use any power source with a "fast-charge" function, to avoid circuit protection.
- Please do not use it in a moist environment such as bathrooms.
- Please do not use it when there is a drastic change in temperature.
- Please do not use it under heavy-dusty or caustic environment.
- Please do not use it where there is not enough space or not good ventilation.
- Please switch off the power and unplug the machine after use.
- Please do not use it when surrounding temperature is up to 40°C.
- Please do not use it when surrounding temperature is below 5°C. The lubricant can thicken quickly and increase the friction between parts.
- Turn off and pull out the plug during thunderstorms.
- Do not use outdoors.
- Do not place it near fire or directly under the sun.
- When moved from a low-temperature storage to a warm environment, the moisture inside the machine may affect, or damage, its functionality. Please rest the product within room temperature for 1 hour before use.



Safety Precautions



3. People who have the below symptoms, please consult a doctor before use

- People who are suffering from skin conditions such as swollen skins, sunburnt or other dermatology diseases.
- People who are suffering from malignant tumors, acute diseases, heart diseases, stroke or severe hypertension.
- People who are suffering from any difficulties related to respiratory organs.
- People who are suffering from any blood circulation disorders such as thrombosis, severe aneurysms and venous tumors; or from any skin infections.
- People who are suffering from perception disorders due to diabetes.
- People who have a fever.
- People who are suffering from injuries or skin diseases.
- Children under 14-year old or people mentally unbalanced without being supervised are not allowed to use this machine.
- People who are told by doctor to have rest or people who feel unwell.
- People whose bodies are wet are forbidden to use this product.
- People who are receiving medical treatment, especially those who are physically unfit.
- People who have open wounds.
- People who aim for recovery.



4. Safety

- Please do not place heavy objects on this product.
- Do not use other therapy devices while using this product.
- This product is not to be used for medical treatment purposes.
- Only one person at a time.
- Recommended usage time of this product is 20 minutes/ session.
- Check if the voltage is suitable to the specification of this product.
- Do not use wet hands to pull the plug.
- Do not allow water to get into the product, to avoid electric shock or damage to this product.
- Do not pull or drag the power cord when unplugging. Please be gentle.
- Do not damage the wires or change the circuit of this product.
- Do not use wet cloths to clean the electric parts such as the switch and plug.
- Be away from this product during power cut to prevent injury in case power resumes suddenly.
- Stop using this product if there appears any irregularity, and please contact us immediately.
- Stop using this product if it causes discomfort, and please consult a doctor.
- This appliance is not intended for people (including children) with reduced physical, sensory or mental capabilities, or a lack of experience and knowledge, unless they are under supervision or instruction concerning the use of the appliance by a person responsible for their safety.
- Children should be supervised to ensure that they do not play with the appliance.
- If the power cord is damaged, it must be replaced by us in order to avoid a hazard.
- Keep the power cord away from hot or moving objects.
- Please cut off power and pull out plug before moving.
- If the product is found damaged, broken, suffering from power leakage, or with exposed parts, please stop using it.
- Do not leave the product unattended when it is operating.



Safety Precautions



5. Maintenance

- This product should be maintained by us, users are prohibited to disassemble or maintain by themselves.
- Please always turn off the main power after use.
- Do not use this product if the socket is loose.
- If this product is to be left idle for long, please tie up the power cords and keep it in dry and dustless environment.
- Do not keep this product at a high temperature or near fire. Do not keep under direct sunlight.
- Please clean this product with dry cloths. Do not use thinner, benzene or alcohol.
- The mechanical components of this product are specially designed and made, no special maintenance required.
- Do not point sharp objects at this product.
- Please use it intermittently; do not keep it continuously running for long.

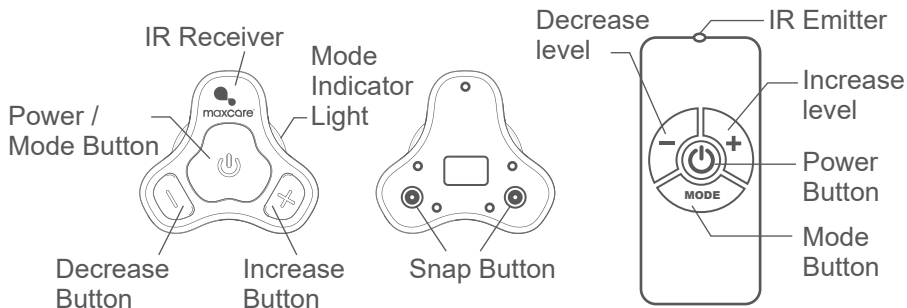
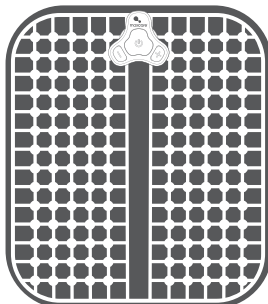





6. Solution of ordinary malfunctions

- It is normal if the motor makes noises during operation.
- If the rated working time is over, the power will automatically go off. If this product works continuously for too long, the temperature control will automatically turn off the machine, if so, please allow a 30-minute break before resuming.

 Name and Function of Components

- 1) Tri-Gear Controller X1
- 2) Foot Pad X1
- 3) USB Cable X1
- 4) Remote Control X1
- 5) CR2025 Button Cell X1
- 6) Drawstring Pocket X1



	<p>Power Button:</p> <ul style="list-style-type: none"> > Press and hold for 3 seconds to turn on the device. > Press and hold for 3 seconds to turn off the device. > Press once to select modes: 5 programs for selection (main unit only). Press "MODE" on remote control to select different programs.
	<p>Intensity Increase Button "+": To increase intensity level. 10 Intensity levels.</p>
	<p>Intensity Decrease Button "-": To decrease intensity level. 10 Intensity levels.</p>

 **How to Use**

Charging

- Plug the micro USB end of the USB cable into the USB charging port at the bottom of the device.
- Plug the other end of the USB cable into a 5V USB plug. Make sure that the input current does not exceed 2A. Too high voltage will damage the micro-processor & circuit of the device.
- The charging indicator light turns red when it is on charge and blue when it is fully charged.
- After being fully charged for 1 hour, the device can work continuously for 1 hour. If the device becomes unresponsive to the power button, the battery may be too low to operate. Please put the device on charge.

Operation

1. Wipe the PU surface of the foot pad with a piece of dry cloth before you start the massage.
2. Clean your feet before starting the massage. Make sure that both feet are completely in touch with the foot pad.

Tips: To increase the conductivity of your sole, you can apply water or water base gel on it.

3. Connect the fully-charged Tri-Gear main device to the foot pad with the metal snap buttons.
4. Press and hold the power button for 2 seconds to turn on the device. The indicator light will green.
5. To use the remote control, please point the IR emitter to the IR receiver on the main unit (within 1m). Adjust your massage by pressing the corresponding buttons. (Please insert the button cell battery before using the remote control)
6. There are 10 EMS intensity levels (1 to 10) in each mode. Press the increase “+” and decrease “-” button to adjust the intensity level of the massage.
7. Each program stops after 20 minutes of operation. Once the intensity is set, the program will start counting down from 20 minutes.
8. After a massage program ends, please make sure the device is off.
9. Please uninstall the main device after use.

Reminders

1. Tri-Gear has a protection function. If the device has been turned on for 2 seconds without skin contact, it will automatically turn off again.
2. When the battery voltage is lower than 2.4V, the red indicator light will start blinking for 5 seconds and then the device will be automatically turned off as a signal to charge battery.

Program	Colour of the Indicator Light (Flashing Light Cycle)	Massage Modes	
1	Green	Relaxation	
2	Purple	Soothing	
3	Yellow	Calming	
4	Blue	Kneading	
5	White	Shiatsu	

3. Begin with the lowest intensity level is recommended. Gradually increase the intensity level to find the most suitable point.
4. Keep the PVC film and put it back on top of the foot pad every time after use to maintain cleanliness.
5. Do NOT peel off the gel pads on the electrode sheets unless replacement is needed.

Cleaning and Maintenance

Main Device & Foot Pad

- Always keep the main device clean by using a piece of dry and soft cloth to wipe the surface.
- Do not spill liquid on the device or immerse it in water.
- Keep it in a cool and dark place.
- Do not clean with chemicals.
- Keep out of the reach of children and pets.

FAQs and Solutions

No.	Symptom	Reasons	Solution
1	No response	Battery low	Put on charge
2	Indicator light blinking	Battery low	Put on charge
3	No feeling after turning on the device	- Poor skin conductance - Too low intensity	- Apply some water-based gel or water on your skin before use - Gradually increase the intensity level
4	Indicator light goes out after 2 seconds	- No skin contact - Poor skin conductance	- Check if the foot pad is directly attached to the skin - Check if the electronic controller is firmly connected with the connectors on the foot pad - Apply some water-based gel or water on your skin before use
5	Normal indicator lighting but no feeling	- Too low intensity - Poor skin conductance - Non-water-based gel is used - Gel pad is worn and greasy - Battery low	- Gradually increase the intensity level - Apply some water-based gel or water on your skin before use - Use water-based gel - Clean the foot pad with a damp cloth or paper towel - Put on charge
6	No buttons work	- Microprocessor hanged - Battery low	- Restart the device - Put on charge
7	No response to the remote control	- Battery low	Change the batteries of the remote control

If the abnormal situations persist after taking the above solutions, please contact us.

 **Technical Specifications**

Product Name	: Mini Foot Pad
Product Model	: MAX-M01MC0512
Rated Voltage (Input)	: 1Hz—33.3Hz
Rated Voltage (Output)	: 86 V (1K load resistance)
Rated Power	: 105mW
Max. Current	: 35 mA (1K load resistance)
Rated time	: 20 minutes
Battery	: Main Unit: Rechargeable 3.7V Lithium Battery Remote Control: CR2025 Button Cell
Product Size	: Main Unit: 65 x 62 x 17 mm Foot Pad: 310 x 280 x 2 mm



商標持有：美斯凱實業有限公司
Trademark holder: Maxcare Industrial Limited

客戶服務熱線 Customer Service Hotline :
香港 Hong Kong - (852) 2950-9383
中國 China - (86) 400 882-1883

www.maxcare.net
www.maxcareonlineshop.com