

ION SONIC TOOTH CLEANER

User Manual

MAX-M01MC0539




maxcare®



COMPANY PROFILE

After two decades of experience in the health industry, **maxcare** founder utilized his knowledge and expertise in massage therapy from an Asian-Pacific perspective to establish our brand - **maxcare**. Our massage therapy business was expanded to Hong Kong in 2012 and is conferred 2017 Hong Kong Top Brand Award. Based in Hong Kong, we are specializing in design & innovation, manufacturing, and retailing of massage and healthcare products.

Believing in the value of 'Technology improves our quality of life' and 'Caring for your family', **maxcare** has been innovating highly qualified products with international patents. Besides, **maxcare** products fulfill various international safety standards, namely EU CE Mark, US UL/ETL Mark, Germany GS Mark, TUV certification, and China CCC Mark.

Motivated by our value, **maxcare** founder has been leading us to be one of the most competitive brands in Hong Kong health industry. **maxcare** applies the contemporary fashion design – neat and modern – on the mechanical engineering of massage chairs. To deliver our value and offer incredible massage experience and healthcare solution to the public, **maxcare** devotes to the innovation, design, manufacturing the top-notch family massage products. By integrating the company philosophy – Spirit Of Healthy Life - into our products, **maxcare** brings health to everyone.

ION SONIC TOOTH CLEANER

Thank you for purchasing the Maxcare products. Before using the product, please read this manual thoroughly, particularly the Safety Precautions section, for the correct use of the product.

This manual should be well kept for easy reference. It enables you to perform a safe and efficient operation.

CONTENT

Safety Precautions	P.01 - P.03
Name and Function of Components	P.04
How to Use	P.04 - P.06
Cleaning and Maintenance	P.07 - P.08
Technical Specifications	P.09



Safety Precautions



1. Important safety warning

- For adults only.
- Do not put this product into water or let water get into it.
- Do not drop or misuse this product.
- Please use well-grounding power supply suited to this product.
- Do not leave this product unattended when working or on charge.
- Keep children and pets away from the movable parts of this product
- When in use, sliding or transporting, please ensure that there are no people or pets nearby to avoid accidents.
- Please use well-grounding power supply suited to this product.
- Please pull out the plug after use or before cleaning to avoid injury or damage to the product.
- Do not use unauthorized accessories.
- Do not use this product outdoors.
- Please read this manual instruction carefully before operation.
- Any other usage that is not listed in this manual instruction is forbidden.
- Recommended usage time of this product is 10 minutes/ session.
- Please check if any part of this products is damaged before use. If so, please do not use it and contact us for repair.
- Please do not use it if the cover is broken or damaged.
- Please do not use it if the shell or surface is damaged.
- Please do not drop anything into this product.
- Please do not fall asleep while using this product.
- Please do not use it while driving.
- Please do not use it if drunk or feeling unwell.
- Please do not use it within one hour after having meal.
- Please do not adjust to the maximum intensity at the beginning of the massage to avoid possible injury.
- Keep the plug clean from needles, rubbish or water.
- If the power cord is damaged, do not attempt to change on your own. Please contact us for repair for safety reasons.



2. Environment for usage

- Please don't use any power source with a "fast-charge" function, to avoid circuit protection.
- Please do not use it when there is a drastic change in temperature.
- Please do not use it under heavy-dusty or caustic environment.
- Please do not use it where there is not enough space or not good ventilation.
- Please switch off the power and unplug the machine after use.
- Please do not use it when surrounding temperature is up to 40°C.
- Please do not use it when surrounding temperature is below 5°C. The lubricant can thicken quickly and increase the friction between parts.



Safety Precautions

- Turn off and pull out the plug during thunderstorms.
- Do not use outdoors.
- Do not place it near fire or directly under the sun.
- When moved from a low-temperature storage to a warm environment, the moisture inside the machine may affect, or damage, its functionality. Please rest the product within room temperature for 1 hour before use.



3. People who have the below symptoms, please consult a doctor before use

- If you have periodontal disease, please consult the dentist before using it.
- People who are suffering from malignant tumors, acute diseases, heart diseases, stroke or severe hypertension.
- People who are suffering from any difficulties related to respiratory organs.
- People who are suffering from any blood circulation disorders such as thrombosis, severe aneurysms and venous tumors; or from any skin infections.
- People who are suffering from perception disorders due to diabetes.
- People who are suffering from osteoporosis.
- People who are suffering from heart disease or carrying electronic medical equipment such as heart pacemaker.
- People who have a fever.
- Pregnant women or women that are in menses period.
- People who are suffering from injuries or skin diseases.
- Children under 14-year old or people mentally unbalanced without being supervised are not allowed to use this machine.
- People who are told by doctor to have rest or people who feel unwell.
- People whose bodies are wet are forbidden to use this product.
- People who are receiving medical treatment, especially those who are physically unfit.
- People who have open wounds.
- People who aim for recovery.



4. Safety

- Please do not place heavy objects on this product.
- Do not use other therapy devices while using this product.
- This product is not to be used for medical treatment purposes.
- Only one person at a time.
- Recommended usage time of this product is 10 minutes/ session.
- Check if the voltage is suitable to the specification of this product.
- Do not use wet hands to pull the plug.
- Do not allow water to get into the product, to avoid electric shock or damage to this product.
- Do not pull or drag the power cord when unplugging. Please be gentle.
- Do not damage the wires or change the circuit of this product.
- Do not use wet cloths to clean the electric parts such as the switch and plug.
- Be away from this product during power cut to prevent injury in case power resumes suddenly.



Safety Precautions

- Stop using this product if there appears any irregularity, and please contact us immediately.
- Stop using this product if it causes discomfort, and please consult a doctor.
- This appliance is not intended for people (including children) with reduced physical, sensory or mental capabilities, or a lack of experience and knowledge, unless they are under supervision or instruction concerning the use of the appliance by a person responsible for their safety.
- Children should be supervised to ensure that they do not play with the appliance.
- If the power cord is damaged, it must be replaced by us in order to avoid a hazard.
- Keep the power cord away from hot or moving objects.
- Please cut off power and pull out plug before moving.
- If the product is found damaged, broken, suffering from power leakage, or with exposed parts, please stop using it.
- Do not leave the product unattended when it is operating.
- Do not turn on the product when it is being charged.



5. Maintenance

- This product should be maintained by us, users are prohibited to disassemble or maintain by themselves.
- Please always turn off the main power after use.
- Do not use this product if the socket is loose.
- If this product is to be left idle for long, please tie up the power cords and keep it in dry and dustless environment.
- Do not keep this product at a high temperature or near fire. Do not keep under direct sunlight.
- Please clean this product with dry cloths. Do not use thinner, benzene or alcohol.
- The mechanical components of this product are specially designed and made, no special maintenance required.
- Do not point sharp objects at this product.
- Please use it intermittently; do not keep it continuously running for long.



6. Solution of ordinary malfunctions

- When the work tips touches the teeth, it makes a buzzing sound, which is the working state.


Product Parts and Function

Product Structure:

1. Work Tips x3
2. Connection ring
3. Power
4. USB port
5. USB charging cable
6. Brush head
7. Small spanner
8. Dental mirror

Product principle

The Dental Calculus Remover Use sound waves to smash stubborn tooth stains and dental calculus. The work tips can be effectively in the deep blind zone in your mouth. Can effectively remove long-term tooth tartar.

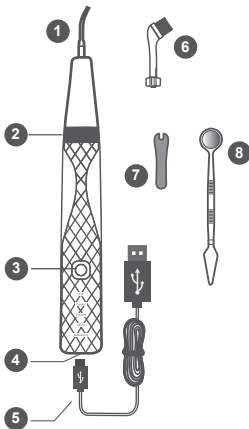

HOW TO USE

1. Make sure the tooth cleanser is fully charged before you use.
2. Press power button for turn on the device and light up to enter Soft mode, make sure you need to touch the connection ring to output the negative ion, negative ion on and ultrasonic wave is Soft mode. Press power button again for Comfort mode, negative ion is unchanged and the ultrasonic wave will change to Comfort mode. And so on, there are 5 modes. You should use Soft mode at first, and then adjust the mode that you need.



Mode:

- | | | |
|------------|------------------|--------------|
| (1)Soft | (2) Comfort | (3) Standard |
| (4) Strong | (5) Super strong | |

3. Select suitable work tips to install to the tooth cleanser with the small spanner. Please do not twist tightly the interface with the small spanner. Choose the suitable mode and keep the tip in a fixed position for 3 to 5 minutes, until the stains are clearly falling off and then move elsewhere. Repeat several times.
4. The brush head is use for cleaning the teeth when you done the treatment, it can not be use for electric toothbrush.
5. If it is low battery, the indicator light will flash. When it is charging, the indicator light will flash green light. When it is fully charged, the indicator light will on green light.
6. Press the power button for few second to turn off the cleanser. The cleanser automatic shut-down when over the 10 minutes preset time.




HOW TO USE
How to change the work tips:

- Please confirm the work tips was locked to the cleanser, and use the small spinner to adjust the work tips.
- Do not use other tools to change the work tips.

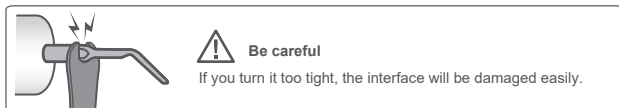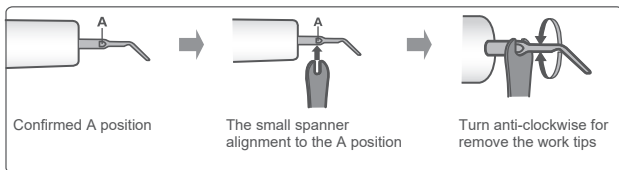
**How to use the work tips**

Slightly press to remove any dental stains. Please adjust elapsed time for different situation.


- ⚠** Do not use the work tips to clean the no stains tarter area.

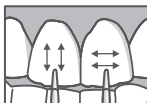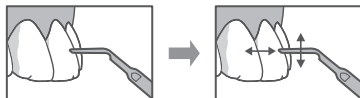


Tooth stains

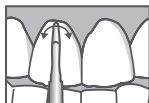


Between gums and tooth

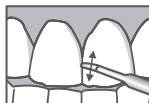


Tooth space

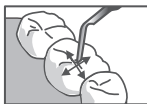

Not smooth tooth surface

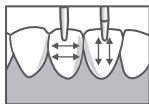


Tooth stains at the back

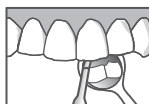


Please use the dental mirror when you cleaning the back side of tooth


 HOW TO USE

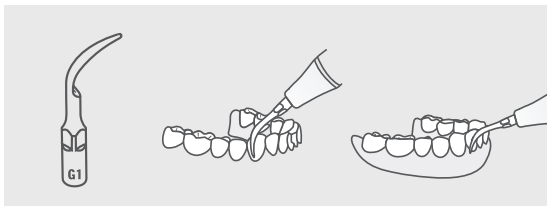

G1 is use for remove dental stains

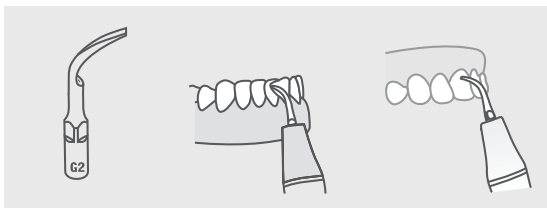

G2 is wider and use for remove the bigger dental stains

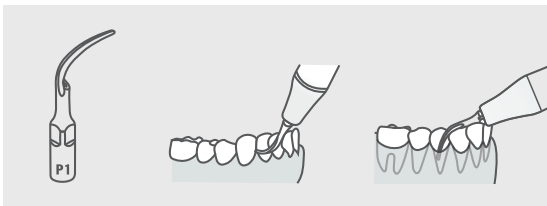

P1 is use for remove dental stains



Cleaning and Maintenance

Cleaning

- Before cleaning, please make sure the switch has been turned off, and the power plug has been pulled out from socket.
- Please use a piece of clean and soft cotton cloth to clean this product. The cloth can be soaked in soap water first, but please dry it to damp before cleaning.
- Work tips and brush head can cleanse.



CAUTION

- Do not pull the plug from socket with wet hands.
- Rinsing this product with water is not allowed.
- Do not use corrosive cleaning agents like benzene or thinner to clean the product, otherwise, the color of product will fade, and the internal parts may be corrupted or broken.
- Do not disassemble the product when it is wet.

Maintenance

- Make sure that the plug has been disconnected from the socket before maintenance.
- Keep the product clean and away from high-temperature/ humid areas, or any environment with volatile gases.
- If this product is to be left idle for long, properly cover it with a dust cover (self-purchase).



CAUTION

- Please keep the product away from direct sunlight, or high temperature, to avoid fading or product damage.


Kindly reminding

- If the product is used through the year, please check conditions as below:
 - * Burnt smell.
 - * Unstable power supply after contact with the power cord.
 - * Power cord heating up.
 - * Power cords damage.
 - * Other abnormalities.





CAUTION

- If above mentioned happens, please stop using the product.
- Switch off the power supply, pull out the plug, and contact us.
- Do not disassembly the product by yourself.

 **Cleaning and Maintenance****How to maintain**

- Do not use it for a long time, please keep it in a dry place.
- Don't wash the device and power adapter with water.
- Pay attention to the failure of the product.
- Please don't use Diluents, Benzene and Solvents ect Kerosene, Chemicals, Scouring powder.
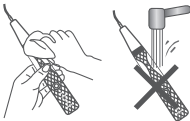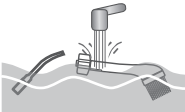
<p>The device can't be cleaned directly with water. And the dirt on the device needs to be wiped with rag.</p>	
<p>Metal heads for tooth cleaning and brush head can be washed with water.</p>	



Technical Specifications

Product Name	: ION SONIC Tooth Cleaner
Product Model	: MAX-M01MC0539
Rated Voltage (input)	: 5V1A
Rated Voltage (output)	: 3.7V/650mA
Rated Power	: 1.85W
Frequency	: 37KHZ+/-3KHZ
Rated Voltage	: 3.7V



商標持有：美斯凱實業有限公司
Trademark holder: Maxcare Industrial Limited

客戶服務熱線 Customer Service Hotline :
香港 Hong Kong - (852) 2950-9383
中國 China - (86) 400 882-1883

www.maxcare.net
www.maxcareonlineshop.com