# SMART NECK

—♥—*Massager*

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MAX-Z05C73G

maxcare®

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MAX-Z05C73G

## COMPANY PROFILE

After two decades of experience in the health industry, **maxcare** founder utilized his knowledge and expertise in massage therapy from an Asian-Pacific perspective to establish our brand - **maxcare**. Our massage therapy business was expanded to Hong Kong in 2012 and is conferred 2017 Hong Kong Top Brand Award. Based in Hong Kong, we are specializing in design & innovation, manufacturing, and retailing of massage and healthcare products.

Believing in the value of 'Technology improves our quality of life' and 'Caring for your family', **maxcare** has been innovating highly qualified products with international patents. Besides, **maxcare** products fulfill various international safety standards, namely EU CE Mark, US UL/ETL Mark, Germany GS Mark, TUV certification, and China CCC Mark.

Motivated by our value, **maxcare** founder has been leading us to be one of the most competitive brands in Hong Kong health industry. **maxcare** applies the contemporary fashion design – neat and modern – on the mechanical engineering of massage chairs. To deliver our value and offer incredible massage experience and healthcare solution to the public, **maxcare** devotes to the innovation, design, manufacturing the top-notch family massage products. By integrating the company philosophy – Spirit Of Healthy Life - into our products, **maxcare** brings health to everyone.

# SMART NECK

——*Massager*

Thank you for purchasing the Maxcare products. Before using the product, please read this manual thoroughly, particularly the Safety Precautions section, for the correct use of the product.

This manual should be well kept for easy reference. It enables you to perform a safe and efficient operation.

## CONTENT

---

Safety Precautions .....	P.01 - P.02
Name and Function of Components .....	P.03
How to Use .....	P.04
Charging Instruction .....	P.05
Cleaning and Maintenance .....	P.06
Troubleshooting .....	P.07
Technical Specifications .....	P.08

## Safety Precautions

### 1. Important safety warning

- Please read and comply with the user manual before operation.
- Do not put any objects into the gap between the massage head and the shell.
- Keep product away from sources of water to prevent short circuit caused by ingress of water into the product.
- Please use a USB cable which is in accordance with the specification of the product.
- Do not leave the product unattended while it is working or charging.
- Please turn off the power and unplug after charging. Do not keep plugging in the charging cable, so that the product is in a charging state for a long time.
- Disconnect the product after use or before cleaning to prevent injury or damage.
- Do not use unauthorized accessories.
- Do not use outdoors.
- Please check if any part of this product is damaged before use. If so, please stop using it immediately and do not dismantle or repair by yourself. You should contact us for repair.
- Please do not use it if drunk or feeling unwell.
- Please do not use it within one hour after having meal.
- Please do not adjust to the maximum intensity at the beginning of the massage to avoid possible injury.
- For those with a neck injury or recent surgery, please consult the doctor before use and follow the instructions carefully.
- Do not wash the product.

### 2. Environment for usage

- Do not use any power source with fast charging function to prevent the circuit protection.
- Do not use it in a humid environment such as a bathroom.
- Do not use it when there is a drastic change in temperature.
- Do not use it under heavy-dusty or caustic environment.
- Do not use it in an environment with poor ventilation.
- Do not use it when the surrounding temperature is up to 40°C.
- Turn off and pull out the plug during thunderstorms.
- Do not place it near fire or direct sunlight.

### 3. People who have the below symptoms, please consult a doctor before use

- People who are suffering from skin conditions such as swollen skins, sunburnt, acute purulent inflammation or other dermatological diseases.
- People who suffer from malignant tumors, acute diseases, heart diseases, stroke or severe hypertension.
- People who suffer from respiratory disorders.
- People who suffer from any bleeding tendency, thrombosis, severe aneurysm or acute venous aneurysm and other symptoms of blood circulation disorder.
- People who suffer from perception disorders caused by diabetes.
- People who suffer from osteoporosis.
- People who suffer from heart disease or carrying electronic medical device such as a heart pacemaker.
- People who have a fever.
- Pregnant or menstruating women.
- People who suffer from injuries or skin diseases.
- Children under 14-year-old or who are delirious without being supervised are not allowed to use this product.
- For those who need to rest or feel unwell after being instructed by a doctor.
- This product has a hot surface. It must be used carefully by those who are sensitive to heat.

## Safety Precautions

### 4. Safety

- Do not use other therapy devices while using this product.
- Do not use for therapeutic purposes.
- Recommended usage time of this product is 15 minutes/ session.
- Check if the voltage in accordance with the specification of this product.
- Do not use wet hands to pull the plug.
- Do not allow liquid to get into the product, to prevent electric shock or damage to this product.
- Do not pull or drag the power cord when unplugging. Please be gentle.
- Do not damage the wires or change the circuit of this product.
- Do not use wet cloths to clean the electric parts such as the switch and USB cable. Please clean this product with dry cloths. Do not use thinner, benzene or alcohol.
- Be away from this product during a power cut to prevent injury in case power resumes suddenly.
- Stop using this product if there is any irregularity, and please contact us immediately.
- Stop using this product if you feel any physical abnormality. Please consult a doctor.
- This product is not intended for people (including children) who are physically handicapped, sensory or neurological deficits, or lack of experience and knowledge, unless they are under supervision or instruction concerning the use of the product by their guardian.
- Children should be supervised to ensure they do not play it as a toy.
- Do not place heavy objects on top of the product and cannot collide with any hard objects.
- This product should be maintained by us, users are prohibited to dismantle, repair or maintain by themselves.
- Keep away from fire, high temperature or flammable, to avoid malfunctions and dangers that may cause deformation, discoloration or damage.
- Do not use the battery after it has been immersed in water.
- Do not dismantle, collide, squeeze the lithium battery or throw it into the fire. If there is serious swelling, please stop using it.
- This product is a personal care massager and cannot replace medicines and medical devices.
- Please be reminded to turn off the power before moving or changing the massage position to avoid strong irritation.
- Please wear the product properly, otherwise an accident or physical discomfort may occur.


### 5. Maintenance

- This product should be maintained by us, users are prohibited to dismantle or repair by themselves.
- If the product will not be used for a long time, please use the attached dust cover to avoid dust.
- Please clean this product with dry cloths. Do not use thinner, benzene or alcohol.
- Keep the product clean and away from the places with high humidity or high temperature and environment with volatile gases.
- Do not place it near fire or direct sunlight, to prevent product color fading.


 Name and Function of Components

## Product Structure

1. Body
2. Speaker
3. Electrode pads
4. Increase + / Decrease – Intensity button
5. Modes selection button
6. USB port
7. Remote control
8. Power button
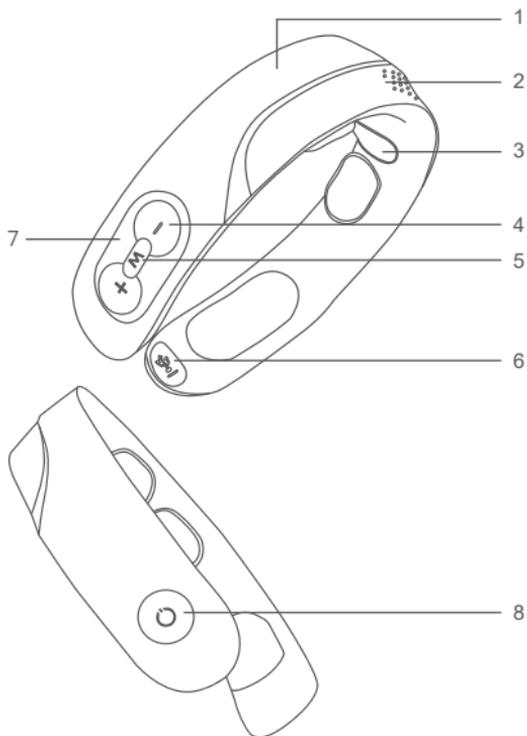
## How to Wear

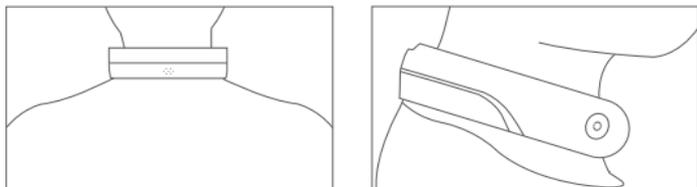Please ensure to wet your neck with a wet towel before use. Then put the product closer to the neck and adjust the position, make sure the electrode pads adhere firmly to the skin, then turn it on and select the suitable mode and intensity.

## HOW TO USE

### 1. Pairing


Press and hold the "Power" Button for 2-3 seconds, until the red light blinks, then press any button on the remote control. The red light will turn off after pairing is successful.

### 2. Operating instruction

1. Press the "⏻" button to turn ON, it will prompt "Power ON" and start with the default mode with intensity 1.
2. When it is operating, press the "⏻" button to turn OFF, it will prompt "Power OFF".
3. Remote control instruction
  - Press "+" button to increase the intensity
  - Press "-" button to decrease the intensity
  - Press "M" button to switch the mode, 6 modes switch cyclically, and the intensity will default as level 1.
    - Mode 1: Gua Sha Scraping mode
    - Mode 2: Tui na mode
    - Mode 3: Acupuncture hammering mode
    - Mode 4: Tapping mode
    - Mode 5: Relax mode
    - Mode 6: Knocking mode
  - Press and hold "M" button to turn ON/OFF the heat function
4. The preset massage time is 10 minutes, it will automatically turn OFF when finished.

### 3. Bluetooth Pairing

1. Turn ON the mobile device's Bluetooth and the product as well.
2. Search "iRest\_C73" on your mobile device and connect. Then you can enjoy the music.

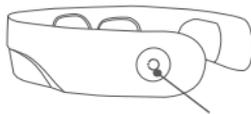1. Turn on the product

2. Search Bluetooth device "iRest\_C73" on your mobile phone.


### 4. Voice Prompt

Power On / Power Off / Pairing successful  
 Mode 1 / Mode 2 / Mode 3 / Mode 4 / Mode 5 / Mode 6 /  
 Intensity 1 / Intensity 2 / Intensity 3 / Intensity 4 / Intensity 5 / Intensity 6 /  
 Intensity 7 / Intensity 8 / Intensity 9 / Intensity 10 / Intensity 11 / Intensity 12 / Intensity  
 13 / Intensity 14 / Intensity 15 / Intensity 16 /  
 Heat on / Heat off / Time is over / Low battery, please charge.

 **Charging Instruction**


1. When the indicator light is dimmed, which means low battery, please charge the product with included USB cable or cables in accordance with the following specification.

- 1 With USB plug
- 2 Input voltage is 5V

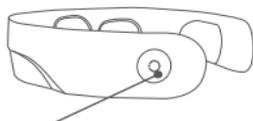Indicated light: The yellow green light will turn on when charging.  
It will turn off when finished.

Plug into  
the power

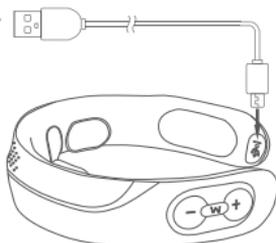2. The product cannot be turned on while charging.
3. If the product is charged during working, it will stop automatically.

**CAUTION**

Do not use the power adapter voltage above 5V,  
otherwise it could damage the product.

 **Cleaning and Maintenance****Cleaning**

- Please ensure to turn off the power before cleaning.
- Use a soft dry cloth or slightly damp cloth to wipe the product surface.

**CAUTION**

- Do not use corrosive liquid such as benzene or thinner to clean the product, to prevent color fading and damage.
- Do not wash or spill any liquid on the product.
- When the product is damp, stop using it and do not dismantle. Leave it to dry naturally.

**Maintenance**

- Keep the product clean and away from the places with high humidity or high temperature and environment with volatile gases.
- If the product will not be used for a long period, please use the included dust cover to prevent dust.

**CAUTION**

- Do not place it near fire or direct sunlight, to prevent product color fading.
- It's recommended to charge the product at least once per month, to preserve the lifespan of the lithium battery
- It's recommended to use it at least every 3 months.
- Keep out of reach of children.
- Do not place heavy objects on top of the product.

**Kindly reminding**

- The product has built-in battery. The indicator red light will blink and the product will turn off automatically when the battery is low. Please charge before use.
- The product will automatically turn off for power saving if it does not adhere to the skin within 20 seconds.

**CAUTION**

- Please ensure that the neck is wet enough to enjoy a better massage experience.
- Please sit properly and wear the product properly, do not swing your head while using.
- Please remove the protective film before use.

## Troubleshooting

The following symptoms do not indicate malfunction. Please check the below phenomenon to find out the solution before contact customer services department. If the problem still cannot be solved, please do not dismantle or repair the product. You should stop using it and contact us immediately.

Issues	Possible Causes	Solutions
Not working	Low battery	Recharge the product
	The product is charging	Use after charging
The charging indicator light does not show or cannot charge	USB charging cable or adapter does not connect properly	Connect again
Display normal, but no sensation is felt in the electrode pads	The electrode pads are not adhering firmly to the skin	The electrode pads must be adhered directly to the skin
	Hairs, clothing or accessories may prevent electrode pads contact with the skin	Remove the obstacle
	Mode error	Switch mode
	Poor conductivity	Apply water or moisturizer to the skin
	The intensity maybe too low	Please increase the intensity
	The electrode pads are dirty	Please clean the surface of electrode pads
Tingling sensation	The electrode pads are not adhered firmly to the neck	Wet the neck with wet towel before putting the electrode pads adhere to skin
	Massage time too long	It's recommended to use 10 minutes at a time, twice a day
	Massage intensity too high	Adjust lower intensity
	Skin too dry	Apply water or moisturizer to the skin
	There is a wound	Do not put the electrode pads on or near the wound

## Technical Specifications

Product Name	: Smart Neck Massager
Product Model	: MAX-Z05C73G
Rated Voltage	: 3.7V 
Rated Power	: 3W
Rated Time	: 10 minutes
Charging Time	: About 2 hours, available for about 1.5 hour
Product Size	: 145 x 152 x 36mm
Packaging Size	: 195 x 195 x 74mm
Net Weight	: 148g
Product Weight	: 480g
Components	: Main product, remote control, USB charging cable

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Trademark holder: Maxcare Industrial Limited