
# PEDICURE INSTRUMENT

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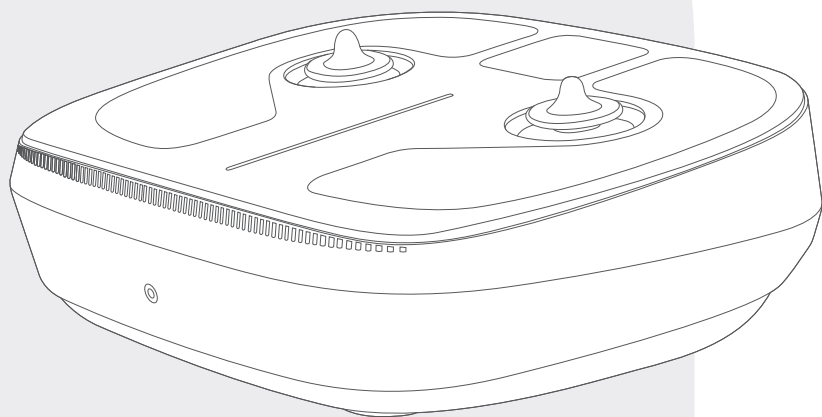MAX-F04G5HFFM

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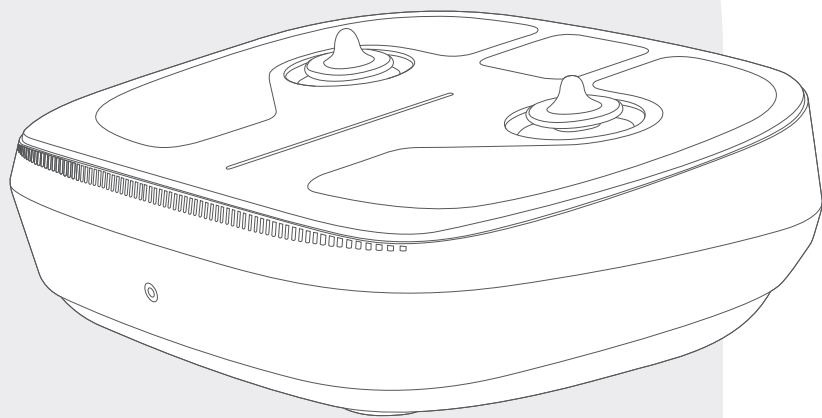
MAX-F04G5HFFM

## COMPANY PROFILE

After two decades of experience in the health industry, **maxcare** founder utilized his knowledge and expertise in massage therapy from an Asian-Pacific perspective to establish our brand - **maxcare**. Our massage therapy business was expanded to Hong Kong in 2012 and is conferred 2017 Hong Kong Top Brand Award. Based in Hong Kong, we are specializing in design & innovation, manufacturing, and retailing of massage and healthcare products.

Believing in the value of 'Technology improves our quality of life' and 'Caring for your family', **maxcare** has been innovating highly qualified products with international patents. Besides, **maxcare** products fulfill various international safety standards, namely EU CE Mark, US UL/ETL Mark, Germany GS Mark, TUV certification, and China CCC Mark.

Motivated by our value, **maxcare** founder has been leading us to be one of the most competitive brands in Hong Kong health industry. **maxcare** applies the contemporary fashion design – neat and modern – on the mechanical engineering of massage chairs. To deliver our value and offer incredible massage experience and healthcare solution to the public, **maxcare** devotes to the innovation, design, manufacturing the top-notch family massage products. By integrating the company philosophy –Spirit of Healthy Life - into our products, **maxcare** brings health to everyone.

# PEDICURE INSTRUMENT

Thank you for purchasing Maxcare's products. Before using the product, please read this manual thoroughly, particularly the Safety Precautions section, for the correct use of the product.

This manual should be well kept for easy reference. It enables you to perform a safe and efficient operation.

## CONTENT

---

Safety Precautions .....	P.14-P.16
Name and Function of Components .....	P.17
How to Use .....	P.18-P.21
Cleaning and Maintenance .....	P.22
FAQs and Solutions .....	P.23
Technical Specifications .....	P.24

## Safety Precautions

### 1. Important safety warning

- Please read and comply with the user manual before operation.
- Keep product away from sources of water to prevent short circuit caused by ingress of water into the product.
- Please use an adapter which is in accordance with the specification of the product.
- Disconnect the product after use or before cleaning to prevent injury or damage.
- Do not use outdoors.
- Be careful when using the heating function, to prevent possible burns.
- Recommended usage time of this product is 15 minutes/ session.
- Please do not keep this product running continuously for over 30 minutes.
- Please check if any part of this product is damaged before use. If so, please stop using it immediately and do not dismantle or repair by yourself. You should contact us for repair.
- Please do not use it if drunk or feeling unwell.
- Please do not use it within one hour after having meal.
- Please do not adjust to the maximum intensity at the beginning of the massage to avoid possible injury.
- Do not wash the product

### 2. Environment for usage

- Do not use it in a humid environment such as a bathroom.
- Do not use it when there is a drastic change in temperature.
- Do not use it under heavy-dusty or caustic environment.
- Do not use it in an environment with poor ventilation.
- Please turn off the power and remove the adapter from socket.
- Do not use it when the surrounding temperature is up to 40°C.
- Please do not use it when surrounding temperature is below 5°C. The lubricant will get thicker quickly and increase the friction between parts.
- Turn off and unplug during thunderstorms.
- Do not place it near fire or direct sunlight.
- When moved from a low-temperature storage to a warm environment, the moisture inside the machine may affect, or damage, its functionality. Please rest the product within room temperature for 1 hour before use.

## Safety Precautions

### 3. People who have the below symptoms, please consult a doctor before use

- People who have rashes or infectious skin diseases.
- People who have diabetes, sensory impairment, malignant tumors, acute illnesses, heart diseases, or severe high blood pressure, etc.
- People who suffer from any bleeding tendency, thrombosis, severe aneurysm or acute venous aneurysm and other symptoms of blood circulation disorder.
- People who are suffering from swollen feet or inflamed feet.
- People who are suffering from phlebitis, thrombosis, joint dysfunction, inflammation, weak bones or osteoporosis, etc.
- People who suffer from heart disease or carrying electronic medical device such as a heart pacemaker.
- People who have pins or screws inside the body, or artificial joints.
- People who have a fever.
- Pregnant or menstruating women.
- Children under 14-year-old or who are delirious without being supervised are not allowed to use this product.
- For those who need to rest or feel unwell after being instructed by a doctor.
- The product has a heating function. It must be used carefully by those who are sensitive to heat.

### 4. Safety

- Do not place heavy objects on top of the product.
- Please place the product on flat ground.
- Do not operate under a blanket or pillow. Overheating can occur causing fire, electric shock, or other injuries.
- Do not stand while using the product.
- The product can only be used by a single user at a time.
- Check if the voltage in accordance with the specification of this product.
- Do not use wet hands to pull the plug or touch the control panel.
- Do not allow liquid to get into the product, to prevent electric shock or damage to this product.
- Do not pull or drag the power cord when unplugging. Please be gentle.
- Do not damage the wires or change the circuit of this product.
- Do not use wet cloths to clean the electric parts such as the switch and plug. Please clean this product with dry cloths. Do not use thinner, benzene or alcohol.
- Be away from this product and turn off the plug during a power cut, to prevent injury in case power resumes suddenly.
- Stop using this product if there is any irregularity, and please contact us immediately.
- Stop using this product if you feel any physical abnormality. Please consult a doctor.

## Safety Precautions

- This product is not intended for people (including children) who are physically handicapped, sensory or neurological deficits, or lack of experience and knowledge, unless they are under supervision or instruction concerning the use of the product by their guardian.
- Children should be supervised to ensure they do not play it as a toy.
- Do not directly pull or drag the power cord when moving the product to avoid damage.
- This product should be maintained by us, users are prohibited to dismantle, repair or maintain by themselves.
- Keep away from fire, high temperature or flammable, to avoid malfunctions and dangers that may cause deformation, discoloration or damage.
- Do not use other therapy devices while using this product.
- This product is a personal care massager and cannot replace medicines and medical devices.
- Please strictly follow the instructions, otherwise an accident or physical discomfort may occur.


## 5. Maintenance


- This product should be maintained by us, users are prohibited to dismantle or repair by themselves.
- Please clean this product with dry cloths. Do not use thinner, benzene or alcohol.
- Keep the product clean and away from the places with high humidity or high temperature and environment with volatile gases.
- Do not place it near fire or direct sunlight, to prevent product color fading.

 Name and Function of Components


Product structure

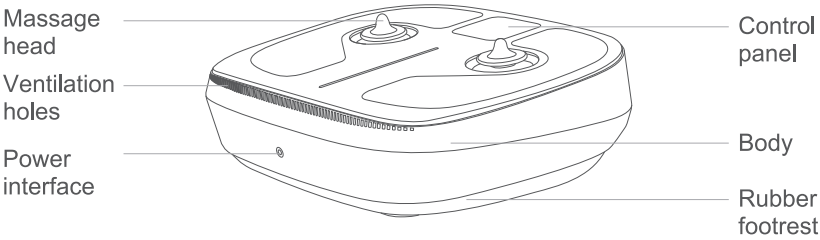
Components

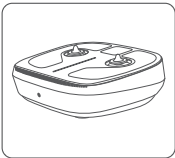Foot massager

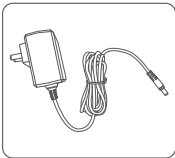Adapter

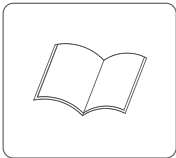User manual

## HOW TO USE

Connect to the power supply, the product turn into standby mode and the indicator light flash.

Function	Icon	Display	Instructions
Power Button			Press and hold this button for two seconds to enter standby mode. Press again to turn off the machine. When the machine is shut down, there is no display on the screen.
Intensity Button			After entering standby mode, press this button to start the massage function with a countdown of 3 seconds on the display screen. The "AU" will be displayed for two seconds.
			There are 3 adjustable intensity levels. The default is level 2 (S2). Press again to switch to level 3 (S3), and press again to switch to level 1 (S1), This cycle repeats.
Heat Button			Press this button to turn on the heat function at level 1 (H1), Press again to switch to level 2 (H2), and press again to turn off the heat (H0). This cycle repeats.
Timer Button			After starting the massage, the default working time is 15 minutes. Press this button to select a massage time of 15 or 30 minutes.

### High Blood Pressure Improvement Plan

Foot reflexology plan: The following are reflex zones/acupressure points (massage in numerical order)

- |                                           |                                                      |
|-------------------------------------------|------------------------------------------------------|
| 1. Head (brain) - kneading for 2 minutes  | 2. Cerebellum and brainstem - pressing for 3 minutes |
| 3. Adrenal gland - pressing for 2 minutes | 4. Kidney - pressing for 2 minutes                   |
| 5. Ureter - pushing for 3 minutes         | 6. Bladder - pressing for 1 minute                   |
| 7. Heart -kneading for 1 minute           | 8. Liver - kneading for 2 minutes                    |
| 9. Neck - pressing for 1 minute           |                                                      |

#### 1.Reasonable diet

Reduce salt, oil, and MSG intake, control staple foods, and limit high-calorie sugary drinks. Consume less preserved foods that are high in salt.

#### 2.Psychological balance

Maintain an optimistic and positive attitude. Avoid negative emotions and actively face life. Avoid anxiety, tension, and staying up late, and establish a regular daily routine.

#### 3.Moderate exercise


Maintain an optimistic and positive attitude. Avoid negative emotions and actively face life. Avoid anxiety, tension, and staying up late, and establish a regular daily routine.

## How to Use

### Tips for Use

#### Step1:

Insert the power cord into the power input interface at the back of the product, and then plug the adapter into the power socket. The display screen will flash for 2 seconds.
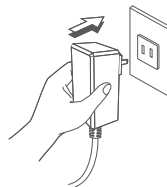
#### Step2:

Press and hold the "power button" for 2 seconds to enter standby mode, and the "power button" will remain on. Press the "intensity button" to start the massage function after a countdown of 3 seconds. The default massage intensity is level 2, and there are a total of 3 adjustable intensity levels.

#### Step3:

Press the "heat button" to choose whether to use the heating function on the massage head. There are two temperature levels that can be adjusted. Press the "timing button" to choose between 15 or 30 minutes of timing. The device will automatically shut down after 15 minutes by default, and there will be a 30-second countdown reminder before shutdown.

#### Step4:


After pressing and holding the "power button" for 2 seconds or automatic shutdown, clean the main unit with a damp cloth or dry cloth. (Remember to turn off the power before cleaning the main unit.)


## How to Use

### Preparation before use

1. Please place this product on a flat surface, 10cm away from the seat.

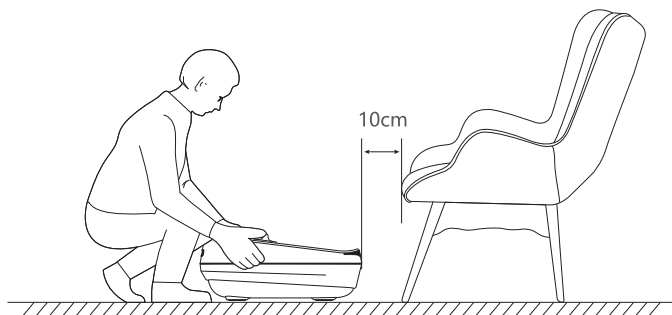2. Please maintain the correct sitting posture, place both feet on the massage heads, and bend both legs in a comfortable position.
3. Keep your feet dry when using.
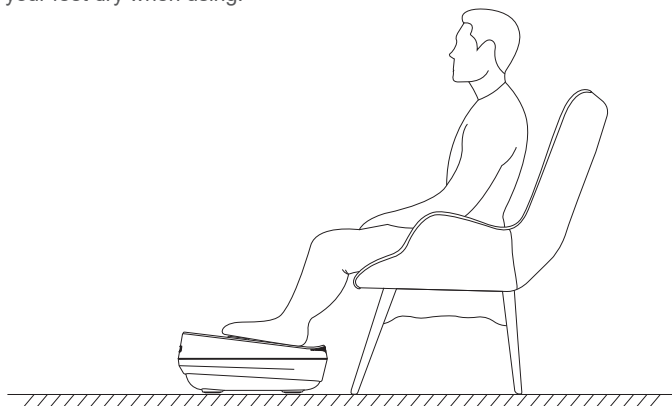
## Cleaning and Maintenance

### Cleaning

- Please ensure to turn off the power and unplug before cleaning.
- Use a soft dry cloth or slightly damp cloth to wipe the product surface.

- Do not touch the plug with wet hands.
- Do not wash or spill any liquid on the product.
- Do not use corrosive liquid such as benzene or thinner to clean the product, to prevent color fading and damage.
- When the product is damp, stop using it and do not dismantle. Leave it to dry naturally.

### Maintenance

- Keep the product clean and away from the places with high humidity or high temperature and environment with volatile gases.
- If the product will not be used for a long period, please use a dust cover (buy separately) to prevent dust, then store in a safe, dry and cool place.

- Do not place it near fire or direct sunlight, to prevent product color fading.
- Do not place heavy objects on top of the product.
- Please ensure to turn off the power and unplug when it's not in use.

### Kindly reminding

- If the product has been used for a long time, please check the conditions as below:
  - \* Burnt smell.
  - \* Unstable power supply after connecting with the power cord.
  - \* Power cord overheat
  - \* Power cords damage
  - \* Do not pull the cloth cover and the zipper forcefully to avoid damage.
  - \* Other abnormalities.

- If the above phenomenon occurs, please stop using the product to prevent malfunctions and accidents.
- Please turn off the power and unplug immediately. Users are prohibited to dismantle and repair by themselves. Please contact us for repair.

## Troubleshooting

The following symptoms do not indicate malfunction. Please check the below phenomenon to find out the solution before contact customer services department. If the problem still cannot be solved, please do not dismantle or repair the product. You should stop using it and contact us immediately.

Issues	Possible Causes	Solutions
"Tik tik" or "Tak tak" sound during operation	Normal mechanical sound during operation or the friction sound between the massage balls and cloth cover	No action required
Not working	The adapter has not been well connected	Connect the adapter properly
	The power switch is off	Switch on the power button
	The massager is used for too long, overheat protection will activated	Allow the massager to cool down before using again
Abrupt stop during operation	Power disconnected	Check if the power cord is connected
	Rated working time (15 minutes) is up and the massager is automatically stop	Press the power button to restart
	If you have been using the massager for a long time, overheat protection may get activated.	Turn the massager off and allow it to cool down before using again
	The massage intensity is too strong, stall protection may get activated	Change the mode and intensity setting, to relax your feet and reduce the machine load

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Trademark holder: Maxcare Industrial Limited