
3D SHOULDER MASSAGE

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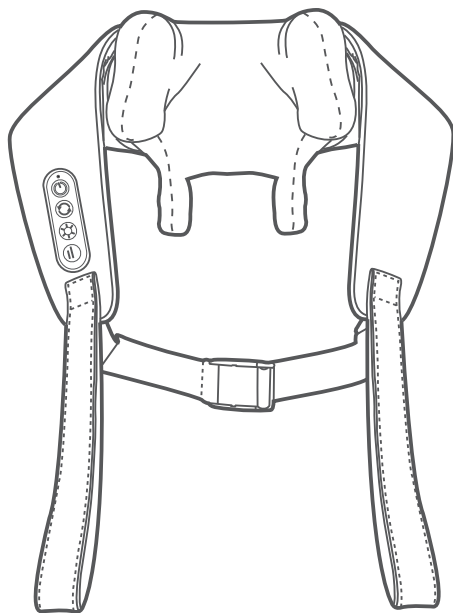MAX-E02EMK158HAB(藍色)

MAX-E02EMK158HAG(綠色)

maxcare®

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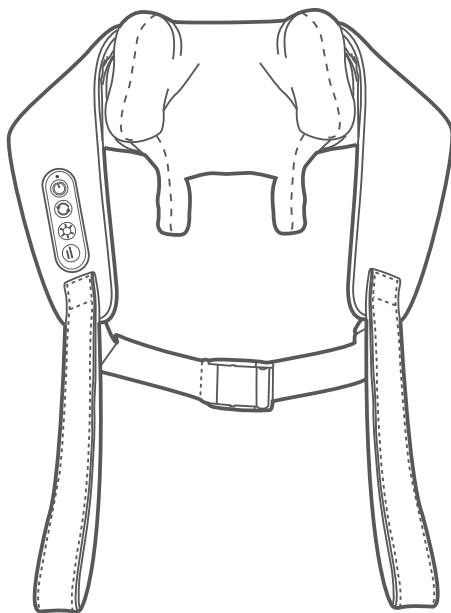
maxcare®

COMPANY PROFILE

After two decades of experience in the health industry, **maxcare** founder utilized his knowledge and expertise in massage therapy from an Asian-Pacific perspective to establish our brand - **maxcare**. Our massage therapy business was expanded to Hong Kong in 2012 and is conferred 2017 Hong Kong Top Brand Award. Based in Hong Kong, we are specializing in design & innovation, manufacturing, and retailing of massage and healthcare products.

Believing in the value of 'Technology improves our quality of life' and 'Caring for your family', **maxcare** has been innovating highly qualified products with international patents. Besides, **maxcare** products fulfill various international safety standards, namely EU CE Mark, US UL/ETL Mark, Germany GS Mark, TUV certification, and China CCC Mark.

Motivated by our value, **maxcare** founder has been leading us to be one of the most competitive brands in Hong Kong health industry. **maxcare** applies the contemporary fashion design – neat and modern – on the mechanical engineering of massage chairs. To deliver our value and offer incredible massage experience and healthcare solution to the public, **maxcare** devotes to the innovation, design, manufacturing the top-notch family massage products. By integrating the company philosophy – Spirit Of Healthy Life - into our products, **maxcare** brings health to everyone.

3D SHOULDER MASSAGE

Thank you for purchasing the Maxcare products. Before using the product, please read this manual thoroughly, particularly the Safety Precautions section, for the correct use of the product.

This manual should be well kept for easy reference. It enables you to perform a safe and efficient operation.

CONTENT

Safety Precautions	P.01 - P.03
Name and Function of Components	P.04
How to Use	P.04-05
Cleaning and Maintenance	P.05 - P.06
Troubleshooting	P.07
Technical Specifications	P.08

Safety Precautions

1. Important safety warning

- For adults only.
- Do not put this product into water or let water get into it.
- Do not drop or misuse this product.
- Please use well-grounding power supply suited to this product.
- Do not leave this product unattended when working or charging.
- Keep children and pets away from the movable parts of this product.
- When in use, please ensure that there are no people or pets nearby to avoid accidents.
- Please pull out the charging cable before cleaning to avoid injury or damage to the product.
- Do not use unauthorized accessories.
- Do not use this product outdoors.
- Please read this manual instruction carefully before operation.
- Any other usage that is not listed in this manual instruction is forbidden.
- Recommended usage time of this product is 10 minutes / session.
- Please check if any part of this products is damaged before use. If so, please do not use it and contact us for repair.
- Please do not use it if the cover is broken or damaged.
- Please do not use it if the shell or cloth cover is damaged.
- Please do not drop anything into this product.
- Please do not fall asleep while using this product.
- Please do not use it while driving
- Please do not use it if drunk or feeling unwell.
- Please do not use it within one hour after having meal.
- Please do not adjust to the maximum intensity at the beginning of the massage to avoid possible injury.
- Keep the plug clean from needles, rubbish or water.
- If the power cord is damaged, do not attempt to change on your own. Please contact us for repair for safety reasons.

2. Environment for usage

- Please do not use it in a moist environment such as bathrooms.
- Please do not use it when there is a drastic change in temperature.
- Please do not use it under heavy-dusty or caustic environment.
- Please do not use it where there is not enough space or not good ventilation.
- Please switch off the power and unplug the machine after use.
- Please do not use it when surrounding temperature is up to 40°C.
- Please do not use it when surrounding temperature is below 5°C. The lubricant can thicken quickly and increase the friction between parts.
- Turn off and pull out the plug during thunderstorms.
- Do not place it near fire or directly under the sun.
- When moved from a low-temperature storage to a warm environment, the moisture inside the machine may affect, or damage, its functionality. Please rest the product within room temperature for 1 hour before use.

Safety Precautions

3. People who have the below symptoms, please consult a doctor before use

- People who are suffering from skin conditions such as swollen skins, sunburnt or other dermatology diseases.
- People who are suffering from malignant tumors, acute diseases, heart diseases, stroke or severe hypertension.
- People who are suffering from any difficulties related to respiratory organs.
- People who are suffering from any blood circulation disorders such as thrombosis, severe aneurysms and venous tumors, or from any skin infections.
- People who are suffering from perception disorders due to diabetes.
- People who are suffering from osteoporosis.
- People who are suffering from heart disease or carrying electronic medical equipment such as heart pacemaker.
- People who have a fever.
- Pregnant women or women that are in menses period.
- People who are suffering from injuries or skin diseases.
- Children under 14-year-old or people mentally unbalanced without being supervised are not allowed to use this machine.
- People who are told by doctor to have rest or people who feel unwell.
- People whose bodies are wet are forbidden to use this product.
- People who are receiving medical treatment, especially those who are physically unfit.

4. Safety

- Please do not place heavy objects on this product and avoid collision with other objects.
- Do not use other therapy devices while using this product.
- This product is not to be used for medical treatment purposes.
- Only one person at a time.
- Recommended usage time of this product is 10 minutes / session.
- Check if the voltage is suitable to the specification of this product.
- Do not use wet hands to pull the plug.
- Do not allow water to get into the product, to avoid electric shock or damage to this product.
- Do not pull or drag the power cord when unplugging. Please be gentle.
- Do not damage the wires or change the circuit of this product.
- Do not use wet cloths to clean the electric parts such as the switch and plug.
- If there is a power outage, the charging of the product should be suspended to avoid damage when the power is suddenly restored.
- Stop using this product if there appears any irregularity, and please contact us immediately.
- Stop using this product if it causes discomfort, and please consult a doctor.
- This appliance is not intended for people (including children) with reduced physical, sensory or mental capabilities, or a lack of experience and knowledge, unless they are under supervision or instruction concerning the use of the appliance by a person responsible for their safety.
- If the power cord is damaged, it must be replaced by us to avoid a hazard.
- Keep the power cord away from hot or moving objects.
- If the product is found damaged, broken, suffering from power leakage, or with exposed parts, please stop using it.
- Do not turn on the product when it is being charged.

Safety Precautions

5. Maintenance


- This product should be maintained by us, users are prohibited to disassemble or maintain by themselves.
- Please always turn off the product after use.
- If this product is to be left idle for long, please tie up the power cords and keep it in dry and dustless environment.
- Do not keep this product at a high temperature or near fire. Do not keep under direct sunlight.
- Please clean this product with dry cloths. Do not use thinner, benzene or alcohol.
- The mechanical components of this product are specially designed and made, no special maintenance required.
- Do not point sharp objects at this product.
- Please use it intermittently; do not keep it continuously running for long.

6. Solution of ordinary malfunctions

- It is normal if the motor makes noises during operation.

⚙️ Name and Function of Components

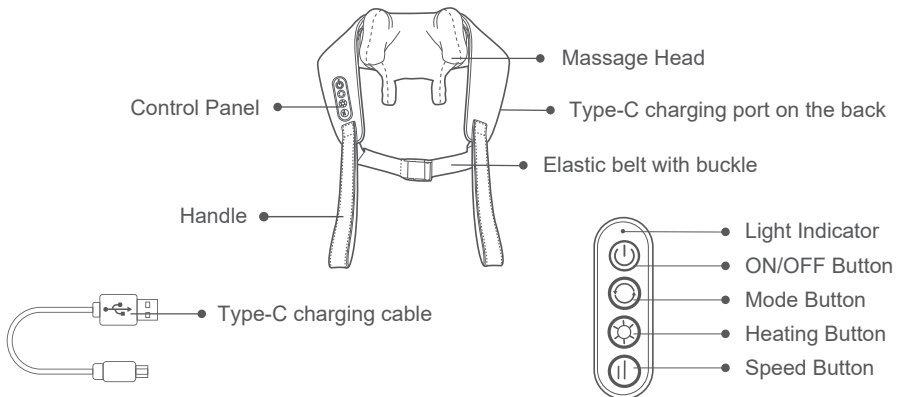⚙️ How to Use

Neck

Back

Abdomen

Thigh

Waist

Soles

🖐️ Operation Instruction

Button	Instruction	Beep sound
ON/OFF	Press and hold this button for 2 seconds to turn on and start the countdown. The default working time is 10 minutes. When the time is up, the product will automatically turn off; you can also press and hold this button for 2 seconds to manually turn it off.	Turn on : Beep once Turn off : Beep once
Mode	The default mode is Shiatsu mode, press this button once to switch to kneading mode, press it again to switch back to Shiatsu mode, and repeat this sequence.	Shiatsu mode: Beep once Kneading mode: Beep twice
Heating	The preset Heating mode is High Heating, press this button once to switch to Low Heating, press it again to turn off heating, and press it again to switch back to High Heating, and repeat this sequence.	Low Heating: Beep once High Heating: Beep twice Turn off the heating function: No beep
Speed	The default speed is Low speed, press this button once to switch to High speed, press it again to switch back to low speed, and repeat this sequence.	Low Speed: Beep once High speed: Beep twice

Charging Instruction

- 1.The product cannot be turned on while charging. Before starting the massager, you need to unplug it.
- 2.When the battery power is more than 10.5V, the indicator light is white.
- 3.When charging, the indicator light flashes red slowly.
- 4.When fully charged, the indicator light will flash white slowly.
- 5.When the battery power is lower than 10.5V, the indicator light will be red,indicating that the product needs to be charged in time.

CAUTION

- Please ensure the product is turned off before charging.
- If the red light is still on (during charging state) after continuous charging for more than 6 hours, please unplug the charger to avoid danger.

Cleaning and Maintenance

Cleaning

- Before cleaning, please ensure the switch has been turned off and the charging cable has been pulled out
- Use soft cloths to clean the surface of product and control panel.
- Use lightly damp cloths to clean and then dry it.

CAUTION

- Do not pull the plug from the socket with wet hands.
- Please do not rinse the product. Avoid water/ liquids from getting into the product.
- Do not use corrosive cleaning agents like benzene or thinner to clean the product, otherwise, the color of product will fade, and the internal parts may be corrupted or broken.
- When the product is wet, do not disassemble or use the product, just let it dry.

 **Cleaning and Maintenance****Maintenance**

- Keep the product clean and away from the places with high humidity or high temperature and environment with volatile gases.
- If the product will not be used for a long period, please use a dust cover (buy separately) to prevent dust, then store in a safe, dry and cool place.

**CAUTION**

- Do not place it near fire or direct sunlight, to prevent product color fading.

Kind reminder

- If the product has been used for a long time, please check the conditions as below:
 - * Burnt smell.
 - * Unstable power supply after connecting with the power cable.
 - * Other abnormalities.

**CAUTION**

- If the above phenomenon occurs, please stop using the product to prevent malfunctions and accidents.
- Please turn off the power and unplug immediately. Users are prohibited to dismantle and repair by themselves. Please contact us for repair.

Troubleshooting

The following symptoms do not indicate malfunction. Please check the below phenomenon to find out the solution before contact customer services team. If the problem still cannot be solved, please do not dismantle or repair the product. You should stop using it and contact us immediately.

Issues	Possible Causes	Solutions
Not working	Low battery	Recharge the product
	Haven't press the power button	Press and hold the ON/OFF button to start
	Rated working time (10 minutes) is up and the massager is automatically stop	Press the power button to restart
Abrupt stop during operation	Low battery	Recharge the product
	If you have been using the massager for a long time, overheat protection may get activated	Turn the massager off and allow it to cool down for 30-50 mins before using it again

 **Technical Specifications**

Product Name	: 3D Shoulder Massage
Product Model	: MAX-E02EMK158HAB (Blue) MAX-E02EMK158HAG (Green)
Rated Voltage	: DC 5V 
Rated Power	: 30W
Battery capacity	: 1800mAh
Charging Time	: 3Hrs
Battery Life	: 90minutes
Heating temperature	: 40-47°C
Rated Time	: 10mins
Net Weight	: 1.25 kg
Gross Weight	: 1.7 kg
Packaging Size	: 245 x 170 x 163 mm
Components	: Massager, charging cable, user manual

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Trademark holder: Maxcare Industrial Limited