# EMS CIRCULATION BOOSTER

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MAX- M01MC0300G

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MAX- M01MC0300G

## COMPANY PROFILE

After two decades of experience in the health industry, **maxcare** founder utilized his knowledge and expertise in massage therapy from an Asian-Pacific perspective to establish our brand - **maxcare** . Our massage therapy business was expanded to Hong Kong in 2012 and is conferred 2017 Hong Kong Top Brand Award. Based in Hong Kong, we are specializing in design & innovation, manufacturing, and retailing of massage and healthcare products.

Believing in the value of 'Technology improves our quality of life' and 'Caring for your family', **maxcare** has been innovating highly qualified products with international patents. Besides, **maxcare** products fulfill various international safety standards, namely EU CE Mark, US UL/ETL Mark, Germany GS Mark, TUV certification, and China CCC Mark.

Motivated by our value, **maxcare** founder has been leading us to be one of the most competitive brands in Hong Kong health industry. **maxcare** applies the contemporary fashion design – neat and modern – on the mechanical engineering of massage chairs. To deliver our value and offer incredible massage experience and healthcare solution to the public, **maxcare** devotes to the innovation, design, manufacturing the top-notch family massage products. By integrating the company philosophy –Spirit of Healthy Life - into our products, **maxcare** brings health to everyone.

Thank you for purchasing Maxcare's products. Before using the product, please read this manual thoroughly, particularly the Safety Precautions section, for the correct use of the product.

This manual should be well kept for easy reference. It enables you to perform a safe and efficient operation.

## CONTENT

---

Safety Precautions .....	P.01 - P.03
Name and Function of Components .....	P.04 - P.05
How to Use .....	P.06 - P.07
Cleaning and Maintenance .....	P.07
Troubleshooting .....	P.08
Technical Specifications .....	P.09

## Safety Precautions

### 1. Important safety warning

- Keep children and pets away from the movable parts of this product.
- Please use well-grounding power supply suited to this product.
- Please pull out the plug after use or before cleaning to avoid injury or damage to the product.
- Do not use unauthorized accessories.
- Do not use this product outdoors.
- Please read this manual instruction carefully before operation.
- Any other usage that is not listed in this manual instruction is forbidden.
- Recommended usage time of this product is 30 minutes/ session.
- Please check if any part of this products is damaged before use. If so, please do not use it, and contact us for repair.
- Please do not use it if the cover, shell or surface is damaged.
- Please do not drop anything into this product.
- Please do not fall asleep while using this product.
- Please do not use it while driving.
- Please do not use it if drunk or feeling unwell.
- Please do not use it within one hour after having meal.
- Please do not adjust to the maximum intensity at the beginning of the massage to avoid possible injury.
- Keep the plug clean from needles, rubbish or water.
- If the power cord is damaged, do not attempt to change on your own. Please contact us for repair for safety reasons.
- Do not stand on this product, only use when seated.
- If you want to move the Electrode Pad to other region during treatment, be sure to turn off the power first, to avoid injury or damage to the product.
- Do not start treatment while wearing an electronic device.

### 2. Environment for usage

- Please do not use it in a moist environment such as bathrooms.
- Please do not use it when there is a drastic change in temperature.
- Please do not use it under heavy-dusty or caustic environment.
- Please do not use it where there is not enough space or not good ventilation.
- Please switch off the power and unplug the machine after use.
- Do not use outdoors.
- Do not place it near fire or direct sunlight.
- When moved from a low-temperature storage to a warm environment, the moisture inside the machine may affect, or damage, its functionality. Please rest the product within room temperature for 1 hour before use.

## Safety Precautions

### 3. People who have the below symptoms, please consult a doctor before use

- People who are suffering from skin conditions such as swollen skins, sunburnt or other dermatology diseases and infectious disease.
  - People who are suffering from malignant tumors, acute diseases, heart diseases, stroke or severe hypertension.
  - People who are suffering from any difficulties related to respiratory organs.
  - People who are suffering from abnormal blood pressure, or any blood circulation disorders such as thrombosis, severe aneurysms and venous tumors, or from any skin infections.
  - People who are suffering from perception disorders due to diabetes.
  - People who have a fever.
  - People who are suffering from injuries or skin diseases.
  - Children under 14-year old or people mentally unbalanced without being supervised are not allowed to use this machine.
  - People who are pregnant.
  - People who are told by doctor to have rest or people who feel unwell.
  - People whose bodies are wet are forbidden to use this product.
  - People who are receiving medical treatment, especially those who are physically unfit.
  - People who have open wounds.
  - People who aim for recovery.
  - This unit must not be used in combination with the following medical devices:
    - Internally transplanted electronic medical devices, e.g., pacemakers.
    - Electronic life support equipment, such as respirators.
    - Electronic medical devices attached to the body, such as electrocardiographs.
- Using this unit with other electronic medical devices may cause erroneous operation of those devices.

### 4. Safety

- Do not use this unit near the heart, above the neck, on the head, around the mouth.
- Please do not place heavy objects on this product.
- Do not use other therapy devices while using this product.
- This product is not to be used for medical treatment purposes.
- Only one person at a time.
- Recommended usage time of this product is 30 minutes/ session.
- Check if the voltage is suitable to the specification of this product.
- Do not use wet hands to pull the plug.
- Do not allow water to get into the product, to avoid electric shock or damage to this product.
- Do not pull or drag the power cord when unplugging. Please be gentle.
- Do not damage the wires or change the circuit of this product.
- Do not use wet cloths to clean the electric parts such as the switch and plug.

## Safety Precautions

- Be away from this product during power cut to prevent injury in case power resumes suddenly.
- Stop using this product if there appears any irregularity, and please contact us immediately.
- Stop using this product if it causes discomfort, and please consult a doctor.
- This appliance is not intended for people (including children) with reduced physical, sensory or mental capabilities, or a lack of experience and knowledge, unless they are under supervision or instruction concerning the use of the appliance by a person responsible for their safety.
- Children should be supervised to ensure that they do not play with the appliance.
- If the power cord is damaged, it must be replaced by us in order to avoid a hazard.
- Keep the power cord away from hot or moving objects.
- Please cut off power and pull out plug before moving.
- If the product is found damaged, broken, suffering from power leakage, or with exposed parts, please stop using it.
- Do not leave the product unattended when it is operating.

## 5. Maintenance

- This product should be maintained by us, users are prohibited to disassemble or maintain by themselves.
- Please always turn off the main power after use.
- Do not use this product if the socket is loose.
- If this product is to be left idle for long, please tie up the power cords and keep it in dry and dustless environment.
- Do not keep this product at a high temperature or near fire. Do not keep under direct sunlight.
- Please clean this product with dry cloths. Do not use thinner, benzene or alcohol.
- The mechanical components of this product are specially designed and made, no special maintenance required.
- Do not point sharp objects at this product.
- Please use it intermittently; do not keep it continuously running for long.

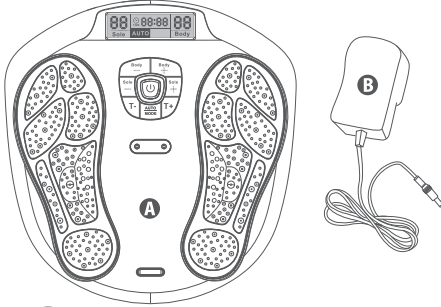
## 6. Solution of ordinary malfunctions

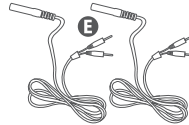
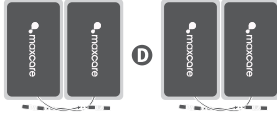
- If the rated working time is over, the power will automatically turn off. If this product works continuously for too long, the temperature control will automatically turn off the machine, if so, please allow a 30-minute break before resuming.

 Name and Function of Components

Product structure

- A. Main Unit
- B. DC Power Adapter
- C. Remote Control
- D. Electrode Gel Pad x 4 pcs.
- E. Wire cable connecting the Electrode gel pads and device
- F. A storage bag

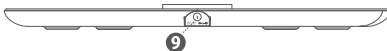
Front View

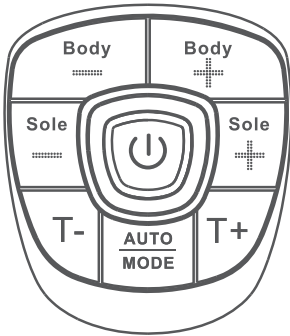
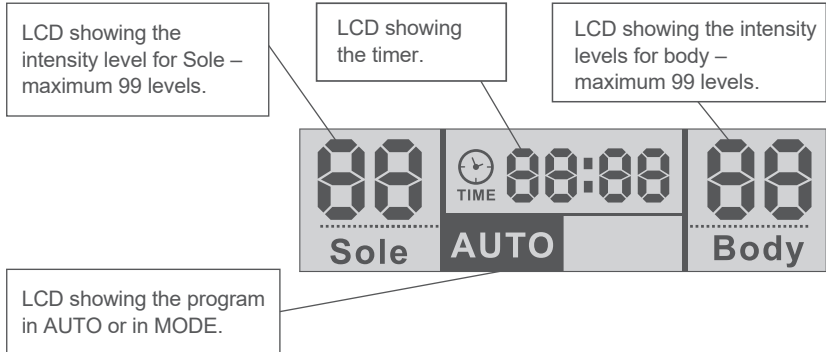
- 1. LCD display screen
- 2. Control Panel
- 3. Electro Pedal for Left Foot
- 4. Electro Pedal for Right Foot
- 5. Gel pad Jack A
- 6. Gel pad Jack B
- 7. Remote control receiver sensor
- 8. Decoration nod
- 9. AC jack for the adapter

Side View

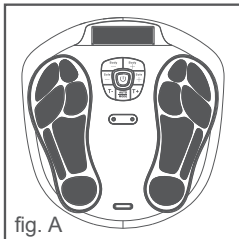
 Name and Function of Components

**Function**

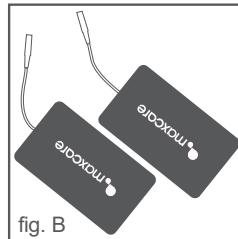
	ON/OFF switch button
<b>BODY -</b>	Decrease the output intensity of body (Available from 1 - 99 levels)
<b>BODY +</b>	Increase the output intensity of body (Available from 1 - 99 levels)
<b>SOLE -</b>	Decrease the output intensity of Sole (Available from 1 - 99 levels)
<b>SOLE +</b>	Increase the output intensity of sole (Available from 1 - 99 levels)
<b>TIME -</b>	Decrease the operation time (available from 1-60 minutes)
<b>TIME +</b>	Increase the operation time (available from 1-60 minutes)
<b>AUTO</b> <b>MODE</b>	Auto - is the preset program with 14 patterns in cycle running for Foot and 10 patterns in cycle running for body. Mode - user can fix the program to the exiting massage pattern on the rest of the time

Electrode area of the unit and gel pad:

On the device the black color area on the device which is the electrode area for the sole.  
(See. the fig. A)

On the gel pad, the black color area on the sticky part is the electrode area for the body, size is 5cm x 9 cm.  
(see. the fig. B)

 How to Use

HOW TO USE

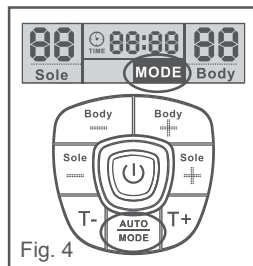
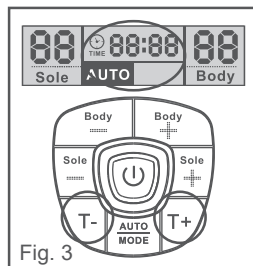
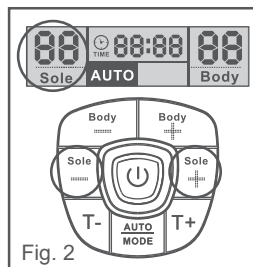
FOR FOOT - SOLE

1. Place both bare feet onto EMS Circulation Booster (do not wear socks).
2. Press the on/off button, the LCD screen will light up in White. And the program will show AUTO and both numbers will show 00, meaning the machine is in standby mode. (See Fig. 1)
3. Gently increase intensity setting by pushing the button of "SOLE +". Or decrease intensity setting by pushing the button of "SOLE-". The intensity level is adjustable between 0 and 99. The LCD display will show the selected level. (See Fig. 2)
4. You can adjust the auto off timer by pressing the "Time -" or "Time +". Timers range from 1-60 minutes. To terminate the massage period, user can turn off the unit anytime by pressing the on/off button once. The timer will start to count down from your chosen time setting. (See Fig.3)
5. If the current massage program you are using is satisfactory, then you can lock the program in place by pressing the Auto/Mode key. This will make the rest of the current massage timer run only on the selected setting. (See. Fig. 4)

 **Attention:**

- A. The aim is not to get up to level '99'.
- B. Chose an intensity level that is comfortable for you!  
This level may vary from day to day.
- C. Drink adequate water to ensure the effectiveness of the device.
- D. If you have dry feet- moisturize them generously with water or moisturizer for optimal results.

EMS Circulation Booster is designed to be used at any time before or after exercise, either to improve circulation before working out, or to aid with the recovery process afterwards. For best results and optimal recovery, use for at least 30 minutes within 8 hours after exercising.

## How to Use

### HOW TO USE

#### FOR BODY

##### USING THE GEL PADS

Wash and dry skin before use.

1. Plug in the 2 cables into the cable jack on the unit. (See Fig.5)
2. Connect pin of the cable to the gel pad properly. (See Fig.6)
3. Remove the protective film on the gel pad. All four gel pads to the area of the body you wish to treat.
4. Press the on/off button to turn on the unit and adjust the EMS output intensity to the desired level for use on the body. Gently increase intensity setting by pushing the button of "Body +". Or decrease intensity setting by pushing the button of "Body -". The display will show the mode and level that you selected and start to count down. (see Fig. 7).

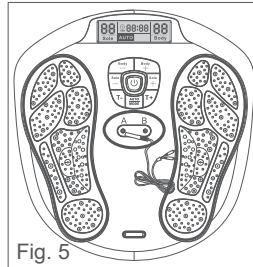
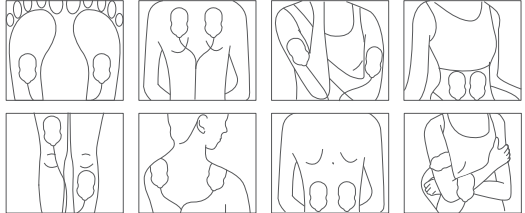
Fig. 5

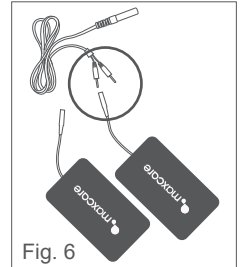
Fig. 6

##### CARE OF YOUR GEL PADS

Never stick two adhesive pads to each other. Keep the adhesive gel pads clean, never expose them to high temperature or direct sunshine. If the electrode gel pads are insufficiently adhesive or dirty, wipe with a wet cloth or change for new ones.

**!** Do not clean the electrode gel pads with any chemical.

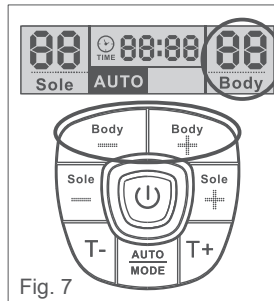
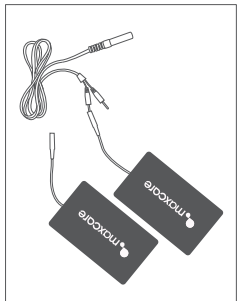
Fig. 7

## Cleaning and Maintenance

##### Main Device & Gel Pad

- Always keep the main device clean by using a piece of dry and soft cloth to wipe the surface.
- Do not spill liquid on the device or immerse it in water.
- Keep it in a cool and dry place.
- Do not clean with chemicals.
- Keep out of the reach of children and pets.

## ? Troubleshooting

The following symptoms do not indicate malfunction. Please check the below phenomenon to find out the solution before contact customer services department. If the problem still cannot be solved, please do not dismantle or repair the product. You should stop using it and contact us immediately.

Issues	Possible Causes	Solutions
Device will not turn on.	<ul style="list-style-type: none"> <li>- Batteries inserted in wrong direction.</li> <li>- The adapter is not plugged well into device properly.</li> </ul>	<ul style="list-style-type: none"> <li>- Insert batteries in correct direction or check the battery is fully powered.</li> <li>- Check the connection of the adapter jack is in well connected.</li> </ul>
Power turns off too soon.	Gel pads not attached correctly to the skin.	Attach Gel pads correctly to the skin.
Power turns off while using massager.	if you are using the batteries operation, then the batteries could be weak/exhausted.	Fit 4 new identical 1.5V alkaline batteries type AA
	Treatment period of 30 minutes is over and power turns off automatically.	Restart treatment
	if you are using the body massage, the electrode gel pad may break.	Replace electrode gel pad
It is difficult to attach Gel pad to the skin.	Transparent film not peeled off	Peel off film on the adhesive surface of Gel pad
	Gel pad applied immediately after washing	Sufficiently dry Gel pad
	Adhesive surface of Gel pad damaged	Replace Gel pad
	The gel pads get dirty and lost their adhesive/stickiness	Replace Gel pad or clean the gel pad with a small drop of water onto the sticky side of the electrode pad and rub into the surface
Adhesive surface of Gel pad is sticky.	Use of Gel pad during perspiring	Leave Gel pad in freezer for overnight
	Gel pad washed too long and/or too frequently	
	Gel pads stored under high temperature, high humidity, direct sunshine	
It is difficult to feel stimulation.	The soles of the feet are too dry and need to be moisturized.	Apply some cream or water to the feet to ensure they are adequately moisturized.
	The sole of the foot is not placed on the foot pad properly	Ensure both of your feet are placed on each pedal properly.
	Gel pads not attached correctly to the skin	Attach Gel pad firmly to the skin
	Gel pads overlap each other	Reattach Long Life pads with no overlap
	Electrode cord not connected correctly	Connect electrode cord correctly
	Applied intensity too weak	Increase the intensity by turning the dial switch
The skin turns red or the skin feels irritated.	Adhesive surface of Gel pads dirty or dry	Wash adhesive surface of Gel pads softly with your fingertips for about 3 seconds under slow running water
	Adhesive surface of Gel pads damaged	Replace Gel pads

 **Technical Specifications**

- Product Name : EMS Circulation Booster
- Model : MAX-M01MC0300G
- Power supply : Main unit: 6V DC or 4x1.5V alkaline batteries type AA\*  
: Remote control: 2x1.5V alkaline batteries type AAA\*
- Adapter Input : AC 100-240V~50-60Hz 0.1A
- Adapter Output : DC 6V 500mA3.0W
- Battery life : >350 minutes
- Frequency generation : Approx. 10 Hz to 55.56 Hz
- Power consumption : 1.05W
- Operating temperature and humidity : -5°C to 50°C, 30%to 90% RH
- Storage temperature and humidity : -10°C to 60°C, 10% to 95% RH
- Main unit dimensions : 338(L) x 324 (W) x 48(H) mm
- Weight Approx. : 950 g
- Package Contents :

Quantity	Parts
1	EMS Circulation Booster
1	AC/DC Adapter
1	Remote Control
2	Cable Wire for Electrode Gel pads
4	Electrode Gel pads
1	Instruction manual

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Trademark holder: Maxcare Industrial Limited