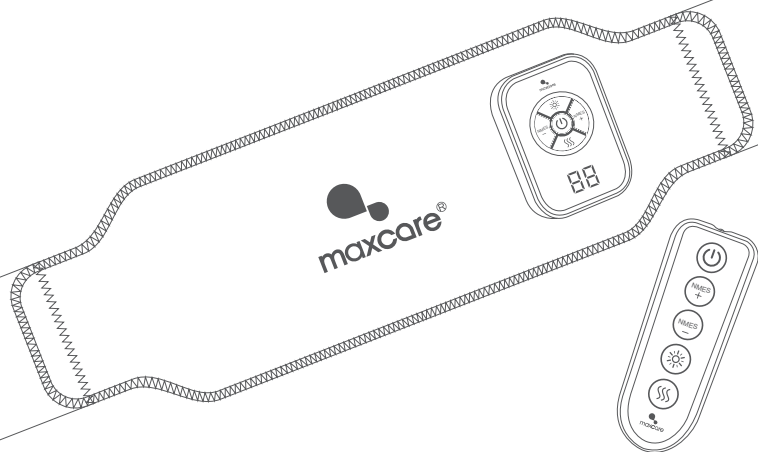
Phototherapy NMES
光療脈衝復康帶

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MAX-M01MC0577B

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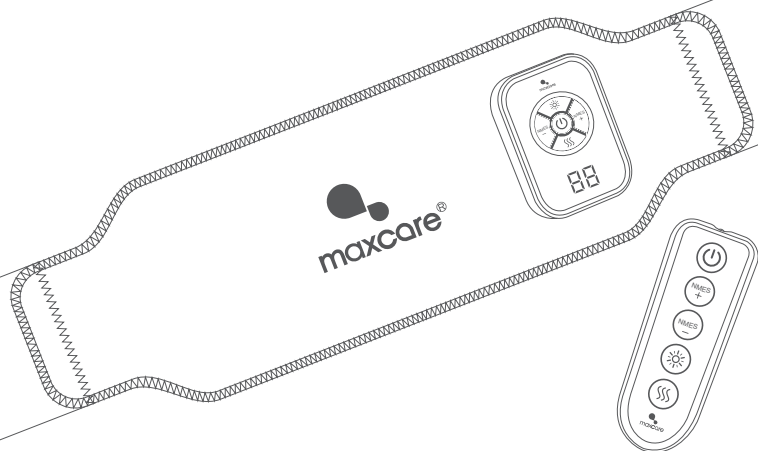
Phototherapy NMES Rehabilitation Belt

User Manual

MAX-M01MC0577B

This product jointly developed by MTCURE (镁天健康) and The Chinese University of Hong Kong, and maxcare provides engineering and manufacturing supports, effectively improves pain conditions while promoting metabolism and blood circulation.

MTCURE (镁天健康), an innovative biotechnology start-up incubated by The Chinese University of Hong Kong.

COMPANY PROFILE

After two decades of experience in the health industry, **maxcare** founder utilized his knowledge and expertise in massage therapy from an Asian-Pacific perspective to establish our brand - **maxcare** . Our massage therapy business was expanded to Hong Kong in 2012 and is conferred 2017 Hong Kong Top Brand Award. Based in Hong Kong, we are specializing in design & innovation, manufacturing, and retailing of massage and healthcare products.

Believing in the value of 'Technology improves our quality of life' and 'Caring for your family', **maxcare** has been innovating highly qualified products with international patents. Besides, **maxcare** products fulfill various international safety standards, namely EU CE Mark, US UL/ETL Mark, Germany GS Mark, TUV certification, and China CCC Mark.

Motivated by our value, **maxcare** founder has been leading us to be one of the most competitive brands in Hong Kong health industry. **maxcare** applies the contemporary fashion design – neat and modern – on the mechanical engineering of massage chairs. To deliver our value and offer incredible massage experience and healthcare solution to the public, **maxcare** devotes to the innovation, design, manufacturing the top-notch family massage products. By integrating the company philosophy –Spirit of Healthy Life - into our products, **maxcare** brings health to everyone.

Phototherapy NMES Rehabilitation Belt

Thank you for purchasing Maxcare's products. Before using the product, please read this manual thoroughly, particularly the Safety Precautions section, for the correct use of the product.

This manual should be well kept for easy reference. It enables you to perform a safe and efficient operation.

CONTENT

Safety Precautions	P.01 - P.03
Product Introduction	P.04
Name and Function of Components	P.05
How to Use	P.06
Cleaning and Maintenance	P.07
Solutions and FAQs	P.08
Technical Specifications	P.09

Safety Precautions

1. Important safety warning

- Please use a well-grounding power supply suited to this product.
- Please pull out the plug after use or before cleaning to avoid injury or damage to the product.
- Do not use unauthorized accessories.
- Please read this manual instruction carefully before operation.
- Any other usages that are not listed in this manual instruction are forbidden.
- Please check if any part of this product is damaged before use. If so, please do not use it and contact us for repair.
- Please do not use it if the product is broken or damaged.
- Please do not use it if the cover or surface is damaged.
- Please do not drop anything into this product.
- Please do not allow water or liquid into this product.
- Please do not fall asleep while using this product.
- Please do not use it if drunk or feeling unwell.
- Please do not overexertion while using the product.
- Keep the plug clean from needles, rubbish, or water.
- If the power cord is damaged, do not attempt to change on your own. Please contact us for repair for safety reasons.

2. Environment for usage

- Please do not use it in a moist environment such as a bathroom.
- Please do not use it when there is a drastic change in temperature.
- Please do not use it under a heavy-dusty or caustic environment.
- Please do not use it where there is not enough space or not a good ventilation.
- Please switch off the power after use.
- Please do not use it when the surrounding temperature is up to 50°C.
- Please do not use it when the surrounding temperature is below -20°C. The lubricant can thicken quickly and increase the friction between parts.
- Turn off and pull out the plug during thunderstorms.
- Do not place it near fire or directly under the sun.
- When moved from a low-temperature storage to a warm environment, the moisture inside the machine may affect, or damage, its functionality. Please rest the product within room temperature for 1 hour before use.

Safety Precautions

3. People who have the below symptoms, please consult a doctor before use

- People who are suffering from heart disease or carrying electronic medical equipment such as heart pacemaker.
- People who have cancers, acute illnesses, heart diseases, or severe high blood pressure, tuberculosis patients, etc.
- People who have acute pain and cause to edema, bacterial infection, inflammation, fracture, sprain, muscle strain.
- People who have skin diseases, such as infections skin diseases, the skin is insensitive to heat and wounds, or patients who are allergic to light and are deemed unfit to use this product by their doctor etc.
- People who are receiving medical treatment, especially those who are physically unfit.
- People who a bleeding tendency.
- People who have a fever.
- Elderly with serious uncomfortable.
- Pregnant women, women who just give birth or are in the menses period, newborn baby.
- Children under 10-year-old or people mentally unbalanced without being supervised are not allowed to use this machine.
- Stop using this product if it causes discomfort, and please consult a doctor.
- We reserve the right of the final interpretation. If the hidden dangers caused by unreasonable use, we would not take all the consequences.

4. Precautions before using the product and some safety instructions

- Do not use the NMES and Heat functions on areas that have implanted metal or other materials, as well as areas with inflammation. It is recommended to use Red Light and Infra-Red Therapy instead and keep the product at a distance of 3-5mm from the skin surface.
- Avoid using the NMES and Heat functions on areas with wounds, pimples, or suppurative diseases. Once the wounds have stopped bleeding or the tissue sprain has healed, it is advisable to use only Red Light and Infra-Red Therapy to promote faster healing. After 72 hours or when the wound has fully healed, it is safe to use the machine.
- Do not use the machine on the heart, head, eyes, or private parts.
- Pregnant women and puerperal are not allowed to use the machine on abdominal.
- Do not directly view the light. (Included wear sunglasses or protective spectacles)
- Please do not place heavy objects on this product.
- Do not use other therapy devices while using this product.
- Please note this is not a medical device, it is not to be used for medical treatment purposes. Only use for improvement and relief purposes.
- Recommended usage time of this product is 10 minutes/ session.
- Check if the voltage is suitable for the specification of this product.
- Do not use wet hands to pull the plug and use the machine.

Safety Precautions

- Do not allow water to get into the product, for the avoidance of electric shock or damage to this product.
- Do not pull or drag the power cord when unplugging. Please be gentle.
- Do not damage the wires or change the circuit and alter the structure or parts of this product.
- Do not use wet cloths to clean the electric parts such as the switch and plug.
- Stop using this product if it appears any irregularity, and please contact us immediately.
- Stop using this product if it causes discomfort, and please consult a doctor.
- This appliance is not intended for people (including children) with reduced physical, sensory or mental capabilities, or a lack of experience and knowledge, unless they are under supervision or instruction concerning the use of the appliance by a person responsible for their safety.
- Children should be supervised to ensure that they do not play with the appliance.
- If the power cord is damaged, it must be replaced by us in order to avoid a hazard.
- Keep the power cord away from hot or moving objects.
- Please cut off power and pull out the plug before moving.
- If the product is found damaged, broken, suffering from power leakage, or with exposed parts, please stop using it.
- This product contains a lithium-ion rechargeable battery.
- Do not push or clamp down on the machine.
- Do not place this product near fire or where it causes fire easily to avoid accidents; the product may also deform, fade, or break.

5. Maintenance

- This product should be maintained by us, users are prohibited to disassemble or maintain by themselves.
- If this product is to be left idle for long, please keep it in a dry and dustless environment.
- Do not keep this product at a high temperature or near a fire. Do not keep under direct sunlight.
- Please clean this product with dry cloths. Do not use thinner, benzene, or alcohol.
- Please use wet cloths to clean the massage head after use.
- Do not point sharp objects at this product.
- Please use it intermittently; do not keep it continuously running for long.

6. Solution of ordinary malfunctions

- You will hear a "sizzling" sound when using it. This is the sound of current flowing. This is not a product fault.

 **Product Introduction**

Product structure

With the increasing pressure of life and work, people’s bodies will have various problems. Phototherapy NMES Rehabilitation Belt is a new treatment that can be used in home or anywhere. It is a safe, efficient, reliable, painless, drug-free, cordless, and natural solution for body pain relief. It can transmit the NMES, Heat and Photobiomodulation Treatment (PBMT) to the human tissue, and directly make muscles, subcutaneous tissue, nerves, etc. It can be used separately or combined at the same time.

*This is not a medical device. This unit is intended to be used only for general physical conditioning.

Usage

Suitable for shoulders, neck, arms, waist, hips, feet, elbows.

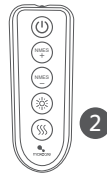
Product Function

<p>NMES</p>	<p>Electronic Muscle Stimulation helps to stimulates muscle blood flow, relief muscle fatigue, neural rehabilitation, enhance muscle strength.</p>
<p>Heat</p>	<p>Hot compress can help to relax tired and sore muscles, available with 3 levels from 38°C > 40°C > 45°C</p>
<p>MT-PBMT (MT-Photobiomodulation)</p>	<p>MT-PBMT combines light in the range of visible light and invisible light, starting at superficial tissue to deeper tissue levels, it promotes wound healing and tissue repair, relieves pain and inflammation.</p>

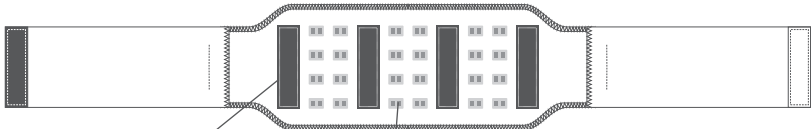
MTCURE (MTCure Biotech Limited) biophysical combined therapy can effectively enhance the effects of various biological treatment signals, thereby avoiding side effects.

 Name and Function of Components

- 1. Phototherapy NMES Rehabilitation Belt
- 2. Remote control
- 3. Extension belt
- 4. CR2032 Battery
- 5. USB-C charging cable

Controller

Heat and NMES
conducting pad

Phototherapy

Velcro Hook

Power
indicator

NMES
intensity
indicator

Charging Jack

How to Use
Operation Instruction

	Switch button - press and hold for 2 seconds to switch on/off the rehabilitation belt, this button on the remote control can only be turned off the device but cannot be turned on.
	NMES function increase intensity button - press to increase the intensity from 0- 20 levels.
	NMES function decrease intensity button - press to decrease the intensity from 20-0 levels.
	Phototherapy button - when selecting Phototherapy, this icon will light up in red.
	Heat button - when selecting heat function, this icon will light up in red for 1st Level 38°C > green 2nd Level 40°C > blue 3rd Level 45°C (+/- 5°C)

How to Use

1. Take out an extension belt and wrap the belt around your body. Then press and hold the switch button for 2 seconds to turn on.
2. The power indicator will light up in blue and the LED show 00 in blue, which means the rehabilitation belt is under the standby mode.
3. After that you can control the belt with remote control or directly press the button of controller.
4. Please choose the corresponding function according to your needs. You can use 3 functions at the same time, or only take 1 function, depending on your needs.
5. If you have selected the NMES function, suggest you apply some water-based gel or water on your skin to enhance the conductivity of your skin. And make sure the rehabilitation belt contacts the skin firmly, Otherwise, the electric impulse cannot be transmitted to your muscle.
6. Please don't use 1 position over 20 mins.
7. Low voltage warning, when the battery voltage is lower than 3.35V, the power indicating light will flash, flashing for 15 seconds to shut down. Then you should charge the rehabilitation belt.

Charging instruction

1. Before the first use, charge the rehabilitation belt for about 5 hours and fully charge it.
2. Plug in the USB-C cable to the charging jack of the device and the other end of the USB plug into the adaptor with 5V/1A output. The power indicator blinking in blue.
3. When all the power indicators stop flashing and keep on, it means the battery is fully charged.
4. For the remote control, put the CR2032 battery by opening the battery compartment and sliding the door located at the back of the remote control. And the remote control can only remote the rehabilitation belt under the standby mode.
5. Charge or discharge batteries to approximately 50% capacity before storage. Recharge the battery to approximately 50% capacity at least every 3-6 months. Store product at a temperature between 5 °C and 25 °C.
6. Do not use while charging.

 **Cleaning and Maintenance****Cleaning**

- Please ensure to turn off the power and unplug before cleaning.
- Use a soft dry cloth or slightly damp cloth to wipe the product surface.

**CAUTION**

- Do not touch the plug with wet hands.
- Do not wash or spill any liquid into the product.
- Do not use corrosive liquid such as benzene or thinner to clean the product, to prevent color fading and damage.
- When the product is damp, stop using it and do not dismantle. Leave it to dry naturally.

Maintenance

- Keep the product clean and away from the places with high humidity or high temperature and environment with volatile gases.
- If the product will not be used for a long period, please use a dust cover (buy separately) to prevent dust, then store in a safe, dry and cool place.

**CAUTION**

- Do not place it near fire or direct sunlight, to prevent product color fading.

Kindly Reminder

If the product has been used for a long time, please check the conditions as below:

- Burnt smell.
- Unstable power supply after connecting with the power cord
- Power cord overheat
- Power cords damage
- Other abnormalities

**CAUTION**

- If the above phenomenon occurs, please stop using the product to prevent malfunctions and accidents.
- Please turn off the power and unplug immediately. Users are prohibited to dismantle and repair by themselves. Please contact us for repair.

? Solutions and FAQs

The following symptoms do not indicate malfunction. Please check the below phenomenon to find out the solution before contact customer services department. If the problem still cannot be solved, please do not dismantle or repair the product. You should stop using it and contact us immediately.

Issues	Possible Causes	Solutions
No Response The device is blinking	- Low battery	- Recharge the product
Can't feel NMES upon selecting NMES function	- Poor skin conductivity - The intensity maybe too low	- Use water-based gel on your skin and make sure the black color NMES conductive pad is in contact with your skin firmly - Increase the intensity level gradually
NMES intensity cannot go up, always drop to 00	- Poor skin contact	- Apply some water-based gel or water on your skin to enhance the sensation of conductivity. - Wrap firmly the belts on the body. Make sure the black color NMES conductive pad is in contact with your skin firmly
All buttons do not work well	- Battery is low	- Recharge the product

- ? Q: When I wipe the electrode area, it becomes black and dirty.
A: If I rub the electrode area forcefully, the color may transfer. However, it will not dirty the skin during the care process, so it is safe.
- ? Q: After use, symptoms such as redness, inflammation, and itching may occur.
A: Individual responses to NMES may vary. If the symptoms do not improve the next day, please consult a healthcare professional.
- ? Q: Feeling discomfort and fatigue after use.
A: Individual responses to NMES may vary. Before the discomfort subsides, please temporarily reduce the frequency of use.

Technical Specifications

Product name	: Phototherapy NMES Rehabilitation Belt
Model no.	: MAX-M01MC0577B
Belts size	: Rehabilitation belt: 140mm x 880mm : Short extension belt: 90mm x 350mm
Charging voltage	: 5V USB 1A
Battery capacity	: 3.7V 4000mAh
Max. Rated power	: 9.8W
Max. output Current	: 2650mA
NMES – Max. Intensity	: 1 mode with 20 intensity
NMES – Frequency Generation	: 1.667 – 16.67Hz
NMES – Max. output voltage	: 13V – 42.1V (under 1KΩ)
Heat	: 3 levels : Red for 1st Level 38°C > Green 2nd Level 40°C > Blue 3rd Level 45°C (+/- 5°C)
Photon	: 72 pcs. red light : 24 pcs. Infra-red light
Charging time	: 5 hours
Work for when full charge	: 90 minutes
Auto Timer	: 10 minutes

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Trademark holder: Maxcare Industrial Limited