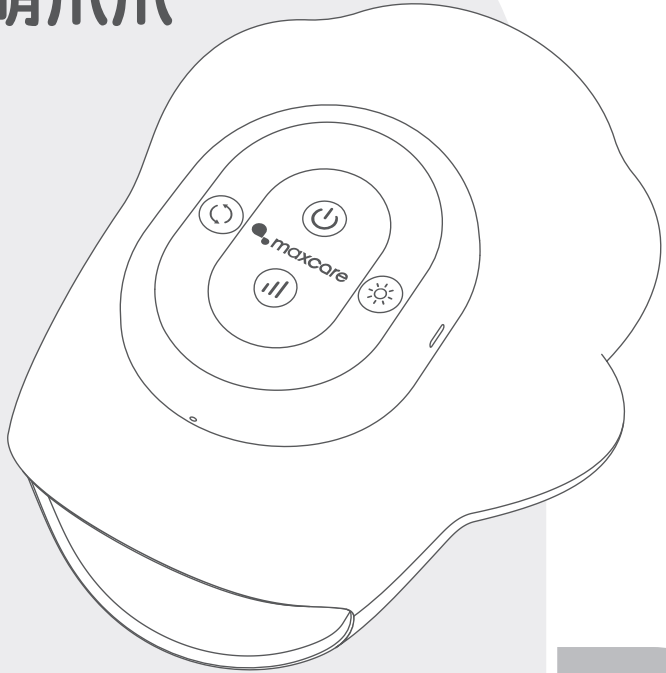
# Adorable Claw

Massager

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MAX-S13H80W

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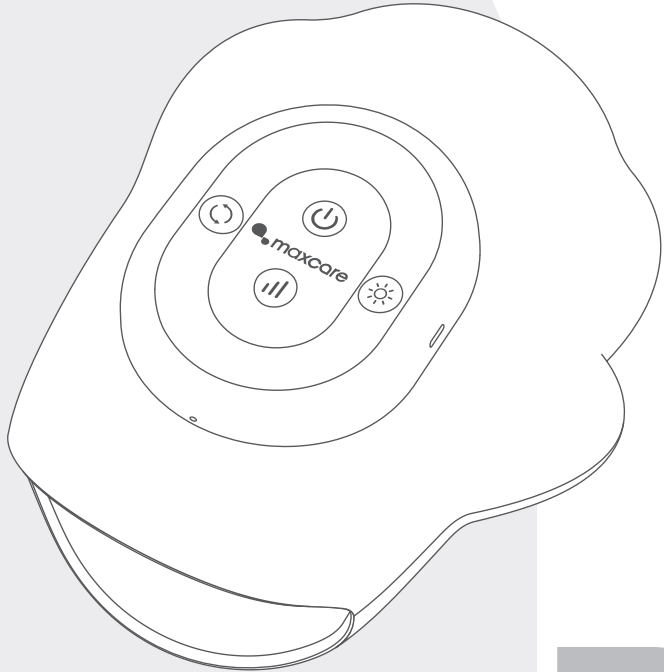
MAX-S13H80W

## COMPANY PROFILE

After two decades of experience in the health industry, **maxcare** founder utilized his knowledge and expertise in massage therapy from an Asian-Pacific perspective to establish our brand - **maxcare** . Our massage therapy business was expanded to Hong Kong in 2012 and is conferred 2017 Hong Kong Top Brand Award. Based in Hong Kong, we are specializing in design & innovation, manufacturing, and retailing of massage and healthcare products.

Believing in the value of 'Technology improves our quality of life' and 'Caring for your family', **maxcare** has been innovating highly qualified products with international patents. Besides, **maxcare** products fulfill various international safety standards, namely EU CE Mark, US UL/ETL Mark, Germany GS Mark, TUV certification, and China CCC Mark.

Motivated by our value, **maxcare** founder has been leading us to be one of the most competitive brands in Hong Kong health industry. **maxcare** applies the contemporary fashion design – neat and modern – on the mechanical engineering of massage chairs. To deliver our value and offer incredible massage experience and healthcare solution to the public, **maxcare** devotes to the innovation, design, manufacturing the top-notch family massage products. By integrating the company philosophy –Spirit of Healthy Life - into our products, **maxcare** brings health to everyone.

# Adorable Claw

## Massager

Thank you for purchasing Maxcare's products. Before using the product, please read this manual thoroughly, particularly the Safety Precautions section, for the correct use of the product.

This manual should be well kept for easy reference. It enables you to perform a safe and efficient operation.

### CONTENT

---

Safety Precautions .....	P.01 - P.03
Name and Function of Components .....	P.04
How to Use .....	P.04 - P.05
Cleaning and Maintenance .....	P.06
Troubleshooting .....	P.07
Technical Specifications .....	P.08

## Safety Precautions

Before use, please read this SAFETY PRECAUTIONS carefully to ensure proper operation.

The following SAFETY PRECAUTIONS help you to use the product correctly. Indicate the risk of danger and damage, the classes of accidents caused by incorrect use are classified into WARNING and CAUTION, hence every user must abide by it.

Indication for WARNING AND CAUTION	
WARNING	Indicating user may be seriously injured once the product is improperly used.
CAUTION	Indicating the user may be injured or the product may be damaged once the product is improperly used.

### 1. Important safety warning

- Please use power supply suited to the product.
- Don't use the product if the power source does not meet the rated voltage to avoid damage or fire hazards.
- Please pull out the plug after use or before cleaning to avoid injury or damage to the product.
- Don't use accessories which are not recommended.
- Don't use the product outdoors.
- Any other usage that is not listed in this manual instruction is forbidden.
- 15 minutes of usage each time is recommended.
- Please don't use it if the case or parts are broken or damaged.
- Please don't use it if the radiator hole is covered or blocked.
- Please don't use it within one hour after having meal.
- Please don't adjust to the greatest intensity at the beginning of the massage to avoid the risk of being injured.
- Keep the plug clean from needles, rubbish or water.
- If the power cord is damaged, don't attempt to change on your own. Please contact us for repair for safety reason.

### 2. Environment for usage

#### CAUTION

- Please don't use it under high temperature and moist environment such as bathrooms.
- Please don't use it immediately while the environmental temperature changes sharply.
- Please don't use it under heavy-dusty or caustic environment.
- Please don't use it where there is not enough space or not good ventilation.
- Please don't use it when surrounding temperature is up to 40°C.
- Please don't use it when surrounding temperature is below 5°C.
- Turn off and pull out the plug during thunderstorms.
- Place the product on a flat surface.
- Don't place it near fire or directly under the sun.
- When moved from a low-temperature storage to a warm environment, the moisture inside the machine may affect, or damage, its functionality. Please rest the product within room temperature for 1 hour before use.

 **Safety Precautions** **3. People who are not suitable to use the product****CAUTION**

- People who are receiving medical treatment, especially physically unfit, are not suitable to use the product.
- People who have cancers, acute illnesses, heart diseases, or severe high blood pressure, etc. are not suitable to use the product.
- Pregnant women or women that are in menses period are not suitable to use the product.
- People who are suffering from fever are not suitable to use the product.
- People who are suffering from heart diseases or carrying electronic medical equipment such as a heart pacemaker are not suitable to use the product.
- People who got injured or skin diseases are not suitable to use the product.
- People who have rheumatoid arthritis, rheumatism, gout are not suitable to use the product.
- People who have osteoporosis are not suitable to use the product.
- People who have blood circulation problems such as thrombosis, severe aneurysm, or acute venous tumor, or infectious skin diseases are not suitable to use the product.
- People who have perception disorders due to diabetes are not suitable to use the product.
- People who have open wounds are not suitable to use the product.
- People whose health condition is poor are not suitable to use the product.
- Patients with tenosynovitis.
- Children under the age of 14 or individuals with mental instability must not use this machine without supervision.
- People who aim for recovery are not suitable to use the product.

 **4. Safety****CAUTION**

- Don't use other therapy devices while using the product.
- The product is not to be used for medical treatment purpose.
- Only one person, i.e. one hand, at a time.
- Recommended usage time of the product 15 minutes/ session.
- Don't pull or drag the power cord when unplugging.
- Don't place heavy objects on the power cord.
- Before use, check if the machine and its parts are damaged. If so, don't use it and contact us for repair.
- Don't tie, bend, pull or twist the power cord to avoid damages.
- Be away from the product during power cut to prevent injury in case power resumes suddenly.
- Stop using the product if there appears any irregularity, and please contact us immediately.
- Stop using the product if it causes discomfort, and please consult a doctor.
- This appliance is not intended for people (including children) with reduced physical, sensory or mental capabilities, or a lack of experience and knowledge, unless they are under supervision or instruction concerning the use of the appliance by a person responsible for their safety.
- Children should be supervised to ensure that they do not play with the appliance.
- Keep the power cord away from hot or moving objects.
- The product is not designed for long-time use. Allow a 30-minute break (power off) for every one hour use. Long-time use could cause overheating and malfunction.

## Safety Precautions

### FORBIDDEN

- It is forbidden to use if there appears dilapidation, crack, leakage or explosion of internal parts on the product.
- It is forbidden to place heavy objects on top of the product.
- It is forbidden to let children use or play with the product.
- It is forbidden to use the product while preparing to sleep or sleeping.
- It is forbidden to use while drunk or sick.
- It is forbidden to put fingers or other foreign objects into the gaps, to prevent harm or malfunction.
- In case of irregular situation (e.g. smoking or smelling), please immediately stop using and cut off power.
- It is forbidden to have water or any liquid on the main parts of the product, to avoid electric shock or fire.
- It is forbidden to operate the console or forcefully drag the power cord with wet hands.
- People whose bodies are wet are forbidden to use the product.
- The surface of the product will be hot, operate carefully.

### DISASSEMBLY FORBIDDEN

Anyone else, except maintenance personnel specified by our company, must not disassemble, maintain, or remodel the product. If any problem, please contact us immediately.

## 5. Maintenance

- The product should be maintained by us, users are prohibited to disassemble or maintain by themselves.
- Please always turn off the main power after use.
- Don't use the product if the socket is loose.
- If the product is to be left idle for long, roll up the wires and keep it in dry and dustless environment.
- Don't keep the product at a high temperature or near fire. Don't keep under direct sunlight.
- Please clean the product with dry cloths. Don't use thinner, benzene or alcohol.
- The mechanical components of the product are specially designed and made, it is not necessary to be specially maintained.
- Don't point sharp objects at the product.
- Please use it intermittently; don't keep it continuously running for long.

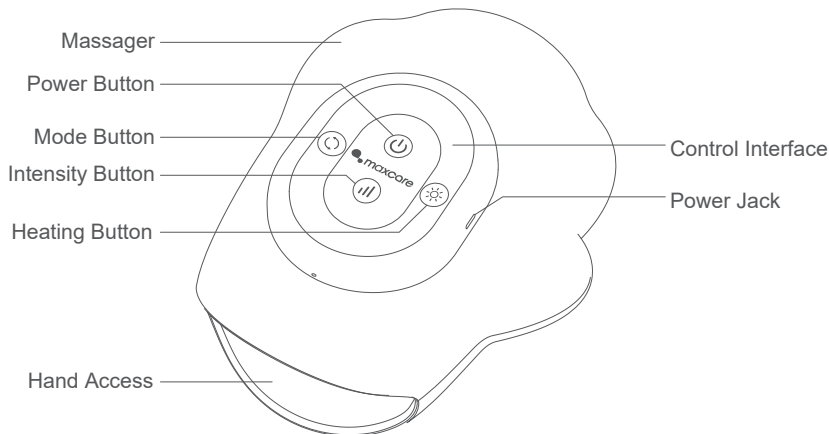
## 6. Solution of ordinary malfunctions

- It is natural if the motors make noises during operation.
- If the remote control can't work normally, please check the plug and socket if they are connected, and make sure the switch is turned on.
- If the rated working time is over, the product will automatically activate the power-off button. If the product works continuously for long, the temperature protector will make the product automatically power-off, it should be used again after half an hour's rest.

## Name and Function of Components

### Product structure

## How to Use

### Instruction

#### Power Button

Press and hold the power button for 1.5 seconds to turn on the device and it will automatically enter the "Vitality Mode". To turn off the device, press and hold the power button again.

#### Mode Button

Press the Mode button to switch between different modes. The mode switching sequence is as follows: "Vitality Mode" - "Soothing Mode" - "Power Mode". The product is set to the default "Vitality Mode" upon startup.

**Mode 1: "Vitality Mode"** - Medium intensity, indicated by a blue light. This mode focuses on massaging the fingers and palm area to activate oxygen in the blood, promote blood circulation, relax muscle tissues, and relieve fatigue. Additionally, it incorporates techniques such as palm kneading, back-of-the-hand pressure, and wrist compression to further enhance the vitality of the hands, relieving stiffness and soreness from daily work and activities.

**Mode 2: "Soothing Mode"** - Medium intensity, indicated by a green light. This mode focuses on massaging the palm and wrist area to effectively promote blood circulation, relieve fatigue, uplift the spirit, enhance alertness, and alleviate tension. It employs a progressive wave-like dynamic kneading technique and gradual meridian massage to soothe hand soreness and fatigue.

**Mode 3: "Power Mode"** - High intensity, indicated by a red light. This mode focuses on massaging the fingers, palm, and wrist areas, providing deep stimulation to the hand meridians. It efficiently relieves hand fatigue and soothes conditions such as tenosynovitis and wrist joint pain.

## How to Use

### Heating Button

Press the Heating button to activate the heating function. Press the Heating button again to deactivate the heating.

### Intensity Button

Press the intensity button to select the desired intensity level. The sequence for switching intensity levels is as follows: "Low Intensity" - "Medium Intensity" - "High Intensity".

Low Intensity: Indicated by blue light.

Medium Intensity: Indicated by green light.

High Intensity: Indicated by red light.

### Charging instruction

1. Connect the USB charging cable to the product's charging port.
2. When the battery level is below 10%, the power indicator light will flash red to indicate low battery warning. After 1 minute, the product will automatically shut down and can only be powered on again after charging.
3. The default working time of the product is 15 minutes. To use it again, you need to turn it on again.
4. The product cannot be turned on and used during the charging process.
5. During charging, the indicator light will show a breathing light pattern and will remain solid green when fully charged.
6. It supports multiple charging methods, including a computer USB interface, power bank, and 5V = 1A power adapter.

### Indicator light description

When the machine is in use, the indicator light displays different colors to indicate the current battery level of the product.

Indicator Light	Status	Description
Red Light	Flash	0 < Current battery level ≤ 9%
Orange Light	Always on	10% ≤ Current battery level ≤ 50%
Green Light	Always on	51% ≤ Current battery level ≤ 100%

When the product is in charging mode, the indicator light will display different colors to indicate the current battery level, as follows:

Indicator Light	Status	Description
Red Light	Breathing light	0 < Current battery level ≤ 9%
Orange Light	Breathing light	10% ≤ Current battery level ≤ 50%
Green Light	Breathing light	51% ≤ Current battery level ≤ 100%
Green Light	Always on	100%

 **Cleaning and Maintenance****Cleaning**

- Before cleaning, please make sure the switch has been turned off, and the power plug has been pulled out from socket.
- Use soft cloths to clean the surface of the product and the control interface.

**CAUTION**

- Do not pull the plug from socket with wet hands.
- Rinsing the product with water is not allowed.
- Do not use corrosive cleaning agents like benzene or thinner to clean the product, otherwise, the color of product will fade, and the internal parts may be corrupted or broken.
- Do not disassemble the product when it is wet.

**Maintenance**

- When not in use, turn off the main power, and remove the plug from socket.
- Keep the product clean and away from high-temperature/ humid areas, or any environment with volatile gases.
- If the product is to be left idle for long, properly cover it with a dust cover (self-purchase).

**CAUTION**

- Please keep the product away from direct sunlight, or high temperature, to avoid fading.

**Kindly Reminder**

If the product has been used for long, please check the conditions as below:

- Burnt smell.
- Unstable power supply after contact with the power cord.
- Power cord heating up.
- Power cords damage.
- Other abnormalities.

**CAUTION**

- If above mentioned happens, please stop using the product.
- Switch off the power supply, pull out the plug, and contact us.
- Don't disassemble the product by self.

## Troubleshooting

The following symptoms do not indicate malfunction. Please check the below phenomenon to find out the solution before contact customer services team. If the problem still cannot be solved, please do not dismantle or repair the product. You should stop using it and contact us immediately.

Issue	Reason	Solution
Not working	The power button is not pressed and held.	Press and hold the power button for 1.5 seconds.
	Low battery level.	Charge the product.
	The 15-minute massage program has been completed, and the product has automatically entered standby mode.	Press the power button to restart.
Unable to charge	The charger is not plugged in properly.	Plug in the charger properly.
	Not using the designated charging cable.	Use the charging cable specified for the product.

## Technical Specifications

Product Name	: Adorable Claw Messenger
Product Model	: MAX-S13H80W
Operating Voltage	: 3.7V
Rated Power	: 6W
Charging Current	: 1A
Charge Time	: 5Hrs
Battery Capacity	: 2200mAh
Battery Life	: Vitality mode available for $\geq 100$ minutes Soothing mode available for $\geq 100$ minutes Power mode available for $\geq 90$ minutes
Heating temperature	: $42 \pm 3^{\circ}\text{C}$
Product Size	: 197x244x111mm
Product Net Weight	: 710g
Interface Type	: Type-C

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Trademark holder: Maxcare Industrial Limited