
INSPIREFOOT

Foot Massager

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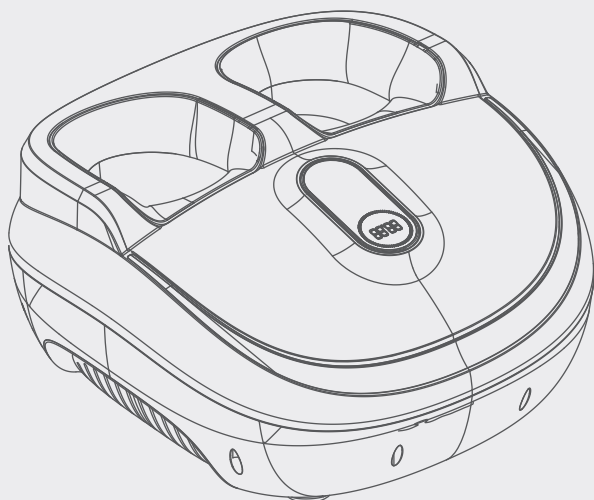MAX-Y03508P

MAX-Y03508B

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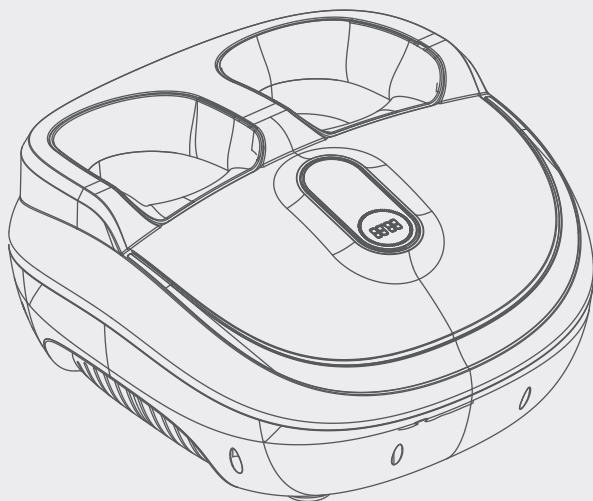
MAX-Y03508P
MAX-Y03508B

COMPANY PROFILE

After two decades of experience in the health industry, **maxcare** founder utilized his knowledge and expertise in massage therapy from an Asian-Pacific perspective to establish our brand - **maxcare** . Our massage therapy business was expanded to Hong Kong in 2012 and is conferred 2017 Hong Kong Top Brand Award. Based in Hong Kong, we are specializing in design & innovation, manufacturing, and retailing of massage and healthcare products.

Believing in the value of 'Technology improves our quality of life' and 'Caring for your family', **maxcare** has been innovating highly qualified products with international patents. Besides, **maxcare** products fulfill various international safety standards, namely EU CE Mark, US UL/ETL Mark, Germany GS Mark, TUV certification, and China CCC Mark.

Motivated by our value, **maxcare** founder has been leading us to be one of the most competitive brands in Hong Kong health industry. **maxcare** applies the contemporary fashion design – neat and modern – on the mechanical engineering of massage chairs. To deliver our value and offer incredible massage experience and healthcare solution to the public, **maxcare** devotes to the innovation, design, manufacturing the top-notch family massage products. By integrating the company philosophy –Spirit of Healthy Life - into our products, **maxcare** brings health to everyone.

INSPIREFOOT

Foot Massager

Thank you for purchasing Maxcare's products. Before using the product, please read this manual thoroughly, particularly the Safety Precautions section, for the correct use of the product.

This manual should be well kept for easy reference. It enables you to perform a safe and efficient operation.

CONTENT

Safety Precautions	P.01 - P.02
Name and Function of Components	P.03
How to Use	P.04 - P.06
Cleaning and Maintenance	P.07
Troubleshooting	P.08
Technical Specifications	P.09

 **Safety Precautions** **1. Important safety warning**

- Please read and comply with the user manual before operation.
- Keep product away from sources of water to prevent short circuit caused by ingress of water into the product.
- Please use an adapter which is in accordance with the specification of the product.
- Disconnect the product after use or before cleaning to prevent injury or damage.
- Do not use outdoors.
- Recommended usage time of this product is 15 minutes/ session.
- Please do not keep this product running continuously for over 30 minutes.
- Please check if any part of this product is damaged before use. If so, please stop using it immediately and do not dismantle or repair by yourself. You should contact us for repair.
- Please do not use it if drunk or feeling unwell.
- Please do not use it within one hour after having meal.
- Do not wash the product

 **2. Environment for usage**

- Do not use it in a humid environment such as a bathroom.
- Do not use it when there is a drastic change in temperature.
- Do not use it under heavy-dusty or caustic environment.
- Do not use it in an environment with poor ventilation.
- Please turn off the power and unplug the power plug after use.
- Do not use it when the surrounding temperature is up to 40°C.
- Please do not use it when surrounding temperature is below 5°C. The lubricant will get thicker quickly and increase the friction between parts.
- Turn off and unplug during thunderstorms.
- Do not place it near fire or direct sunlight, to prevent product color fading.
- When moved from a low-temperature storage to a warm environment, the moisture inside the machine may affect, or damage, its functionality. Please rest the product within room temperature for 1 hour before use.

 **3. People who are not suitable to use the product**

- People who have rashes or infectious skin diseases.
- People who have diabetes, sensory impairment, malignant tumors, acute illnesses, heart diseases, or severe high blood pressure, etc.
- People who suffer from any bleeding tendency, thrombosis, severe aneurysm or acute venous aneurysm and other symptoms of blood circulation disorder.
- People who are suffering from swollen feet or inflamed feet.
- People who are suffering from phlebitis, thrombosis, joint dysfunction, inflammation, weak bones or osteoporosis, etc.
- People who suffer from heart disease or carrying electronic medical device such as a heart pacemaker.
- People who have pins or screws inside the body, or artificial joints.
- People who have a fever.
- Pregnant or menstruating women.
- Children under 14-year-old or who are delirious without being supervised are not allowed to use this product.
- For those who need to rest or feel unwell after being instructed by a doctor.
- The product has a heating function. It must be used carefully by those who are sensitive to heat.

Safety Precautions





4. Safety

- Do not place heavy objects on top of the product.
- Please place the product on flat ground.
- Do not operate under a blanket or pillow. Overheating can occur causing fire, electric shock, or other injuries.
- Do not stand while using the product.
- The product can only be used by a single user at a time.
- Check if the voltage in accordance with the specification of this product.
- Do not use wet hands to pull the plug or touch the control panel.
- Do not allow liquid to get into the product, to prevent electric shock or damage to this product.
- Do not pull or drag the power cord when unplugging. Please be gentle.
- Do not damage the wires or change the circuit of this product.
- Do not use wet cloths to clean the electric parts such as the switch and plug. Please clean this product with dry cloths. Do not use thinner, benzene or alcohol.
- Be away from this product and turn off the plug during a power cut, to prevent injury in case power resumes suddenly.
- Stop using this product if there is any irregularity, and please contact us immediately.
- Stop using this product if you feel any physical abnormality. Please consult a doctor.
- This product is not intended for people (including children) who are physically handicapped, sensory or neurological deficits, or lack of experience and knowledge, unless they are under supervision or instruction concerning the use of the product by their guardian.
- Children should be supervised to ensure they do not play it as a toy.
- Do not directly pull or drag the power cord when moving the product to avoid damage.
- This product should be maintained by us, users are prohibited to dismantle, repair or maintain by themselves.
- Keep away from fire, high temperature or flammable, to avoid malfunctions and dangers that may cause deformation, discoloration or damage.
- Do not use other therapy devices while using this product.
- This product is a personal care massager and cannot replace medicines and medical devices.
- Please strictly follow the instructions, otherwise an accident or physical discomfort may occur.

⚙️ Name and Function of Components

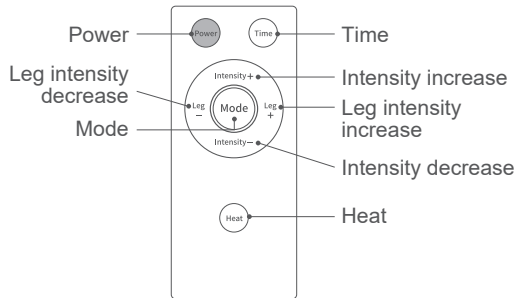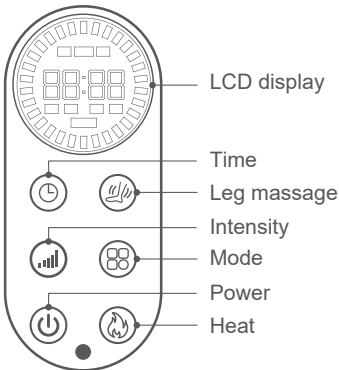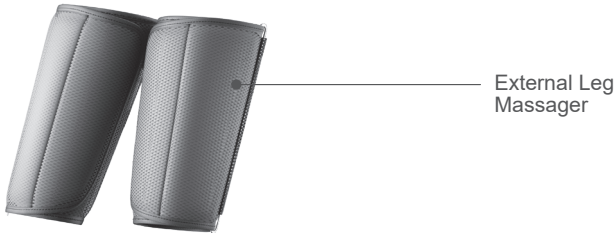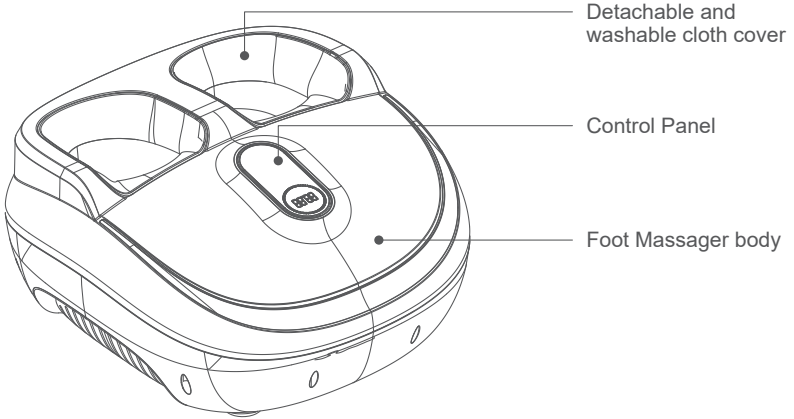Product structure

 How to Use

How to Use

When connected to the power supply, the product will enter standby mode, indicated by the light displaying "OFF."

Buttons	Description	Display									
ON/OFF  	Press the power button once to turn on the device and enter the default working mode: P1 mode, with a working time of 15 minutes, air pressure intensity at L5, and heat level at H1.	Lights On									
	Press the power button again to stop the device from operating and enter standby mode.	OFF									
Time  	The default working time is 15 minutes. Each press increases the time by 5 minutes, up to a maximum of 30 minutes.	Countdown									
Intensity 	There are 9 levels of air pressure intensity. When turned on, the default setting is L5. Press the button once to adjust the intensity to L6, then to L7, L8, L9, L1, L2, L3, L4, and finally back to L5 in a sequential loop. Weakest <table border="1" style="display: inline-table; vertical-align: middle;"> <tr> <td>L1</td><td>L2</td><td>L3</td><td>L4</td><td>L5</td><td>L6</td><td>L7</td><td>L8</td><td>L9</td> </tr> </table> Strogest	L1	L2	L3	L4	L5	L6	L7	L8	L9	Intensity L1-L9
	L1	L2	L3	L4	L5	L6	L7	L8	L9		
	Intensity increase Intensity+	Press this button to increase the air pressure intensity.									
Intensity decrease Intensity-	Press this button to decrease the air pressure intensity.										
Mode  Mode	There are 5 preset modes to choose from, with P1 mode being the default upon startup. P1: Airbag kneading works + roller rotates and massages the feet for 60 seconds with a slow frequency. It intermittently holds the feet, and the roller switches between forward and reverse cycles every 60 seconds.	Mode P1									
	P2: Airbag kneading works + roller rotates and massages the feet for 30 seconds with a moderate frequency. It intermittently holds the feet, and the roller switches between forward and reverse cycles every 30 seconds.	Mode P2									
	P3: Airbag kneading works + roller rotates and massages the feet for 10 seconds with a moderate frequency. It intermittently holds the feet, and the roller switches between forward and reverse cycles every 10 seconds.	Mode P3									
	P4: Airbag kneading is turned off + roller rotates and massages the feet for 60 seconds with a slow frequency. The roller switches between forward and reverse cycles every 60 seconds.	Mode P4									
	P5: Airbag kneading works + roller stops working.	Mode P5									

 How to Use

How to Use


Buttons	Description	Display
Leg massage 	Default setting for the 3-level calf air pressure mode, with air pressure turned off at startup (L0).	L0
	Press once to adjust to L1 (weak), then to L2 (medium) > L3 (strong) > L0 (Turn off) in a sequential loop.	Leg L1-L3
Leg intensity increase Leg +	Press this button to increase the intensity of calf air pressure.	
Leg intensity decrease Leg -	Press this button to decrease the intensity of calf air pressure.	
Heat  	Default setting for the 2-level heating temperature selection, with H1 mode being the default at startup. H1 mode: Heating temperature 45°C	Heat H1
	Press once to enter H2 mode: Heating temperature 50°C	Heat H2
	Press again to turn off the heating function.	Heat H0

 How to Use**How to Use**


1. Please check if the product is damaged before use. Place the product on a stable surface and insert the plug into a power outlet.
2. The display screen will show "OFF" indicating that the product is in standby mode.
3. Adjust the foot to a comfortable position and press the power button to turn on. The default massage mode is P1, with 15 minutes massage time, air pressure intensity at L5, and heat level at H1.
4. During use, press  button to select the massage mode. Press  button to adjust the air pressure intensity. Press  button to turn on the leg air pressure massage mode. Press  button to turn on or off the heating mode. Press  button to adjust the massage time.
5. When the massage time is up, the product will automatically shut down. Alternatively, you can press  button directly to turn it off during use. The display screen will show "OFF" to indicate that it has entered standby mode.

Note: Before using the external leg massager, please insert the connector into the socket on the back of the machine as shown in the diagram below.

Before Installation:

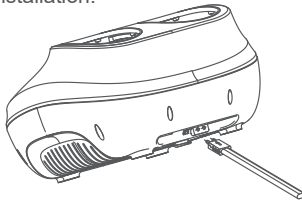
After Installation:

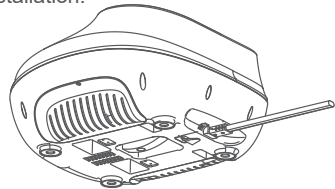 **Cleaning and Maintenance****Cleaning**

- Please ensure to turn off the power and unplug before cleaning.
- Use a soft dry cloth or slightly damp cloth to wipe the product surface.
- The cloth cover is removable and washable. Please refer to the instructions in the diagram to unzip and detach it.

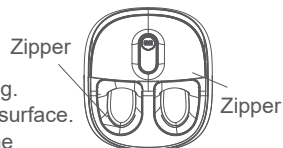- Do not touch the plug with wet hands.
- Do not wash or spill any liquid on the product.
- Do not use corrosive liquid such as benzene or thinner to clean the product, to prevent color fading and damage.
- When the product is damp, stop using it and do not dismantle. Leave it to dry naturally.

Maintenance

- Keep the product clean and away from the places with high humidity or high temperature and environment with volatile gases.
- If the product will not be used for a long period, please use a dust cover (buy separately) to prevent dust, then store in a safe, dry and cool place.

- Please ensure to turn off the power and unplug when it's not in use.

Kindly Reminder

If the product has been used for a long time, please check the conditions as below:

- Burnt smell.
- Unstable power supply after connecting with the power cord.
- Power cord overheat.
- Power cords damage.
- Do not pull the cloth cover and the zipper forcefully to avoid damage.
- Other abnormalities.

- If the above phenomenon occurs, please stop using the product to prevent malfunctions and accidents.
- Please turn off the power and unplug immediately. Users are prohibited to dismantle and repair by themselves. Please contact us for repair.

Troubleshooting

The following symptoms do not indicate malfunction. Please check the below phenomenon to find out the solution before contact customer services team. If the problem still cannot be solved, please do not dismantle or repair the product. You should stop using it and contact us immediately.

No.	Issues	Possible Causes	Solutions
1	“Tik tik” or “Tak tak” sound during operation	Normal mechanical sound during operation or the friction sound between the massage balls and cloth cover	No action required
2	Not working	The adapter has not been well connected	Connect the adapter properly
		The power switch is off	Switch on the power button
		The massager is used for too long, overheat protection will be activated	Allow the massager to cool down before using again
3	Abrupt stop during operation	Power disconnected	Check if the power cord is connected
		Rated working time is up and the massager is automatically stop	Press the power button to restart
		If you have been using the massager for a long time, overheat protection may get activated.	Turn the massager off and allow it to cool down before using again

Technical Specifications

Product Name	: InspireFoot
Product Model	: MAX-Y03508B (Silver) MAX-Y03508P (Purple)
Rated Voltage	: 220V~ 50Hz
Rated Power	: 60W
Rated Time	: 15 minutes
Product Size	: 47.5 x 27 x 43cm
Heating temperature	: H1 45°C H2 50°C
Packaging Size	: 47.5 x 27 x 43cm
Net Weight	: 5.9kg
Product Weight	: 7.2kg
Components	: Foot massager, external leg massager, remote control (batteries excluded), adapter, user manual

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Trademark holder: Maxcare Industrial Limited