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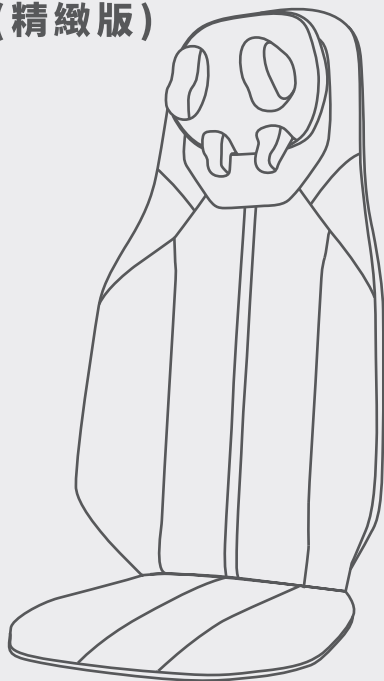
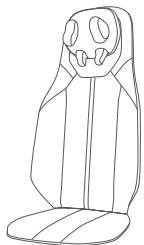
**Compact**  
**GRANDMASTER**  
BACK MASSAGER

## **User Manual**

MAX-S03FRC432

## COMPANY PROFILE

After two decades of experience in the health industry, **maxcare** founder utilized his knowledge and expertise in massage therapy from an Asian-Pacific perspective to establish our brand - **maxcare** . Our massage therapy business was expanded to Hong Kong in 2012 and is conferred 2017 Hong Kong Top Brand Award. Based in Hong Kong, we are specializing in design & innovation, manufacturing, and retailing of massage and healthcare products.

Believing in the value of 'Technology improves our quality of life' and 'Caring for your family', **maxcare** has been innovating highly qualified products with international patents. Besides, **maxcare** products fulfill various international safety standards, namely EU CE Mark, US UL/ETL Mark, Germany GS Mark, TUV certification, and China CCC Mark.

Motivated by our value, **maxcare** founder has been leading us to be one of the most competitive brands in Hong Kong health industry. **maxcare** applies the contemporary fashion design – neat and modern – on the mechanical engineering of massage chairs. To deliver our value and offer incredible massage experience and healthcare solution to the public, **maxcare** devotes to the innovation, design, manufacturing the top-notch family massage products. By integrating the company philosophy –Spirit of Healthy Life - into our products, **maxcare** brings health to everyone.

# **GRANDMASTER** BACK MASSAGER

Compact

Thank you for purchasing Maxcare's products. Before using the product, please read this manual thoroughly, particularly the Safety Precautions section, for the correct use of the product.

This manual should be well kept for easy reference. It enables you to perform a safe and efficient operation.

## CONTENT

---

Safety Precautions .....	P.01 - P.03
Name and Function of Components .....	P.04 - P.05
How to Use .....	P.06
Cleaning and Maintenance .....	P.07
Troubleshooting .....	P.08
Technical Specifications .....	P.09

## Safety Precautions

### 1. Important safety warning

- Please read and comply with the user manual before operation.
- Keep product away from sources of water to prevent short circuit caused by ingress of water into the product.
- Please use an adapter which is in accordance with the specification of the product.
- Disconnect the power supply before cleaning or after use to prevent injury or damage.
- Do not use outdoors.
- Recommended usage time of this product is 15 minutes/ session.
- Please check if any part of this product is damaged before use. If so, please stop using it immediately and do not dismantle or repair by yourself. You should contact us for repair.
- Please do not use it if the cover is broken or damaged.
- Please do not use it if the shell or surface is damaged.
- Please do not drop anything into this product.
- Please do not fall asleep while using this product.
- Please do not use it if drunk or feeling unwell.
- Please do not use it within one hour after having meal.
- Please do not adjust to the maximum intensity at the beginning of the massage to avoid possible injury.
- Keep the plug clean from needles, rubbish or water.
- If the power cord is damaged, do not attempt to change on your own. Please contact us for repair for safety reasons.

### 2. Environment for usage

- Do not use it in a humid environment such as a bathroom.
- Do not use it when there is a drastic change in temperature.
- Do not use it under heavy-dusty or caustic environment.
- Do not use it in an environment with poor ventilation.
- Please switch off the power and unplug the machine after use.
- Do not use it when the surrounding temperature is up to 40°C.
- Please do not use it when surrounding temperature is below 5°C. The lubricant will get thicker quickly and increase the friction between parts.
- Turn off and unplug during thunderstorms.
- Place the product on a flat surface.
- Do not place it near fire or direct sunlight.
- When moved from a low-temperature storage to a warm environment, the moisture inside the machine may affect or damage its functionality. Please rest the product within room temperature for 1 hour before use.

## Safety Precautions

### 3. People who have the below symptoms, please consult a doctor before use

- People who are suffering from osteoporosis.
- People who suffer from heart disease or carrying electronic medical device such as a heart pacemaker.
- People who have a fever.
- Pregnant or menstruating women.
- People who got injured or skin diseases.
- Children under 14-year-old or who are delirious without being supervised are not allowed to use this product.
- For those who need to rest or feel unwell after being instructed by a doctor.
- People whose bodies are wet are forbidden to use this product.
- The product has a heating function. It must be used carefully by those who are sensitive to heat.
- People who are receiving medical treatment, especially those who are physically unfit.
- People who have open wounds.
- People who aim for recovery.

### 4. Safety

- Please do not place heavy objects on this product.
- The product can only be used by a single user at a time.
- Check if the voltage in accordance with the specification of this product.
- Do not use wet hands to pull the plug.
- Do not allow liquid to get into the product, to prevent electric shock or damage to this product.
- Do not pull or drag the power cord when unplugging. Please be gentle.
- Do not damage the wires or change the circuit of this product.
- Do not use wet cloths to clean the electric parts such as the switch and plug.
- Be away from this product during power cut to prevent injury in case power resumes suddenly.
- Stop using this product if there is any irregularity, and please contact us immediately.
- Stop using this product if you feel any physical abnormality. Please consult a doctor.
- This product is not intended for people (including children) who are physically handicapped, sensory or neurological deficits, or lack of experience and knowledge, unless they are under supervision or instruction concerning the use of the product by their guardian.
- Children should be supervised to ensure they do not play it as a toy.
- This product should be maintained by us, users are prohibited to dismantle, repair or maintain by themselves.
- Keep the power cord is damaged, it must be replaced by us to avoid a hazard.
- Please cut off power and pull-out plug before moving.
- If the product is found damaged, broken, suffering from power leakage, or with exposed parts, please stop using it.

## Safety Precautions

### 5. Maintenance

- This product should be maintained by us, users are prohibited to dismantle or repair by themselves.
- Please always turn off the main power after use.
- Do not use this product if the socket is loose.
- If this product is to be left idle for long, please tie up the power cords and keep it in dry and dustless environment.
- Do not keep this product at a high temperature or near fire. Do not keep under direct sunlight.
- Please clean this product with dry cloths. Do not use thinner, benzene or alcohol.
- Do not point sharp objects at this product.
- Please use it intermittently. Do not keep it continuously running for long.


### 6. Solution of ordinary malfunctions


- It is normal if the motor makes noises during operation.
- If the remote control can't work normally, please check the plug and socket if they are connected, and make sure the switch is turned on.
- If the rated working times is over, the power will automatically turn off. If thus product works continuously for too long, the temperature control will automatically turn off the machine, if so, please allow a 30-minute break before resuming.

 Name and Function of Components

Product structure

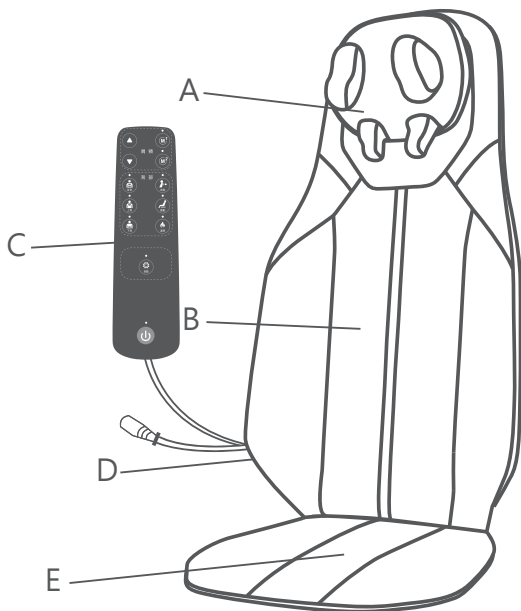- A: Neck massage ball
- B: Back massage ball
- C: Controller
- D: Controller bag
- E: The seat cushion
- F: AC Adaptor plug
- G: DC Adaptor cable
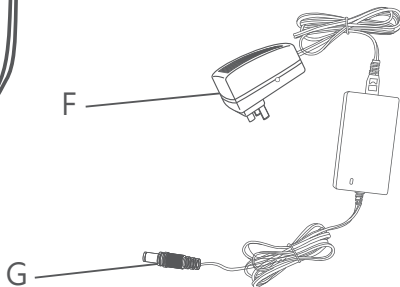
## Name and Function of Components

### How To Use

#### Product operation procedure

1. Place the product upright on a sofa or chair. (use the elastic belt at the back of the massage seat to fix the product on the chair)
2. Connect the power cord with the adapter, insert the round head of the DC cable into the DC connector of the massager, and then insert the AC plug into the socket. Sit comfortably on the massage cushion with your back and neck resting on the product.
3. Start the hand controller to enjoy whole-body massage.

#### Operation instructions


- Neck**
1. Neck adjustment: Press and hold this button to adjust the height of neck massage ball to meet the requirement of different people.
  - 2-3. Neck massage mode: Press this button to turn on/off the neck massage Mode 1: Continuous Kneading or Mode 2: Intermittent Kneading.

- Back**
4. Back massage area: full back, upper back and lower back massage.
  5. Target: kneading massage nodes on a fixed point.
  6. Heat: turn on/off the heating function, the indicator light lights up or goes out accordingly.
  7. Auto mode: turn on the full massage functions, the default mode is vibration low level, heating on, neck massage Mode 1 Continuous Kneading, full back massage, indicator light stays on.  
Press it the second time switches to automatic mode 2: Focus on upper back + lower back massage, indicated by a fast flashing light.  
Press it the third time switches to automatic mode 3: Focus on lower back+ upper back massage, indicated by a slow flashing light.  
Press it the fourth time, switched off the automatic mode and it will cycle through sequentially.

- The seat cushion**
8. Vibration: turn on/off vibration function of the seat cushion, three levels of the intensity can be adjusted, including low level (indicated by slow flashing light), medium level (indicated by fast flashing light), and high level (indicated by steady light).

- ON/OFF**
9. Power button: turn on the product, the power indicator light will light up and the product will enter standby mode.  
Automatic shutdown: Automatically shuts down after 15 minutes of use.  
Manual shutdown: When the product is on, press the power button to turn it off, and the indicator light will go out.

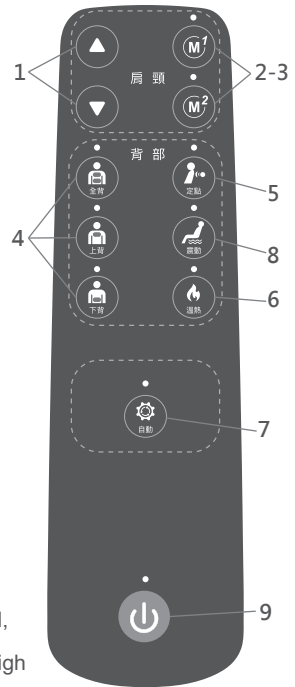※ After the product is shut down, the massage ball will automatically return to its original position and then shut down.

 How to Use

How To Use

Button	Description
	Full Back: Massage full back
	Upper Back: Massage upper back
	Lower Back: Massage lower back
	Target: Fix at a point to specially massage
	Neck position adjustment: Press and hold to adjust the height of the neck massage ball
	Neck massage mode: turn on/off Mode 1: Continuous Kneading / Mode 2: Intermittent Kneading
	Vibration: turn on/off vibration function of the seat cushion, three levels of the intensity can be adjusted, including low level (indicated by slow flashing light), medium level (indicated by fast flashing light), and high level (indicated by steady light).
	Heat: turn on/off the heating function
	Auto mode: turn on the default massage mode, which is vibration low level, heating on, neck massage Mode 1 Continuous Kneading, full back massage, indicator light stays on. Press it the second time switches to automatic mode 2: Focus on upper back + lower back massage, indicated by a fast flashing light. Press it the third time switches to automatic mode 3: Focus on lower back+ upper back massage, indicated by a slow flashing light. Press it the fourth time, switched off the automatic mode and it will cycle through sequentially.
	Power button: turn on/off the product

## Cleaning and Maintenance

### Cleaning

- Please ensure to turn off the power and unplug before cleaning.
- Use a soft dry cloth or slightly damp cloth to wipe the product surface.

### CAUTION

- Do not touch the plug with wet hands.
- Do not wash or spill any liquid on the product.
- Do not use corrosive liquid such as benzene or thinner to clean the product, to prevent color fading and damage.
- When the product is damp, stop using it and do not dismantle. Leave it to dry naturally.

### Maintenance

- Make sure the plug has been discounted from the socket before maintenance.
- Keep the product clean and away from the places with high humidity or high temperature and environment with volatile gases.
- If the product will not be used for a long period, please use a dust cover (buy separately) to prevent dust, then store in a safe, dry and cool place.

### CAUTION

- Do not place it near fire or direct sunlight, to prevent product color fading.

### Kindly Reminder

If the product has been used for a long time, please check the conditions as below:

- Burnt smell.
- The power cords do not work every time.
- The power cords overheat.
- The power cords damaged.
- Other abnormalities.

### CAUTION

- If the above phenomenon occurs, please stop using the product to prevent malfunctions and accidents.
- Please turn off the power and unplug immediately. Users are prohibited to dismantle and repair by themselves. Please contact us for repair.

## ? Troubleshooting


The following symptoms do not indicate malfunction. Please check the below phenomenon to find out the solution before contact customer services team. If the problem still cannot be solved, please do not dismantle or repair the product. You should stop using it and contact us immediately.

No.	Issues	Possible Causes	Solutions
1	“Tik tik” or “Tak tak” sound during operation	Normal mechanical sound during operation or the friction sound between the massage head and cloth cover	No action required
2	Not working	The adapter has not been well connected	Connect the adapter properly
		The power switch is off	Press the ON/OFF button to start
3	Abrupt stop during operation	Power disconnected	Check if the power cord is connected
		Rated working time (15 minutes) is up and the massager is automatically stop	Press the power button to restart
		If you have been using the massager for a long time, overheat protection may get activated	Turn the massager off and allow it to cool down before using again
		The massage balls maybe stuck	Reduce the pressure on the product

## Technical Specifications

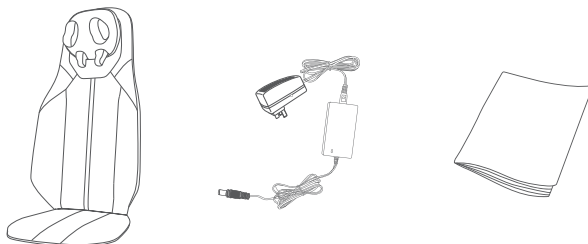Product Name	: Grandmaster Back Massager (Compact)
Product Model	: MAX-S03FRC432G (Green) MAX-S03FRC432B (Beige)
Rated Voltage	: 12V
Rated Power	: 48W
Rated Time	: 15 minutes
Heating temperature	: 45±5°C
Product Size	: 52.5x44.5x78cm
Packaging Size	: 81.5x24.5x47cm
Product Weight	: 6kg
Product Standard	: GB4706.1-2005 GB4706.10-2008 GB4343.1-2018
Components	: Massager, adapter, user manual

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Trademark holder: Maxcare Industrial Limited