POCKET EMS

Massager

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MAX-M01MC0608

POCKET EMS

Massager

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POCKET EMS Massager

Thank you for purchasing our products. Before use, please read this instruction manual carefully to ensure correct operation. Please keep this instruction manual for future reference.

PRODUCT INTRODUCTION

The Pocket EMS Massager is a mini pocket massager that provides EMS (Electronic Muscle Stimulator) to effectively help relax, firming and tonic muscle groups. Comes with 2 electronic devices and 2 gel pads to use on different parts of your body. Cordless and compact size can be used anywhere and anytime.


PRECAUTIONS FOR USE




Never use damaged, deteriorated, or dirty gel pads, otherwise it may make you feel pain/ burned.

- Always make sure gel pads well contact with the electrodes when using. If the electrode touches the skin directly, it may be painful or get burned.
- Attach the product to the area it will be used and then turn the power on.
- When you are done using it, turn the power off before removing it from your body.
- If you sense something abnormal with your skin or body, stop use immediately and consult your doctor. Not doing so may result in an accident or poor health.
- Using the product at a very high level or for a long time may cause muscle fatigue, tingling and/or redness on your skin.
- Do not heat or cool the product before use, as doing so may deform it. Also do not use in a manner that subjects the product to sudden changes in temperature, such as pouring hot water on it, using it in a sauna, etc.
- Clean the area where it will be used.
- The product is not waterproof. Do not use in a location that is high humid or where water may get on it.

PRODUCT STRUCTURE

- A. 2 pcs Gel pad
- B. 2 pcs Electronic Devices
- C. Charging base with cover
- D. USB/C charging cable

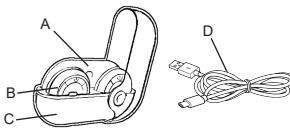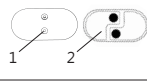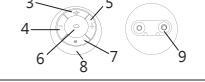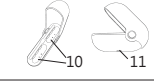Gel Pad	Device	Charging Base
		
1. Snap buttons 2. Electrode area	3. Power button 4. Intensity decrease 5. Intensity increase 6. Power indicator light 7. Mode button 8. Charging contact point 9. Snap button	10. Charging contact points used to connect to the device during charging 11. USB TYPE-C charging socket

CHARGING INSTRUCTION

- Before first use, charge your massager for about 1 hour till it's fully charged.
- Put the 2 devices into the charging base and make sure the devices are sit well on the charging base. The charging contact points on the device are well connected to the charging contact points of the charging base.
- Connect USB/C charging cable to the charging base and the other end of the USB should plug into the 5V 1A USB adapter.
- The power indicator light of device flashes in red when charging and stays on red when fully charged.
- Disconnect the product from the USB-C cable once it is fully charged.
- If the device does not respond to the power button, the battery may be too low. Please charge device.

HOW TO USE

- Connect the gel pads to the devices with the snap buttons.
- Peel the PVC protective film off the gel pads and keep the protective film for gel pad protection after use.
- Then place them on your body where you want to have treatment.
- Press and hold the power button for 2 seconds to begin massage, and the green light will be on. The massage function will activate at the lowest level massage intensity.
- 10 intensity levels available: press "+" or "-" button to increase or decrease the intensity.
- There are 5 massage modes available. Press the "M" button to switch the massage mode. (Note that this is a cycle key).

Here below is the color of the mode that the indicating light will light up.

Mode	Mode 1	Mode 2	Mode 3	Mode 4	Mode 5
Color of LED	Green	Pink	Yellow	Blue	White
Frequency	62.5Hz	62.5Hz	10-16.7 Hz	3.3-33.3 Hz	1-33 Hz

- After use, press and hold the power button for 3 seconds to turn off and the working indicator light will be off.
- The auto timer is 15 minutes, and the red-light flashes for 5 seconds to shut down.
- Please be reminded to stick the protective film onto the patches after use.
- If the intensity or massage mode has been changed during use, the auto timer will count from 15 mins again.

Note

- You can only adjust the mode and intensity level when the device is on and the patch is in contact with your skin. If the patch is not placed correctly on the skin, the massager will automatically turn off.
- Mild discomfort may occur if your skin is too oily or moist, or if the patch is too damaged to function.
- The electrode gel pad can be used about 30 times. Do not place it on or remove the protective film if you do not intend to use it.
- To make it easier for users to identify the buttons during use, there is one dot on the intensity increase button and two dots on the power button.

TROUBLE SHOOTING

The following symptoms do not indicate malfunction. Please check the below phenomenon to find out the solution before contact customer services department. If the problem still cannot be solved, please do not dismantle or repair the product. You should stop using it and contact us immediately.

Problem	Cause	Solution
No Led Light Up	- Low Battery	- Charge it.
LED light up but why can't feel	- No skin contact. - Unit not snapped on gel pad correctly. - Intensity level at 0 - Microprocessors hang up	- Make sure the skin is well touch with the electrode pad. - Check that device is connected firmly with both snap button of the electrode pad. - Make sure the gel pad and device connected well and stick firm on your skin - Reset device by turn it off for 1 minute.
All Button Does Not Work Well	- Microprocessor hang up - Battery is low	- Reset the device by turn it off for 1 minute - Charge the battery
Gel pad is not sticky	- Adhesive surface (hydrogel surface) is dirty or exhausted	- Gel pad is a consuming part, need to change a new gel pad with your distributor.
Skin redness and pricking	- Skin allergy	- There are individual differences in response to the use of EMS. If symptoms do not improve the next day, please consult your doctor.
the LED on and off automatically	- Low battery	- Charge the battery.

MAINTENCE

- Do not immerse in liquid or allow liquid to come into contact the device.
- Never use abrasive or chemical cleaners. Use a soft and dry cloth to clean the product.
- Charge or discharge batteries to approximately 50% capacity before storage. Recharge the battery to approximately 50% capacity at least every 3-6 months.
- Store product in dry, well-ventilated, non-corrosive gas environment. At a temperature between 5°C and 25°C.

TECHNICAL SPECIFICATION

Product Name	: Pocket EMS Massager	EMS Pulse Frequency	: 1Hz - 62.5Hz
Model Number	: MAX-M01MC0608B (Black)	Rated Power	: 105MW
	: MAX-M01MC0608W (Wood)	Massage mode	: 5
Input voltage	: USB 5V 1A	Intensity	: 10
Max. current	: 35mA	Battery capacity	: 3.7V 60mAh
		Charging & Duration	: Charge 1 hour can use for 6 hours
		Auto Timer	: 15 minutes