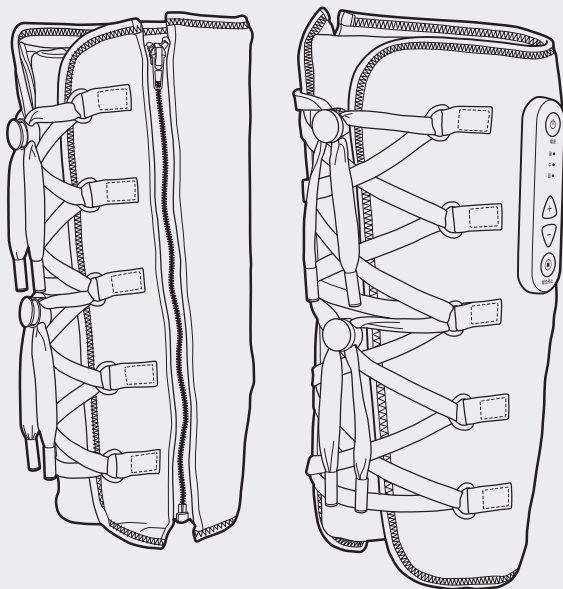
LEG SHAPING FAIRY

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MAX-M01MC0606A

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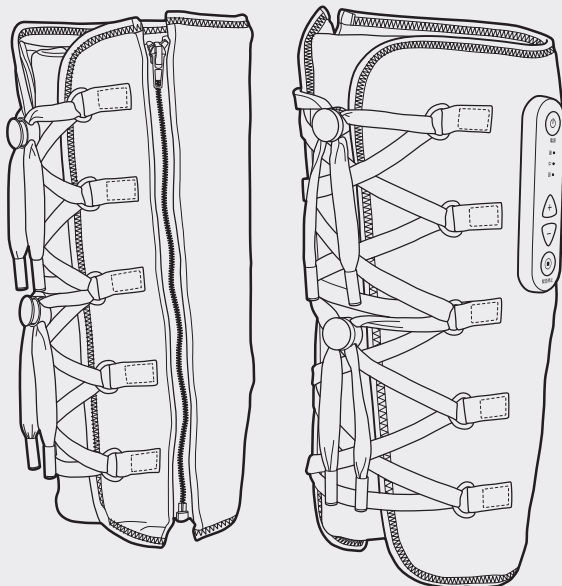
MAX-M01MC0606A

COMPANY PROFILE

After two decades of experience in the health industry, **maxcare** founder utilized his knowledge and expertise in massage therapy from an Asian-Pacific perspective to establish our brand - **maxcare** . Our massage therapy business was expanded to Hong Kong in 2012 and is conferred 2017 Hong Kong Top Brand Award. Based in Hong Kong, we are specializing in design & innovation, manufacturing, and retailing of massage and healthcare products.

Believing in the value of 'Technology improves our quality of life' and 'Caring for your family', **maxcare** has been innovating highly qualified products with international patents. Besides, **maxcare** products fulfill various international safety standards, namely EU CE Mark, US UL/ETL Mark, Germany GS Mark, TUV certification, and China CCC Mark.

Motivated by our value, **maxcare** founder has been leading us to be one of the most competitive brands in Hong Kong health industry. **maxcare** applies the contemporary fashion design – neat and modern – on the mechanical engineering of massage chairs. To deliver our value and offer incredible massage experience and healthcare solution to the public, **maxcare** devotes to the innovation, design, manufacturing the top-notch family massage products. By integrating the company philosophy –Spirit of Healthy Life - into our products, **maxcare** brings health to everyone.

LEG SHAPING FAIRY

Thank you for purchasing Maxcare's products. Before using the product, please read this manual thoroughly, particularly the Safety Precautions section, for the correct use of the product.

This manual should be well kept for easy reference. It enables you to perform a safe and efficient operation.

CONTENT

Safety Precautions	P.01 - P.02
Name and Function of Components	P.03
How to Use	P.04
Cleaning and Maintenance	P.05
Troubleshooting	P.06
Technical Specifications	P.07

Safety Precautions

1. Important safety warning

- Please read and comply with the user manual before operation.
- Keep product away from sources of water to prevent short circuit caused by ingress of water into the product.
- Please use an adapter which is in accordance with the specification of the product.
- Disconnect the product after use or before cleaning to prevent injury or damage.
- Do not use outdoors.
- Recommended usage time of this product is 15 minutes/ session.
- Please do not keep this product running continuously for over 30 minutes.
- Please check if any part of this product is damaged before use. If so, please stop using it immediately and do not dismantle or repair by yourself. You should contact us for repair.
- Please do not use it if drunk or feeling unwell.
- Please do not use it within one hour after having meal.
- Do not wash the product
- Please do not keep this product running continuously for over 30 minutes.

2. Environment for usage

- Do not use it in a humid environment such as a bathroom.
- Do not use it when there is a drastic change in temperature.
- Do not use it under heavy-dusty or caustic environment.
- Do not use it in an environment with poor ventilation.
- Please turn off the power and unplug the power plug after use.
- Do not use it when the surrounding temperature is up to 40°C.
- Turn off and unplug during thunderstorms.
- Do not place it near fire or direct sunlight, to prevent product color fading.

3. People who are not suitable to use the product

- People who have rashes or infectious skin diseases.
- People who have diabetes, sensory impairment, malignant tumors, acute illnesses, heart diseases, or severe high blood pressure, etc.
- People who suffer from any bleeding tendency, thrombosis, severe aneurysm or acute venous aneurysm and other symptoms of blood circulation disorder.
- People who are suffering from swollen feet or inflamed feet.
- People who are suffering from phlebitis, thrombosis, joint dysfunction, inflammation, weak bones or osteoporosis, etc.
- People who suffer from heart disease or carrying electronic medical device such as a heart pacemaker.
- People who have pins or screws inside the body, or artificial joints.
- People who have a fever.
- Pregnant or women who just gave birth.
- Children under 14-year-old or who are delirious without being supervised are not allowed to use this product.
- For those who need to rest or feel unwell after being instructed by a doctor.
- This product is not intended for medical use and is only intended for calf massage purposes.

Safety Precautions

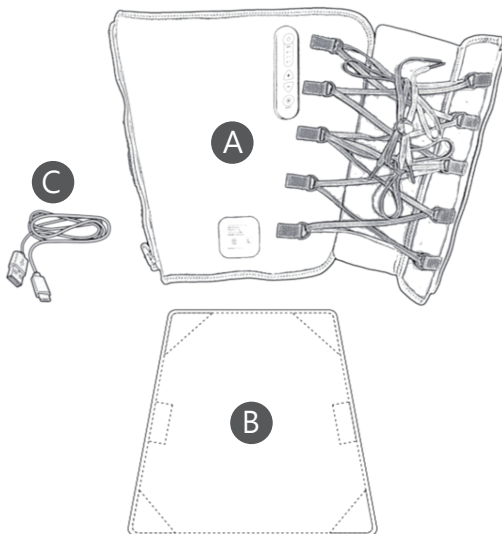
4. Safety

- Do not place heavy objects on top of the product.
- Do not operate under a blanket or pillow. Overheating can occur causing fire, electric shock, or other injuries.
- The product can only be used by a single user at a time.
- Check if the voltage in accordance with the specification of this product.
- Do not use wet hands to pull the plug or touch the control panel.
- Do not allow liquid to get into the product, to prevent electric shock or damage to this product.
- Do not pull or drag the power cord when unplugging. Please be gentle.
- Do not damage the wires or change the circuit of this product.
- Do not use wet cloths to clean the electric parts such as the switch and plug. Please clean this product with dry cloths. Do not use thinner, benzene or alcohol.
- Be away from this product and turn off the plug during a power cut, to prevent injury in case power resumes suddenly.
- Stop using this product if there is any irregularity, and please contact us immediately.
- Stop using this product if you feel any physical abnormality. Please consult a doctor.
- This product is not intended for people (including children) who are physically handicapped, sensory or neurological deficits, or lack of experience and knowledge, unless they are under supervision or instruction concerning the use of the product by their guardian.
- Children should be supervised to ensure they do not play it as a toy.
- Do not directly pull or drag the power cord when moving the product to avoid damage.
- This product should be maintained by us, users are prohibited to dismantle, repair or maintain by themselves.
- Keep away from fire, high temperature or flammable, to avoid malfunctions and dangers that may cause deformation, discoloration or damage.
- Do not use other therapy devices while using this product.
- This product is a personal care massager and cannot replace medicines and medical devices.
- Please strictly follow the instructions, otherwise an accident or physical discomfort may occur.

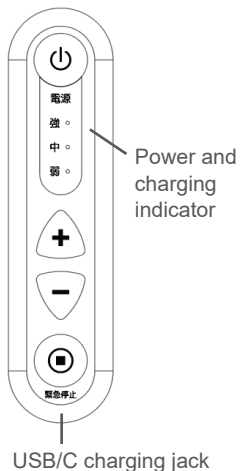
Name and Function of Components

Product structure

- A. 2 x Air Leg Compressor
- B. 2 x Detachable Pad
- C. 1 x USB /C charging cable

Control panel and functions

Power button: Press and hold for two seconds to switch on/off the device.
Short press to stop the product immediately.

Increase level button,
3 levels available – Low > Medium > High

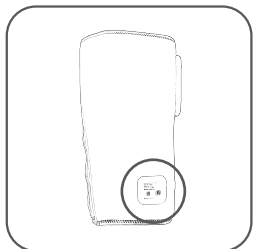
Decrease level button,
3 levels available – High > Medium > Low

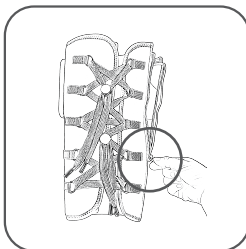
Emergency button:
Short press to stop the product immediately

 **How to Use****Charging instructions**

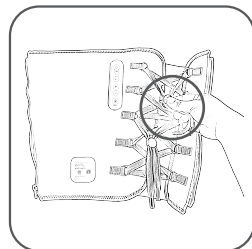
1. Connect the USB/C cable to the charging jack of both control panel.
2. Then connect the other end USB head to charging adaptor and plug into the mains supply.
3. The charging indicator will be blinking in blue while charging and stay on blue when fully charged.
4. The product takes 5-6 hours to fully charge and can be used approximately 9 times. Each use lasts for 10 minutes.
5. After the charging adapter is fully charged, be sure to disconnect and unplug the power

How to use

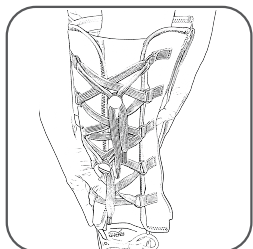
1. Choose the right sleeve first

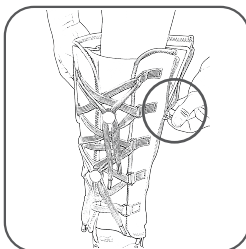
2. Unzip the sleeve

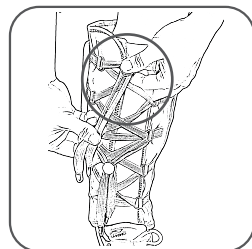
3. Loose the string

4. Wrap the lower leg

5. Zip the sleeve

6. Adjust the string to get comfortable massage

Cleaning and Maintenance

Cleaning

- Please ensure to turn off the power and unplug before cleaning.
- Use a soft dry cloth or slightly damp cloth to wipe the product surface.
- Only Detachable pads can be machine washed.

CAUTION

- Do not touch the plug with wet hands.
- Do not wash or spill any liquid on the product.
- Do not use corrosive liquid such as benzene or thinner to clean the product, to prevent color fading and damage.
- When the product is damp, stop using it and do not dismantle. Leave it to dry naturally.
- Do not iron/Do not use bleach/Do not tumble dry/Do not dry clean.
- Detachable pads can be machine washed or hand washed below 30°C.

Maintenance

- Keep the product clean and away from the places with high humidity or high temperature and environment with volatile gases.
- If the product will not be used for a long period, please use the storage Bag to prevent dust, then store in a safe, dry and cool place.

CAUTION

- Please ensure to turn off the power and unplug when it's not in use.

Kindly Reminder

If the product has been used for a long time, please check the conditions as below:

- Burnt smell.
- Unstable power supply after connecting with the power cord.
- Power cord overheat
- Power cords damage
- Other abnormalities.

CAUTION

- If the above phenomenon occurs, please stop using the product to prevent malfunctions and accidents.
- Please turn off the power and unplug immediately. Users are prohibited to dismantle and repair by themselves. Please contact us for repair.

Troubleshooting

The following symptoms do not indicate malfunction. Please check the below phenomenon to find out the solution before contact customer services department. If the problem still cannot be solved, please do not dismantle or repair the product. You should stop using it and contact us immediately.

Issues	Possible Causes	Solutions
Not working	Low battery	Recharge the product
	Haven't press the power button	Press and hold the ON/OFF button to start
	Rated working time (10 minutes) is up and the massager is automatically stop	Press the power button to restart
Abrupt stop during operation	Low battery	Recharge the product
	If you have been using the massager for a long time, overheat protection may get activated.	Turn the massager off and allow it to cool down for 30-50 mins before using it again.

Technical Specifications

Product Name	: Leg Shaping Fairy
Product Model	: MAX-M01MC0606A
Rated Voltage	: 5V 2A
Rated Power	: 2.35W (Max.)
Battery Capacity	: 1800 mAh
Rated Time	: 10 minutes
Charging time and battery life	: Charging time of 3.5-4.5 hours; Battery life of 90 minutes
Packaging Size	: 36.5 x 23 x 9.2 cm
Product Weight	: 1.34 kg

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Trademark holder: Maxcare Industrial Limited