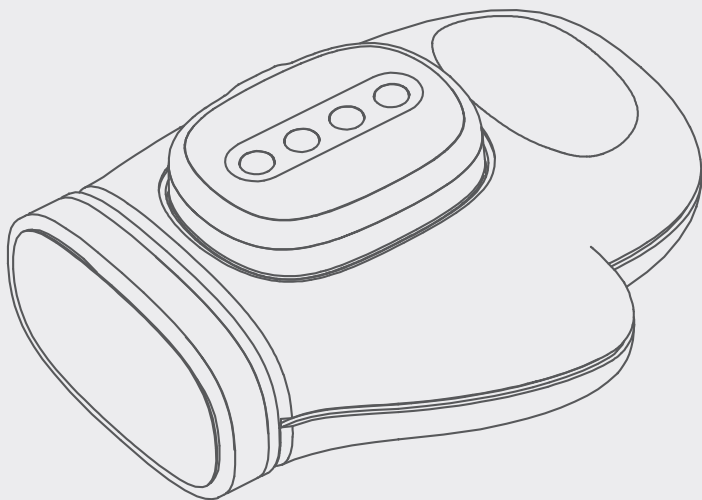
Cat Style

Hand massager

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MAX-S15ST1802A

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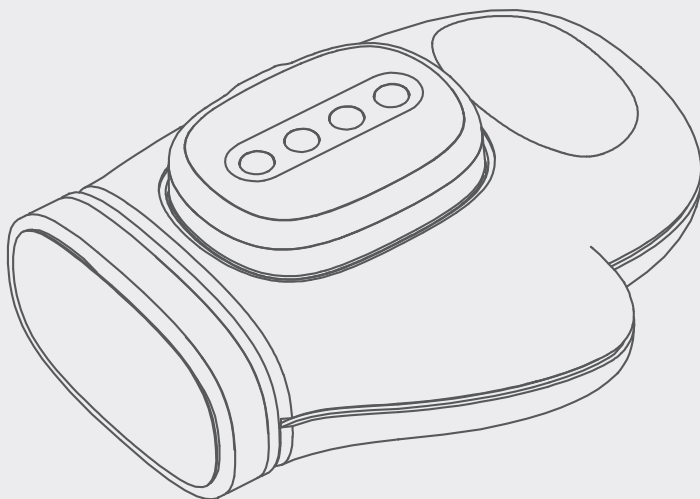
MAX-S15ST1802A

COMPANY PROFILE

After two decades of experience in the health industry, **maxcare** founder utilized his knowledge and expertise in massage therapy from an Asian-Pacific perspective to establish our brand - **maxcare** . Our massage therapy business was expanded to Hong Kong in 2012 and is conferred 2017 Hong Kong Top Brand Award. Based in Hong Kong, we are specializing in design & innovation, manufacturing, and retailing of massage and healthcare products.

Believing in the value of 'Technology improves our quality of life' and 'Caring for your family', **maxcare** has been innovating highly qualified products with international patents. Besides, **maxcare** products fulfill various international safety standards, namely EU CE Mark, US UL/ETL Mark, Germany GS Mark, TUV certification, and China CCC Mark.

Motivated by our value, **maxcare** founder has been leading us to be one of the most competitive brands in Hong Kong health industry. **maxcare** applies the contemporary fashion design – neat and modern – on the mechanical engineering of massage chairs. To deliver our value and offer incredible massage experience and healthcare solution to the public, **maxcare** devotes to the innovation, design, manufacturing the top-notch family massage products. By integrating the company philosophy –Spirit of Healthy Life - into our products, **maxcare** brings health to everyone.

Cat Style

Hand massager

Thank you for purchasing Maxcare's products. Before using the product, please read this manual thoroughly, particularly the Safety Precautions section, for the correct use of the product.

This manual should be well kept for easy reference. It enables you to perform a safe and efficient operation.

CONTENT

| | |
|---------------------------------------|-------------|
| Safety Precautions | P.01 - P.02 |
| Name and Function of Components | P.03 |
| Before use | P.04 |
| How to Use | P.05 - P.06 |
| Cleaning and Maintenance | P.07 |
| Troubleshooting | P.08 |
| Technical Specifications | P.09 |

Safety Precautions

Before use, please read this SAFETY PRECAUTIONS carefully to ensure proper operation. The following SAFETY PRECAUTIONS help you to use the product correctly.

1. Important safety warning

- Please use power supply suited to the product.
- Don't use the product if the power source does not meet the rated voltage to avoid damage or fire hazards.
- Please pull out the plug after use or before cleaning to avoid injury or damage to the product.
- Don't use accessories which are not recommended.
- Any other usage that is not listed in this manual instruction is forbidden.
- Please don't use it if the case or parts are broken or damaged.
- When in use, do not wear hard objects such as jewelry or nail decorations around the massage area.

2. Environment for usage

- Please don't use it under high temperature and moist environment such as bathrooms.
- Please don't use it under heavy-dusty or caustic environment.
- Please don't use it where there is not enough space or not good ventilation.
- Please don't use it when surrounding temperature is up to 40°C.
- Please don't use it when surrounding temperature is below 5°C.
- When moved from a low-temperature storage to a warm environment, the moisture inside the machine may affect, or damage, its functionality. Please rest the product within room temperature for 1 hour before use.

3. People who are not suitable to use the product

- People who are receiving medical treatment, especially physically unfit, are not suitable to use the product.
- People who have cancers, acute illnesses, heart diseases, or severe high blood pressure, etc. are not suitable to use the product.
- People who have rheumatoid arthritis, rheumatism, gout are not suitable to use the product.
- People who have osteoporosis are not suitable to use the product.
- People who have perception disorders due to diabetes are not suitable to use the product.
- People who have open wounds are not suitable to use the product.
- People whose health condition is poor are not suitable to use the product.
- Patients with tenosynovitis.
- Children under the age of 14 or individuals with mental instability must not use this machine without supervision.
- People who aim for recovery are not suitable to use the product.

Safety Precautions

4. Safety

- Don't use other therapy devices while using the product.
- Please do not use for purposes other than hand massage.
- The product is not to be used for medical treatment purpose.
- This product is for single-handed use only, please do not use both hands simultaneously.
- Recommended usage time of the product 15 minutes/ session.
- Before use, check if the machine and its parts are damaged. If so, don't use it and contact us for repair.
- Stop using the product if there appears any irregularity, and please contact us immediately.
- Stop using the product if it causes discomfort, and please consult a doctor.
- This appliance is not intended for people (including children) with reduced physical, sensory or mental capabilities, or a lack of experience and knowledge, unless they are under supervision or instruction concerning the use of the appliance by a person responsible for their safety.
- Children should be supervised to ensure that they do not play with the appliance.
- It is forbidden to use the product while preparing to sleep or sleeping.
- Please do not use with wet hands.
- This appliance has a heating surface. Caution is required when used by individuals who are not sensitive to heat.

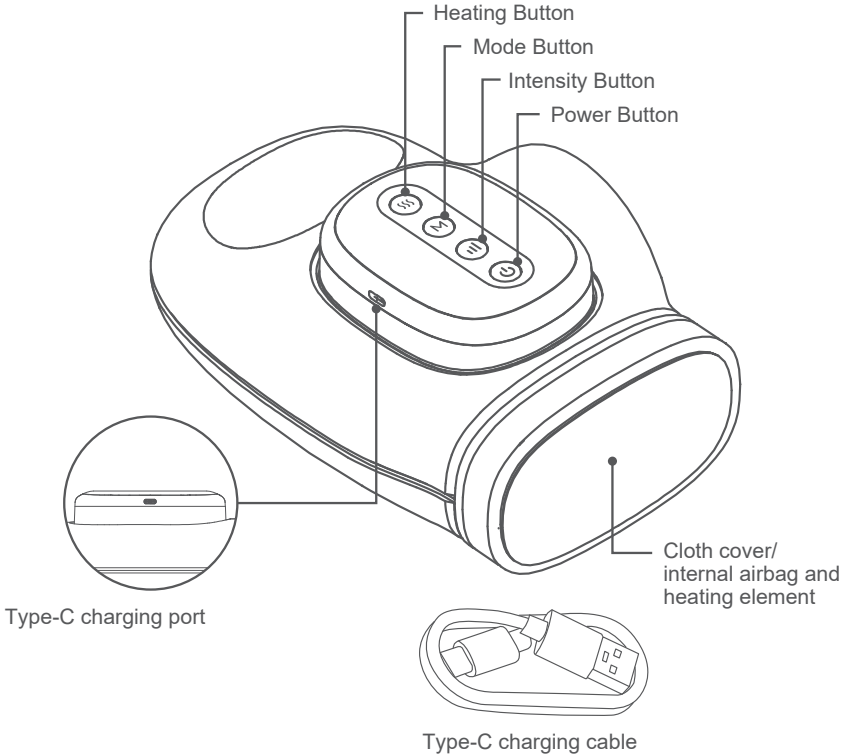
5. Maintenance

- The product should be maintained by us, users are prohibited to disassemble or maintain by themselves.
- Don't keep the product at a high temperature or near fire. Don't keep under direct sunlight.
- Please clean the product with dry cloths. Don't use thinner, benzene or alcohol.
- Please use it intermittently; don't keep it continuously running for long.

 Name and Function of Components

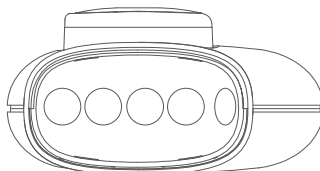
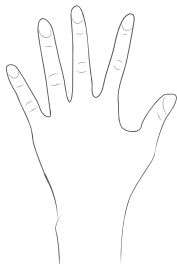
Product structure

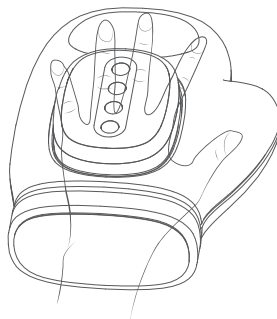
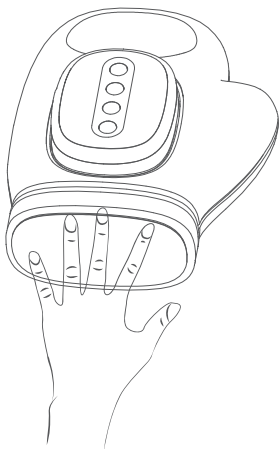
Before use

Please check if the product has power. If the battery is low and the product shuts down automatically or cannot be turned on, please recharge it before using it again.

1. Please ensure that your hands are clean. Do not use with wet hands.

2. Insert your fingers into the correct slots in sequence.

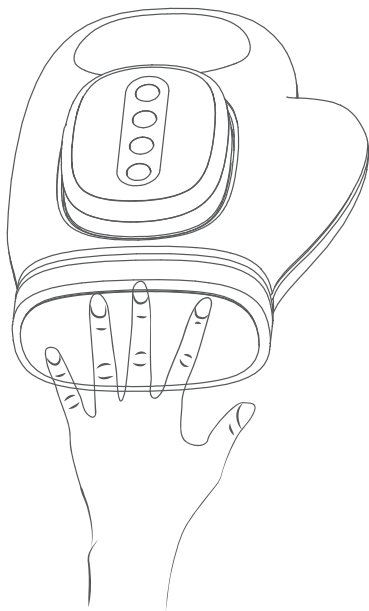
Note that inserting and removing your hand must be done while the machine is turned off.

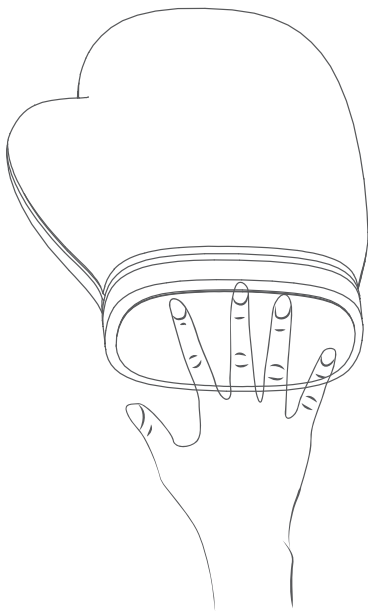
 How to Use

Introduction to using the left and right hands:

This product is designed for single-hand use. Please do not use both hands simultaneously.

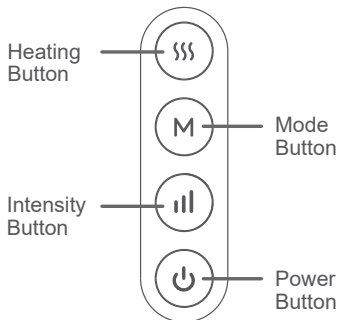
Left hand

Right hand

 How to Use

Button instructions

Power Button

- Press and hold the power button to turn on/off.
- Sufficient battery: green light
- Low battery: flashing red light
- Charging: slow flashing red light
- Fully charged: green light stays on

Heating Button

- Press the heating button to switch between low/medium/high/off heating functions. The default setting when turned on is without heating.
- Low setting: green light
- Medium setting: yellow light
- High setting: red light
- Turn Off: light off

Mode Button

- The default mode is Mode 1.
- Press the mode button to switch between Mode 1/ Mode 2 / Mode 3.
- Mode 1: green light
- Mode 2: blue light
- Mode 3: red light

Intensity Button

- Press the intensity button to switch between levels 1-5, with the default setting at level 1.
- Level 1: green light
- Level 2: blue light
- Level 3: yellow light
- Level 4: purple light
- Level 5: bright red light

| Function | Air Pressure | Time |
|----------|---------------------------------|------------|
| Mode 1 | Finger-Palm-Wrist (Overall) | 15 minutes |
| Mode 2 | Finger-Palm-Wrist (Progressive) | 15 minutes |
| Mode 3 | Finger-Palm-Wrist (Alternating) | 15 minutes |

Charging instructions

1. Connect the USB charging cable to the product's charging port.
2. The default working time for the product is 15 minutes. It needs to be restarted for subsequent use.
3. During the charging process, the product cannot be turned on.
4. For prolonged periods of non-use, it is recommended to recharge the product approximately every 30 days to prevent the battery from over-discharging protection, which may lead to charging abnormalities or inability to charge.

Cleaning and Maintenance

Cleaning

- Before cleaning, please make sure the switch has been turned off.
- Use soft cloths to clean the surface of the product and the control interface.

CAUTION

- Do not pull the plug from socket with wet hands.
- Rinsing the product with water is not allowed.
- Do not use corrosive cleaning agents like benzene or thinner to clean the product, otherwise, the color of product will fade, and the internal parts may be corrupted or broken.
- Do not disassemble the product when it is wet.
- When not in use, turn off the main power.
- Keep the product clean and away from high-temperature/ humid areas, or any environment with volatile gases.
- If the product is to be left idle for long, properly cover it with a dust cover (self-purchase).

Kindly Reminder

If the product has been used for long, please check the conditions as below:

- Burnt smell.
- Unstable power supply after contact with the power cord.
- Power cord heating up.
- Power cords damage.
- Other abnormalities.

CAUTION

- If above mentioned happens, please stop using the product.
- Switch off the power supply, pull out the plug, and contact us.
- Don't disassemble the product by self.

? Troubleshooting

The following symptoms do not indicate malfunction. Please check the below phenomenon to find out the solution before contact customer services team. If the problem still cannot be solved, please do not dismantle or repair the product. You should stop using it and contact us immediately.

| Issue | Reason | Solution |
|------------------|---|---|
| Not working | The power button is not pressed and held. | Press and hold the power button for 1.5 seconds. |
| | Low battery level. | Charge the product. |
| | The 15-minute massage program has been completed, and the product has automatically entered standby mode. | Press the power button to restart. |
| Unable to charge | The charger is not plugged in properly. | Plug in the charger properly. |
| | Not using the designated charging cable. | Use the charging cable specified for the product. |

Technical Specifications

| | |
|----------------------|-----------------------------|
| Product Name: | : Hand massager - Cat style |
| Product Model: | : MAX-S15ST1802A |
| Operating Voltage: | : 3.7V |
| Rated Power: | : 5W |
| Charging voltage: | : 5V |
| Charging Current: | : 1A_MAX |
| Charge Time: | : 2Hrs |
| Battery Capacity: | : 2000mAh |
| Battery Life: | : 1.5-2 Hrs |
| Heating temperature: | : 40-50°C |
| Product Size: | : 19*24.5*10CM |
| Product Net Weight: | : 685g |

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Trademark holder: Maxcare Industrial Limited