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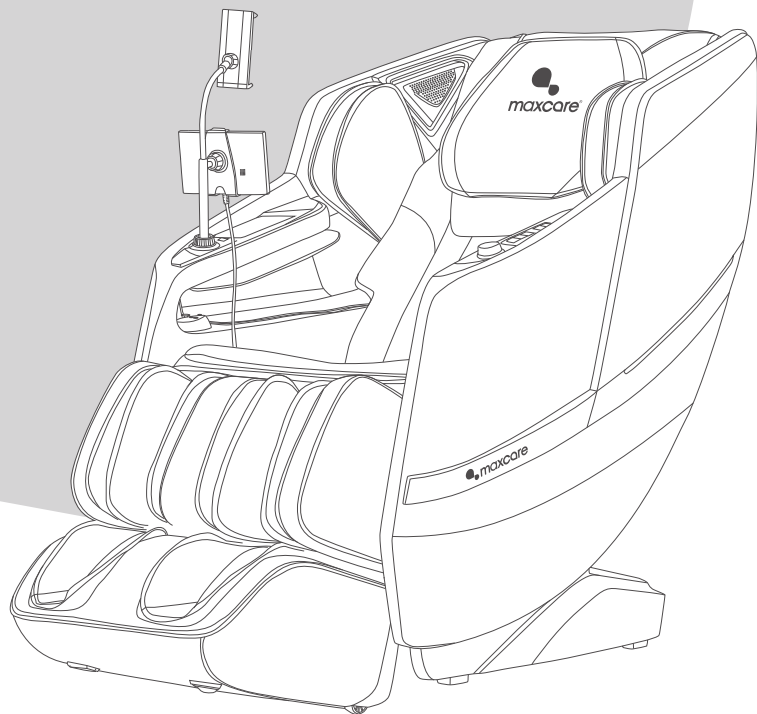
— Lovey Pro —

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MAX-Z05A369

maxcare®

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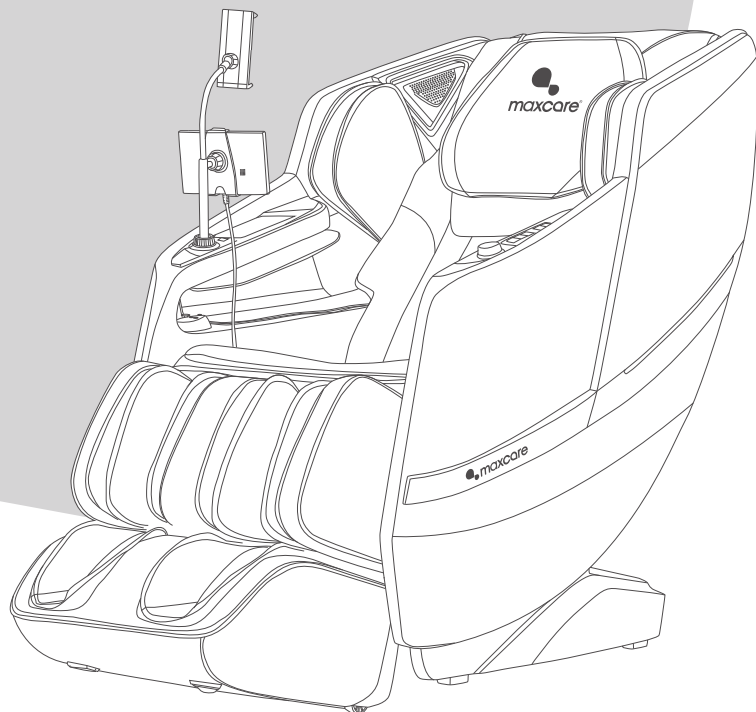
maxcare®

COMPANY PROFILE

After two decades of experience in the health industry, **maxcare** founder utilized his knowledge and expertise in massage therapy from an Asian-Pacific perspective to establish our brand - **maxcare**. Our massage therapy business was expanded to Hong Kong in 2012 and is conferred 2017 Hong Kong Top Brand Award. Based in Hong Kong, we are specializing in design & innovation, manufacturing, and retailing of massage and healthcare products.

Believing in the value of 'Technology improves our quality of life' and 'Caring for your family', **maxcare** has been innovating highly qualified products with international patents. Besides, **maxcare** products fulfill various international safety standards, namely EU CE Mark, US UL/ETL Mark, Germany GS Mark, TUV certification, and China CCC Mark.

Motivated by our value, **maxcare** founder has been leading us to be one of the most competitive brands in Hong Kong health industry. **maxcare** applies the contemporary fashion design – neat and modern – on the mechanical engineering of massage chairs. To deliver our value and offer incredible massage experience and healthcare solution to the public, **maxcare** devotes to the innovation, design, manufacturing the top-notch family massage products. By integrating the company philosophy –Spirit of Healthy Life - into our products, **maxcare** brings health to everyone.

— Lovey Pro —

Thank you for purchasing **maxcare**'s products. Before using the product, please read this manual thoroughly, particularly the Safety Precautions section, for the correct use of the product.

This manual should be well kept for easy reference. It enables you to perform a safe and efficient operation.

CONTENT

Safety Precautions	P.01 - P.06
Name and Function of Components	P.07 - P.13
How to Use	P.14 - P.35
Cleaning and Maintenance	P.36 - P.37
FAQs and Solutions	P.38
Technical Specifications	P.39

Safety Precautions

1. Important safety warning

- Please read this manual instruction carefully before operation.
- Any other usage that is not listed in this manual instruction is forbidden.
- Keep children and pets away from the movable parts of this product to avoid accidents caused by the up-and-down-ward movement of the calf rest and back cover.
- When in use, sliding or transporting, please ensure that there are no people or pets nearby (the back, bottom, front and sides of the product), and that no foreign objects or pets inside the product, to avoid accidents.
- Please use well-grounding power supply suited to this product.
- Please pull out the plug after use or before cleaning to avoid injury or damage to the product.
- Do not use unauthorized accessories.
- Do not use this product outdoors.
- People whose bodies are wet are forbidden to use this product.
- Recommended usage time of this product is 20 minutes/ session.
- Please check if any part of this products is damaged before use. If so, please do not use it and contact us for repair.
- Please do not use it if the cover is broken or damaged.
- Please do not use it if the shell or surface is damaged.
- Please do not drop anything into this product.
- Please do not fall asleep while using this product.
- Please do not use it if drunk or feeling unwell.
- Please do not use it within one hour after having meal.
- Please do not adjust to the maximum intensity at the beginning of the massage to avoid possible injury.
- Please do not keep sharp objects in pockets when use.
- Keep the plug clean from needles, rubbish or water.
- If the power cord is damaged, do not attempt to change on your own. Please contact us for repair for safety reasons.

2. Environment for usage

- Please do not use it in a moist environment such as bathrooms.
- Please do not use it when there is a drastic change in temperature.
- Please do not use it under heavy-dusty or caustic environment.
- Please do not use it where there is not enough space or not good ventilation.
- Please switch off the power and unplug the machine after use.
- Please do not use it when surrounding temperature is up to 40°C.
- Please do not use it when surrounding temperature is below 5°C. The lubricant can thicken quickly and increase the friction between parts.
- Turn off and pull out the plug during thunderstorms.
- Place the product on a flat surface.
- Do not use outdoors.
- Do not place it near fire or directly under the sun.
- When moved from a low-temperature storage to a warm environment, the moisture inside the machine may affect, or damage, its functionality. Please rest the product within room temperature for 1 hour before use.

3. People who have the below symptoms, please consult a doctor before use

- People who are suffering from osteoporosis.
- People who have cancers, acute illnesses, heart diseases, or severe high blood pressure, etc. are not suitable to use this product.
- People who have/ had back pain, or legs, waist, neck numbness/ injuries (disc herniation, spondylosis and cervical protrusion, etc. patients) are not suitable to use this product.
- People who have rheumatoid arthritis, rheumatism, gout are not suitable to use this product.
- People who have respiratory difficulties are not suitable to use this product.
- People who have blood circulation problems such as thrombosis, severe aneurysm, or acute venous tumor, or infectious skin diseases are not suitable to use this product.
- People who have perception disorders due to diabetes are not suitable to use this product.
- People who have spinal abnormalities or spine curvature disorders are not suitable to use this product.
- People whose health condition is poor are not suitable to use this product.
- People who are suffering from heart disease or carrying electronic medical equipment such as heart pacemaker.
- People who have a fever.

Safety Precautions

- Pregnant women or women that are in menses period.
- People who got injured or skin diseases.
- Children under 14-year old or people mentally unbalanced without being supervised are not allowed to use this machine.
- People who are told by doctor to have rest or people who feel unwell.
- The appliance has heating function. People who are insensitive to heat must be careful when using the appliance.
- People who are receiving medical treatment, especially those who are physically unfit.
- People who have open wounds.
- People who aim for recovery.

4. Safety

- Please do not place heavy objects on this product.
- Do not use other therapy devices while using this product.
- This product is not to be used for medical treatment purpose.
- Only one person at a time.
- Recommended usage time of this product is 20 minutes/ session.
- It is possible to place a mat underneath the machine to avoid damage to the floor.
- For safety reason, maintain a minimum of 0.5m between the product and others when use.
- Check if the voltage is suitable to the specification of this product.
- Do not use wet hands to pull the plug.
- Do not allow water to get into the product, to avoid electric shock or damage to this product.
- Do not pull or drag the power cord when unplugging. Please be gentle.
- Do not damage the wires or change the circuit of this product.
- Do not use wet cloths to clean the electric parts such as the switch and plug.
- Be away from this product during power cut to prevent injury in case power resumes suddenly.
- Stop using this product if there appears any irregularity, and please contact us immediately.
- Stop using this product if it causes discomfort, and please consult a doctor.
- This appliance is not intended for people (including children) with reduced physical, sensory or mental capabilities, or a lack of experience and knowledge, unless they are under supervision or instruction concerning the use of the appliance by a person responsible for their safety.
- Children should be supervised to ensure that they do not play with the appliance.
- If the power cord is damaged, it must be replaced by us in order to avoid a hazard.
- Do not ride on the handle to avoid injury or malfunction.
- Keep the power cord away from hot or moving objects.
- Please cut off power and pull out plug before moving.
- If the product is found damaged, broken, suffering from power leakage, or with exposed parts, please stop using it.

5. Maintenance

- This product should be maintained by us, users are prohibited to disassemble or maintain by themselves.
- Please always turn off the main power after use.
- Do not use this product if the socket is loose.
- If this product is to be left idle for long, please tie up the power cords and keep it in dry and dustless environment.
- Do not keep this product at a high temperature or near fire. Do not keep under direct sunlight.
- Please clean this product with dry cloths. Do not use thinner, benzene or alcohol.
- The mechanical components of this product are specially designed and made, no special maintenance required.
- Do not point sharp objects at this product.
- Do not roll or pull this product above uneven ground, it should be lifted before moving.
- Please use it intermittently; do not keep it continuously running for long.

6. Solution of ordinary malfunctions

- It is normal if the motor makes noises during operation.
- If the remote control can't work normally, please check the plug and socket if they are connected, and make sure the switch is turned on.
- If the rated working time is over, the power will automatically go off. If this product works continuously for too long, the temperature control will automatically turn off the machine, if so, please allow a 30-minute break before resuming.

**CAUTION**

Please make sure all the power is cut off before moving.
The power cord and all other cables should be disconnected.

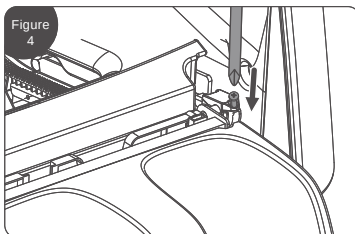
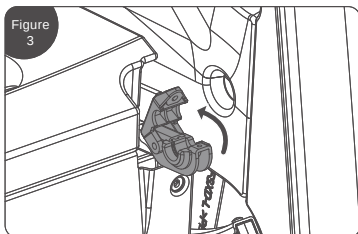
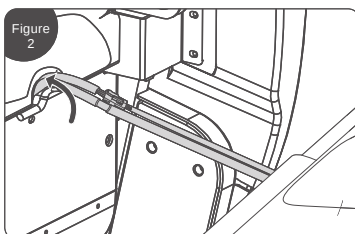
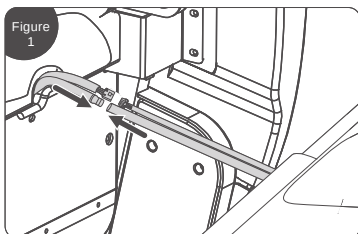
Safety Precautions

Installation parts

Calf installation instructions

1. Insert the wire of the calf and the air pipe into the front socket of the steel frame assembly (Figure 1).
2. Then push the connected wire and air pipe into the hole in the backrest frame (Figure 2).
3. After the wire docking is completed, open the calf fixing covers on both sides of the steel frame assembly (Figure 3).
4. Lift the calf into the notch, close the fixing covers, and secure it in place with 2 pieces screws (Figure 4).

Note: this step can be ignored for all-in-one installation.

Safety Precautions

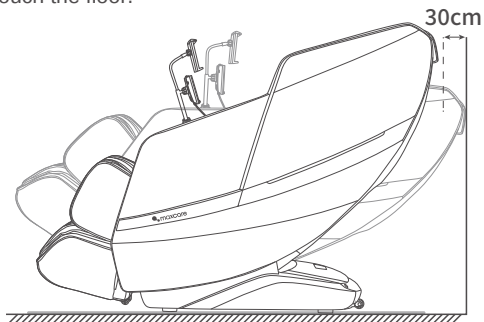
Place massage chair

Installation placement

Make sure there is enough space for the massage chair can be tilted at least 30 cm.

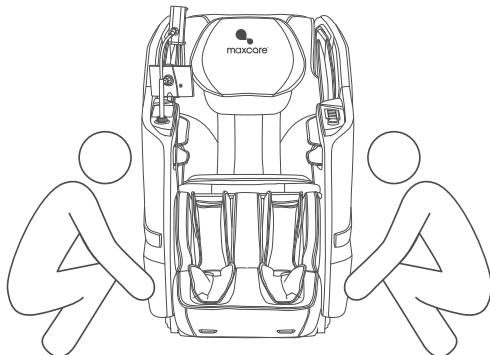
Do not expose the massage chair to direct sunlight or high-temperature environments (e.g., in front of heating devices). Place it in an area where it can avoid heat and direct sunlight. For long-term use in a fixed location, it is recommended to place a soft cushion underneath. The cushion should be large enough to cover the area where the chair and the footrest touch the floor.

Product lift instruction

Two-person lifting required due to product weight.

1. Please lift the massage chair when moving it on delicate surfaces, such as wooden floors.
2. When moving on similar surfaces, lifting is recommended to avoid damaging the casters. Be cautious when lowering the chair into the designated position to avoid hand or foot injuries. Do not release the chair until it is completely flat. To prevent floor damage, it is recommended to place a mat or protective covering on the floor.

Safety Precautions

Place massage chair

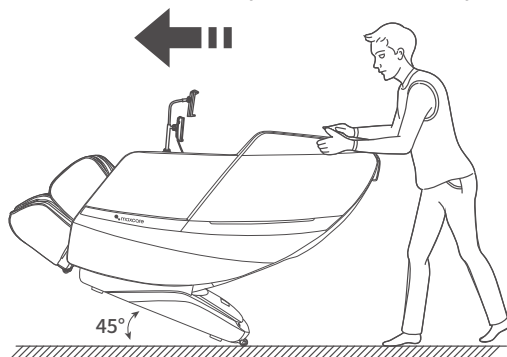
Move method

Use the rear casters to move the massage chair. Tilt the chair backward approximately 45° as shown in the diagram to move it.

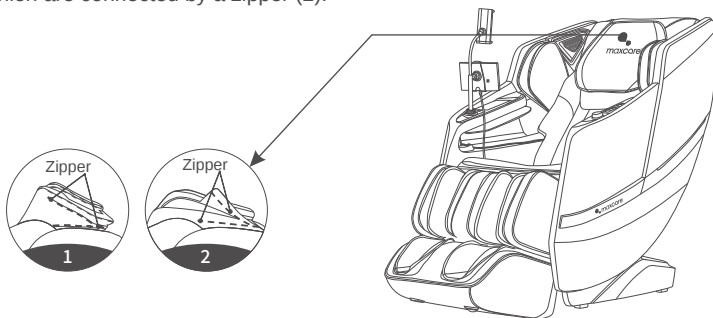
NOTE: Excessive force to tilt will cause the chair to tip over completely and product damage.

1. Before moving, ensure the massage chair is properly powered off, with both the backrest and calf sections fully reset, and the power supply disconnected.
2. The chair should remain empty during movement; no people, pets, or objects should be on it.
3. Do not move the chair over obstacles higher than 2mm or gaps wider than 5mm.
4. When moving the chair using the casters, avoid continuous movement over distances greater than 50 meters. If necessary, move the chair slowly and at a steady pace.

How to disassemble and use headrest cushion and backrest cushion

Using a headrest cushion can reduce the intensity of kneading massage on the neck and shoulders, and you can decide whether to use a headrest cushion according to your needs (recommended to use a headrest cushion). The backrest cushion and the backrest are connected by a zipper (1), and the headrest cushion and the backrest cushion are connected by a zipper (2).

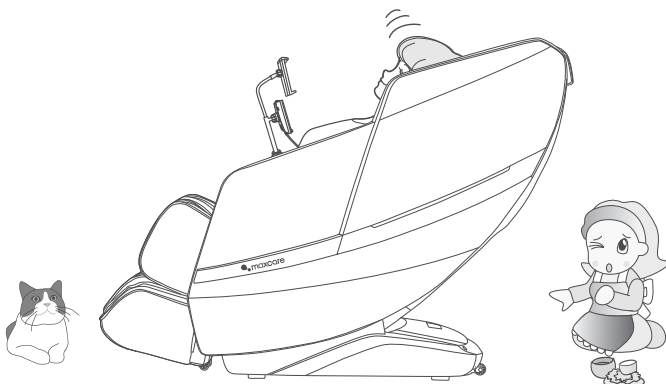
Safety Precautions

Before use

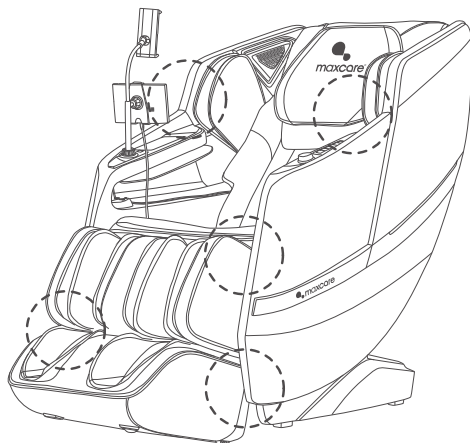
Check around

Make sure that there are no people, pets and other objects around the unit.

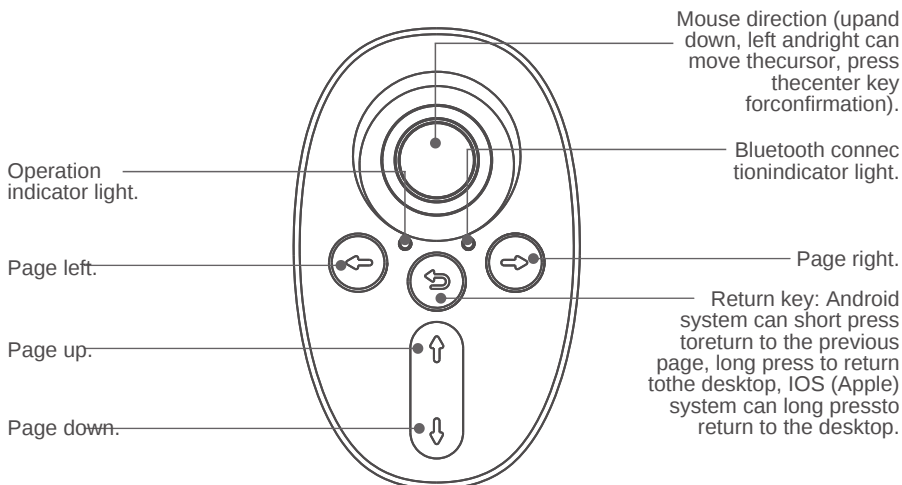
Make sure there are no foreign objects in the gaps of the machine

Before massaging, make sure that no foreign objects are caught in the legs, feet, etc. Under no circumstances allow hands, feet, and head to penetrate into the area between the lower leg mechanism and main body, or the area between the armrest and the capsule. Under no circumstances should hands, feet, or head sink under the fabric or leather covering of the massage area.

Name and Function of Components

Mobile phone bluetoothcontroller operation instructions

Supported systems: Android, IOS (users who want to display the cursor need to turn on Assistive Touch in System Settings-Accessibility-Touch-Assistive Touch).

- 1.Apple users recommend mode 1.
- 2.Android users recommend mode 2.

Combination key usage

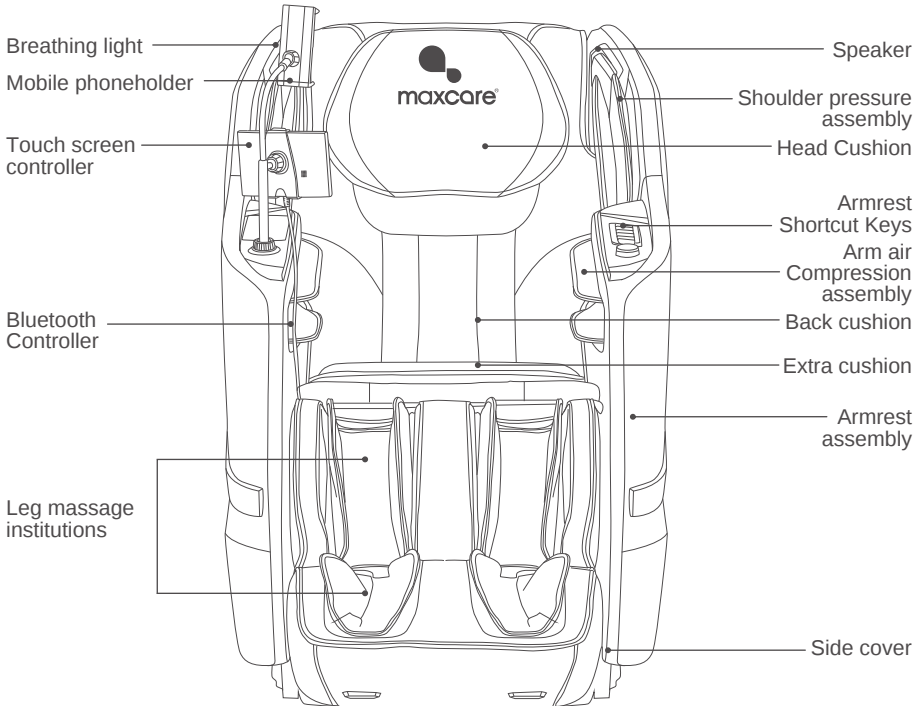
- Cursor speed +: long press the return key + page up (support power outage save).
- Cursor speed -: long press the return key + page down (support power outage save).
- Switch Apple system user: Return key + left key.
- Switch Android system user: Return key + right key.
- Disconnect the currently connected device: Up key + Down key, (the previously paired device is invalid).

Note: This mobile phone Bluetooth controller needs to be connected with mobile 10

Name and Function of Components

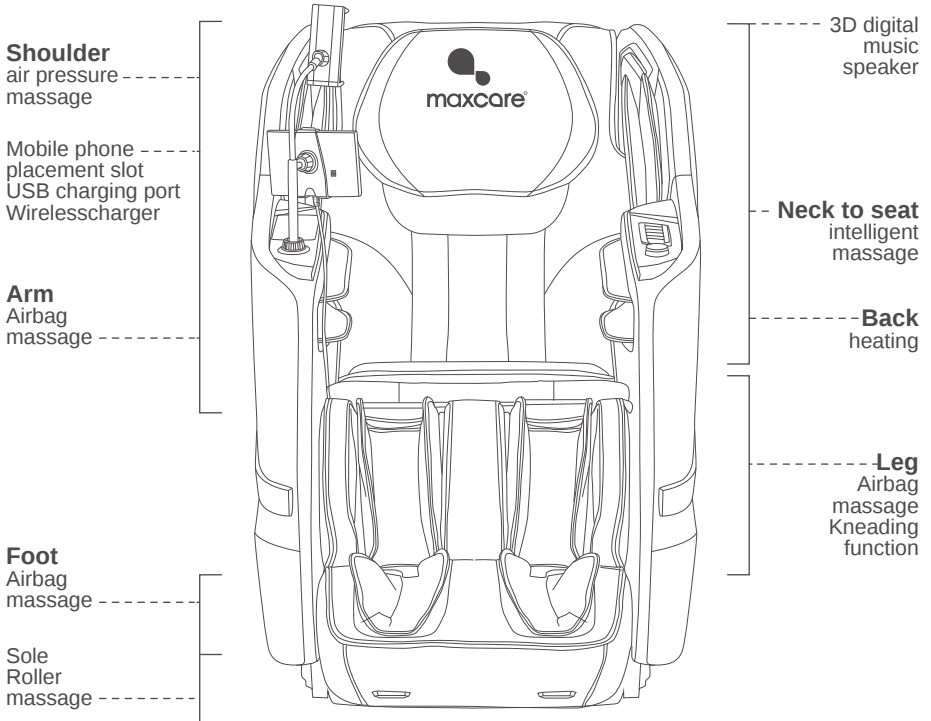
Parts description

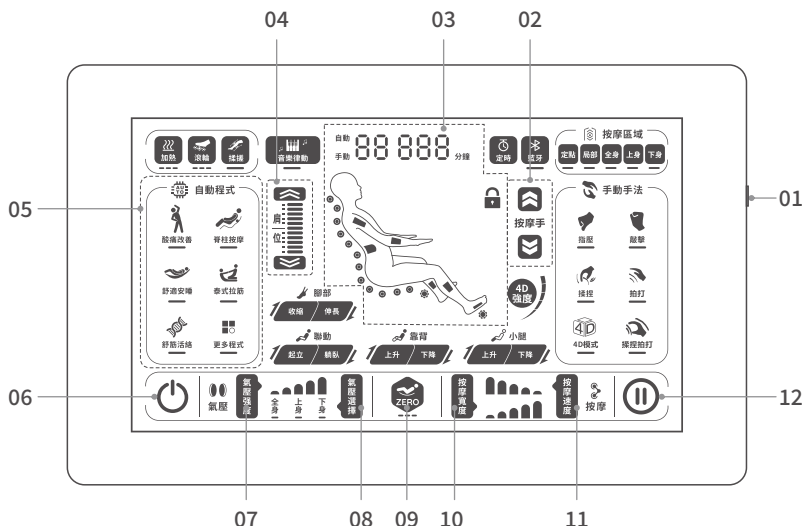
Name and Function of Components

Function distribution

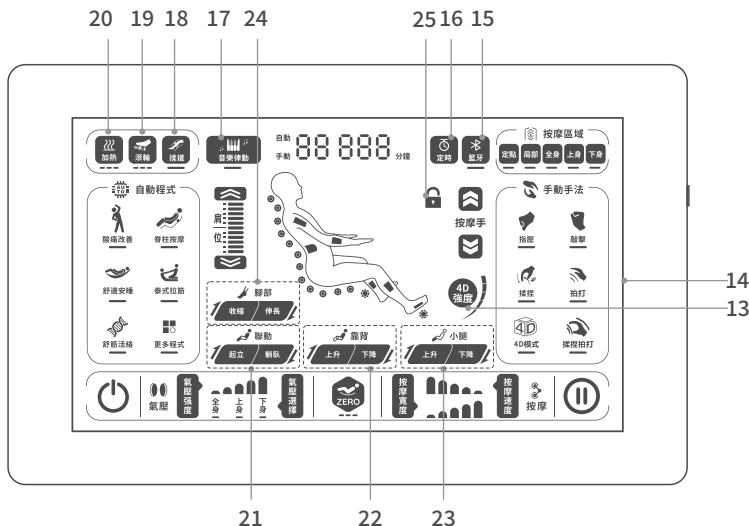
Name and Function of Components

I.Introduction of the remote controller

01. Power Switch Button
Turn the massage chair on or off.
02. Mechanism Up/Down Adjustment Key
Adjust the mechanism position upward or downward.
03. Information Display Screen
Displays the status of the massage chair (program in use / remaining time / mechanism position / air bag / heating / rollers, etc.).
04. Shoulder Position Adjustment Key
Adjust the shoulder position to set the massage reference.
05. Automatic Program Selection Key
Function for selecting automatic massage programs.
06. Touch Power Switch Key
Turn the massage chair functions off.
07. Air Pressure Intensity Key
Adjust the air pressure intensity (5 levels available).
09. Zero Gravity Key
Adjust to Zero Gravity (4 levels available).
10. Massage Speed Adjustment Key
Adjust the massage speed (5 levels available).
11. Massage Area Key
Select massage area function (5 area modes available).
12. Pause Key
Short press to pause the massage or exit the paused state. Press and hold to lock the screen function or exit the locked state.

Name and Function of Components

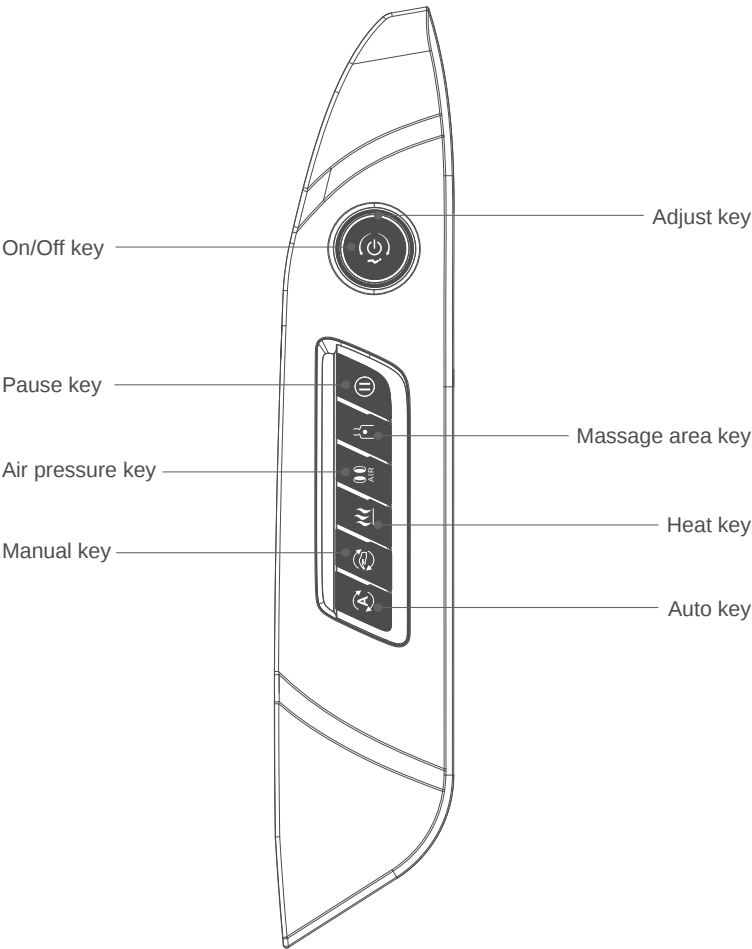
13. 4D Strength Key
Adjust the 4D intensity.
14. Massage Technique Selection Key
Select from different massage techniques (9 options available)
15. Bluetooth turn on/off key
Turn Bluetooth on or off.
16. Time Adjustment Key
Adjust the massage time.
17. Music rhythm keys
Turn on Off the unique automatic program Music rhythm function.
18. Claf Knead Key
Turn the foot calf knead function on or off.
19. Foot Roller Key
Turn the foot roller function on (3 speed levels available) or off.
20. Heating Function Key
Turn the heating function on or off.
21. Sitting Position Adjustment Key
Manually adjust the risen or fall of the backrest and calf support.
22. Backrest Adjustment Key
Manually adjust the rise or fall of the backrest.
23. Calf Adjustment Key
Manually adjust the risen or fall of the leg support.
24. Foot Adjustment Key
Manually adjust the shorten or extend of the foot.
25. Lock Display
Display whether the Screen function is locked or not.

Name and Function of Components

II. Introduction to Armrest Shortcut Keys

Name and Function of Components

Button	Description
	Power key: Short press this key to turn ON/OFF the massage.
	Linkage lying knob: By rotating this knob clockwise, the backrest will move downward slowly and the legrest will rise slowly; by releasing it, such DOWN/UP action of the backrest and the legrest will be stopped.
	Linkage sitting knob: By rotating this knob counterclockwise, the backrest will rise slowly and the legrest will move downward slowly; by releasing it, such UP/DOWN action of the backrest and the legrest will be stopped.
	Press this key during massage, and all the massage actions will be suspended at once. (i.e. the massage chair will turn off after 20min). Then press this key again to restore all the previous massage actions.
	Massage area key: press this button to switch the back stretch mode.
	Air pressure key: Each time the button is pressed, it switches to an air pressure mode, with a total of four modes (including turning off air pressure), which can be switched cyclically.
	Heat treatment key: Click the key to enable or disable the heat treatment function.
	Manual key: each time the key is short pressed, it will be switched to next technique massage function.
	Auto mode key: each time the key is short pressed, it will be switched to next automatic massage function.

How to Use

III. Start massage

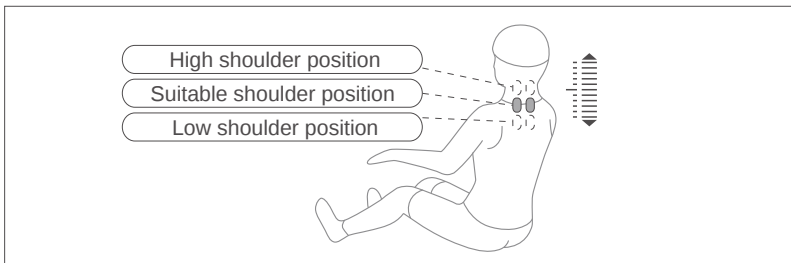
Power Switch Button
(Press this button to
power on/off)

(Touch screen)

1. Press the power button to start the massage, which defaults to Relax massage
2. Automatic angle adjustment (you can press the Handheld controller display “” key/ Armrest shortcut “” key)
3. When detecting the shoulder position, please maintain a proper sitting posture (you can press the Handheld controller display “”/ Armrest short cut key “” key to skip).
4. Fine-tune the shoulder position; if there is a discrepancy between the detected position and the actual position, you can adjust it with the up and down keys. There are a total of 11 shoulder positions, and if no adjustment is made within 10 seconds, it will automatically start the massage function.

Note: You can directly select your preferred automatic mode; pressing the same automatic mode twice will skip steps 2, 3, and 4 and directly start the massage function.

Button	Description	Display
	Adjust the position of the massage ball upwards to reaches the proper shoulder position.	
	Adjust the position of the massage ball downwards to reaches the proper shoulder position.	

How to Use

IV. Pause/Lock Function (please use this function after the chair is powered on)

1. After turning on the massage chair, pressing the pause key  will enter pause mode, stopping all massage functions. Pressing the pause key  again will resume the previous massage state.

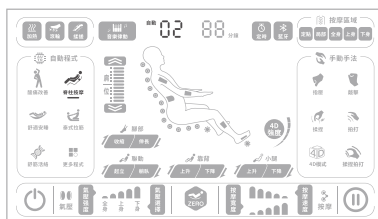
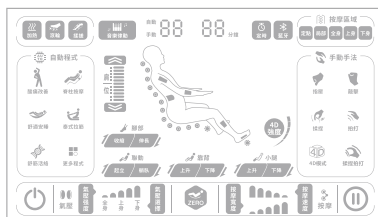
Note: In pause mode, all functions except for turning off the chair are unavailable.

2. Press and hold the pause key  for 3 seconds to lock the screen function. The massage function will continue. Press and hold the key  again to unlock the screen, restoring all massage function keys, and the lock icon  will disappear from the screen.

Note: In lock mode, only the armrest shortcut keys and the handheld remote button can be used; others cannot be operated.

V. Automatic Program Function

1. Select your desired automatic program from the left automatic program menu bar. Once selected, the red box under the automatic program will light up, indicating that a specific automatic program has been selected, while the corresponding number on the information display screen will show which automatic program is running.

2. In the automatic program selection interface, there are a total of 26 automatic programs to choose from: Pain Relief, Spine, Comfort, Thai Massage, Meridian Treatment, Neck And Shoulder, Waist, Bum Toning, Leg&foot Relieve, Executive, Office Regimen, President Cosy, Full-body Stretch, Body Shaping, Circulation Boost, Lactic Acid Breakdown, Sport Recovery, Refresh, Noon Recharge, Vitality Recovery, Relax, Quick Refresh, Tui-Na and Heat, Brain Refresh, Gentle waves knead, Musical Rhythm etc.

How to Use

Function	Key	Display 1	Description	Display 2
自動程式	 酸痛改善	 酸痛改善	Pain Relief: Relieve fatigue and soreness all over the body, relax the limbs and improve micro-circulation, relieve joint pain and restore vitality.	自動 01
	 脊柱按摩	 脊柱按摩	Spine: Care for the health of the spine and improve the lack of exercise.	自動 02
	 舒適安睡	 舒適安睡	Comfort: With kneading, shiatsu massage techniques, through the 12 health acupoints massage, help to regulate and balance the human body; At the same time, the massage chair automatically swings back and forth to form a zero-gravity cradle massage state so that people can fall asleep more quickly, improving sleep quality.	自動 03
	 泰式拉筋	 泰式拉筋	Thai Massage: It mainly uses electric push rods on the calf and backrest to stretch and massage your legs and shoulders with air pressure. The 3D massage hand moves along the back and waist to create deep pressure, imitating being in a Thai massage. It is suitable for people with stress and can improve hunchback.	自動 04
	 舒筋活絡	 舒筋活絡	Meridian Treatment: Mainly massage the shoulders, neck, back and waist with kneading, tapping and acupressure to dredge the meridians.	自動 05
	 更多程式	 更多程式	Neck And Shoulder: The unique massage techniques for the shoulders and neck can help promote blood circulation, relieve neck and shoulder pain, and improve blood supply to the brain.	自動 06
			Waist: Mainly massage the waist with kneading, kneading and tapping techniques, and then do 3D massage on the waist with Tui-na massage techniques.	自動 07
			Bum Toning: Massage focus on the buttocks and thighs and aims to eliminate fat by body sculpting.	自動 08
			Leg&foot Relieve: Mainly massage the shoulders and waist with kneading, tapping, kneading and tapping techniques. Then massage and soothe the legs and feet with air pressure. It is suitable for sedentary people, it can relieve cervical and lumbar pain, as well as legs and feet.	自動 09

How to Use

Function	Key	Display 1	Description	Display 2
 自動程式	 更多程式	 更多程式	Executive:Mainly massage the shoulders, neck, back and waist with kneading,tapping,kneading and tapping techniques. Itcan relieve shoulder and neck pain through deep massage,improve waist muscle soreness, and protectthe lumbar spine.It is suitable for sedentary people.	自動 10
			Office Regimen:Designed for sedentary people. It helps to improve and increase physical activity and improve health.	自動 11
			President Cosy:Mainly massage the shoulders, neck and back with kneading,tapping,kneading and tapping techniques.It helps to relieve fatigue, promote blood circulation and release stress through deep tissue massage on shoulders.neck and back.	自動 12
			Full-body Stretch:Turn on the shoulder airbags, and slowly slide the backrest. At the same time, turn on the leg air bags.Lower the calf slowly,and massage hands massage to the waist. The airbags slowly inflate and deflate, and they massage to the waist. The airbags in the arms will turn on, the massage hands will massage with tapping from the waist to the thighs, and then massage back from the thighs to the shoulders and necks. Loop action.	自動 13
			Body Shaping:Mainly massage the shoulders, neck, waist and buttocks with kneading and shiatsu techniques. To betreated like a queen, to relax your body and mind with a gentle massage, which is specifically designed for women.	自動 14
			Circulation Boost:Dredge the meridians to improve the body'smicrocirculation and enhance immunity.	自動 15
			Lactic Acid Breakdown:Mainly massage the back and waist with tapping,kneading and tapping techniques. It helps to promote blood circulation, eliminate muscle soreness, and accelerate the lactic acid breakdown. Suitable for peoplewith muscle strain and manual labor.	自動 16
			Sport Recovery:Pressing and rubbing with the airbags, and rolling to relieve the soreness of shoulders, arms, knees, calves,and other parts after exercise or fitness.	自動 17
			Refresh:The massage hand moves to and fore from the shoulders and neck to the waist with kneading and tappingtechniques.Then it changes to tapping and knocking massage techniques.	自動 18

How to Use

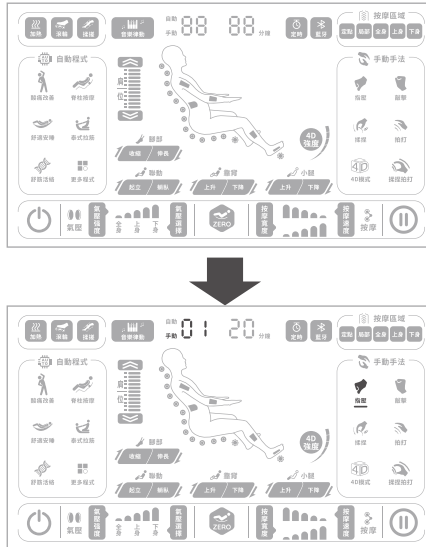
Function	Key	Display 1	Description	Display 2
 自動程式	 更多模式	 更多模式	Noon Recharge: Massage hands,airbags, and rollers are implemented at the same time, you may relax and enjoy with Bluetoothmusic.	自動 19
			Vitality Recovery: Mainly massage the waist, shoulders, and neck with kneading, tapping, kneading and tapping techniques, the deep tissue massages the whole body. This massage is more powerful and can wake up the vitalityof the body and stimulate physical potential. It issuitable for athletes.	自動 20
			Relax: Gentle and comfortable deep tissue massage. It can soothe the tired body and relax.	自動 21
			Quick Refresh: This is a relaxing and most comfortable massageexperience.	自動 22
			Tui-Naand Heat: Kneading massage for the shoulder and neck with airbags, kneading and tapping massage for the back and waist, kneading for the buttocks, and kneading and tapping massage for the thighs. At the same time, thefoot roller, the leg and arm airbags, and the back heat therapy turn on.	自動 23
			Brain Refresh: Mainly massage the shoulders, neck, back, and waist with kneading, 3D kneading, and shiatsu techniques. Toachieve harmony between mind and body through deeptissue massage.Suitable for people who are depressedand have brain overload.	自動 24
			Gentle waves knead: Mainly based on full air pressure massage, using theairbag kneading and pressing in a cycle. To promoteblood circulation, relieve pain and muscle fatigue.	自動 25
	 音樂律動	 音樂律動	Musical Rhythm: The massage techniques mainly focus on kneading, tapping, and a combination of both. By incorporating changes in the rhythm of the music, the intensity and frequency of the techniques are also adjusted accordingly to achieve the best relaxation and soothing effects, enhancing the comfort and effectiveness of the massage experience.	自動 26

How to Use

VI. Manual Program Function

1. Manual Program Selection

- Select your desired manual program from the right manual program menu bar.
Once selected, the red box under the manual program will light up, indicating that a specific manual technique has been chosend.

- In the manual program selection interface, there are 9 manual programs available: Shiatsu, Tapping, Kneading, knocking, Kneading & Tapping, 4D Tapping, 4D Kneading, and 4D Kneading & Tapping.

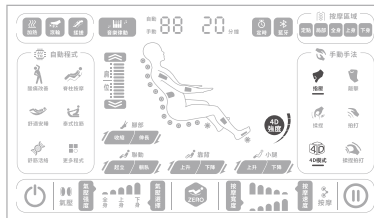
Function	Key	Display 1	Description	Display 2
手動手法	指壓	指壓	Shiatsu: Width.	手動 01
	敲擊	敲擊	Knocking: Speed, Width.	手動 02
	揉捏	揉捏	Kneading: Speed.	手動 03
	拍打	拍打	Tapping: Speed, Width.	手動 04

How to Use

Function	Key	Display 1	Description	Display 2
 手動手法	 揉捏拍打	 揉捏拍打	Kneading & Tapping: Speed.	手動 05
	 指壓  敲擊	 4D模式  指壓	4D Shiatsu: Width.	手動 06
		 4D模式  敲擊	4D Tapping: Speed, Width.	手動 07
		 4D模式  揉捏	4D Grasp: Speed.	手動 08
		 4D模式  拍打	4D Tap: Speed, Width.	手動 09
		 4D模式  揉捏拍打	4D Kneading & Tapping: Speed.	手動 10

2. 4D Intensity Adjustment

The 4D intensity key can adjust the 4D intensity, with a total of 5 levels adjustable; level 1 represents the weakest intensity, and level 5 represents the strongest. The higher the 4D intensity, the higher the massage rollers protrude and the greater the massage force. Pressing the key once changes the 4D intensity from weak to strong, cycling through the 5 levels.

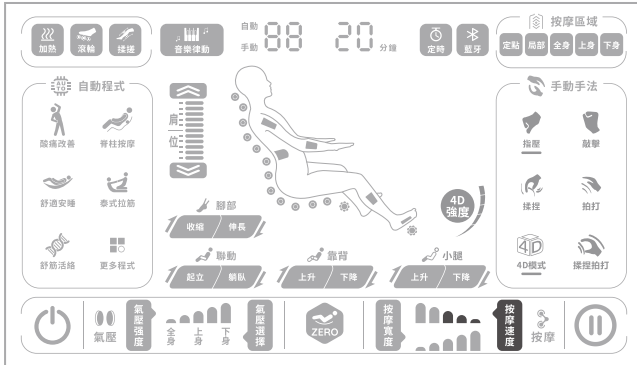
Button	Description				
	Each time a Button is pressed, the 4D intensity changes from weak to strong. There are a total of 5 levels, cycling through them.				
 1 Level	 2 Level	 3 Level	 4 Level	 5 Level	

How to Use

3. Massage Speed

- (1) The user must select any manual program before adjusting the massage speed.
- (2) The massage speed can be adjusted across 5 levels, with level 5 being the fastest. Speed adjustments are not available when the massage techniques are paused or during shiatsu. Each press changes the massage speed from slow to fast, cycling through the 5 levels.

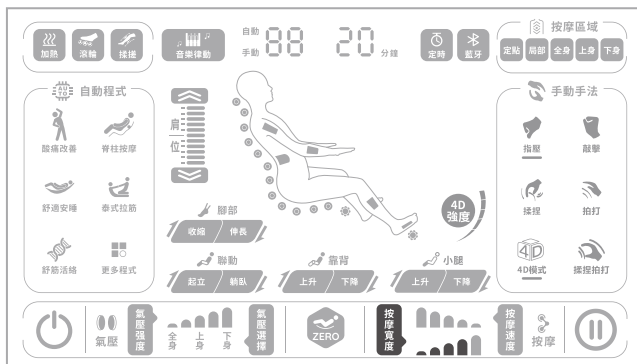
Button	Description				
	Each time key is pressed, the speed changes from slow to quick. There are a total of 5 levels, cycling through them.				
	1 Level		2 Level		3 Level
			4 Level		5 Level

How to Use

4. Massage Width

- (1) The user must select any manual program before adjusting the massage width.
- (2) The massage width has 5 adjustable levels, with level 5 representing the widest. Please refer to page 21-22 for specific manual programs that allow width adjustment.

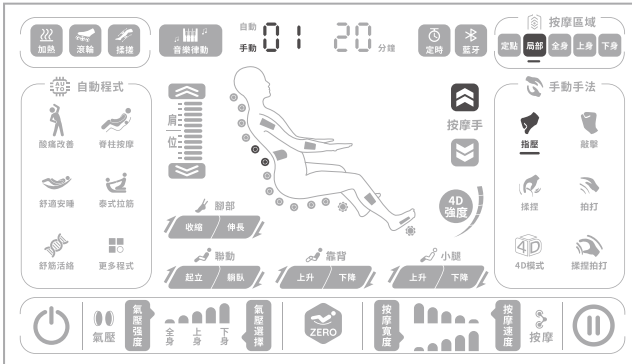
Button	Description				
	Each time key is pressed, the width changes from narrow to wide. There are a total of 5 levels, cycling through them.				
	1 Level		2 Level		3 Level
	4 Level		5 Level		

How to Use

5. Massage Area

There are 5 selectable massage areas: Fixed Point, Local, Upper Body, Lower Body, and Full Body. Each press of the Button switches the area massage mode.

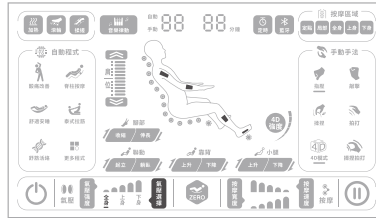
Function	Key	Display 1	Description	
按摩區域	定點	定點	Fixed Point: Targeted massage at a specific point.	
	局部	局部	Partial: Small area back-and-forth massage.	
	全身	全身	Full Body: Full back back-and-forth massage.	
	上身	上身	Upper Body: Upper back back-and-forth massage.	
	下身	下身	Lower Body: Lower back back-and-forth massage.	
按摩手		↑	Manipulator uplink	Massage position adjustment (in fixed point and local mode, massage position can be adjusted).
		↓	Manipulator descending	

How to Use

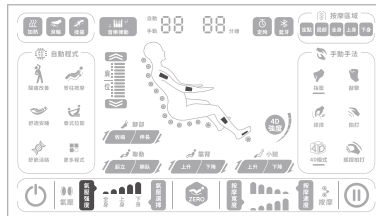
VII. Air Pressure Massage

- (1) Air Pressure Selection: Options include upper body air pressure, lower body air pressure, full body air pressure, and turning off the function. Each press key the air pressure mode (including turning off air pressure, with a total of 4 modes cycling). Once selected, the red box under the mode will light up, indicating a selected massage area as following:

Button	Description	Display
氣壓選擇	Full Body Air Pressure: Air pressure massage for arms, shoulders, legs, and feet.	全身 上身 下身
	Upper Body Air Pressure: Focused air pressure massage for arms and shoulders.	全身 上身 下身
	Lower Body Air Pressure: Focused air pressure massage for legs and feet.	全身 上身 下身

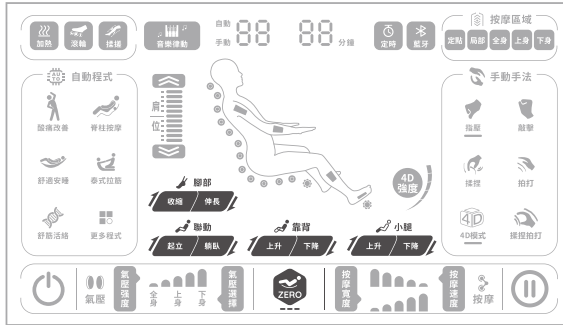
- (2) Air Pressure Intensity Adjustment: After activating the air pressure function, the intensity can be adjusted with 5 levels available; level 1 represents the weakest intensity and level 5 the strongest.

Button	Description				
	Each press changes the air pressure intensity from weak to strong, cycling through the 5 levels.				
					
1 Level	2 Level	3 Level	4 Level	5 Level	

How to Use

VIII. Sitting Position Adjustment

1. Zero Gravity Adjustment

Button	Description	Level	Display
	Each press key to a different zero gravity mode. With a total of 4 modes (including reset standing mode), cycling through them.	Reset Stand	
		1 Level	
		2 Level	
		3 Level	

2. Sitting Position Adjustment

Button	Description
	Linked Standing Key: Pressing this key slowly raises the backrest while lowering the leg support. Release to stop.
	Linked Recline Key: Pressing this key slowly lowers the backrest while raising the leg support. Release to stop.

How to Use

3. Backrest Angle Adjustment

Button	Description
	Backrest Up Key: Pressing this key slowly raises the backrest; release to stop.
	Backrest Down Key: Pressing this key slowly lowers the backrest; release to stop.

4. Leg Support Angle Adjustment:

Button	Description
	Leg Up Key: Pressing this key slowly raises the leg support; release to stop.
	Leg Down Key: Pressing this key slowly lowers the leg support; release to stop.

5. Foot shorten/Foot extend adjustment.

Button	Description
	Remotely adjust shortening of the leg unit; Press and hold this button to move the leg unit up slowly and stop when released.
	Remotely adjust extension of the leg unit; press and hold this button to move the leg unit down slowly and stop when released.

Note:

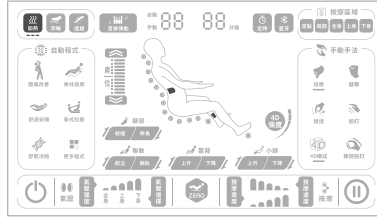
After the Zero gravity adjustment/Leg angle adjustment /Backrest unit angular adjustment/Lying angle adjustment action is completed, the leg unit will automatically detect the foot length.

How to Use

IX. Other Functions (please use after the chair is powered on)

1. Heating Function

- (1). Click the “ 加熱” button · turn on the heating function or switch the heating area. The heating function is already open (After turning on the heating function for 3 minutes, it will be warm). As shown in the following figure:

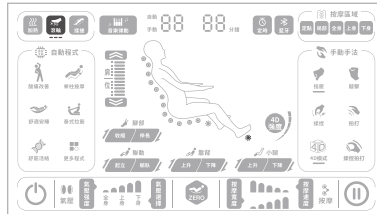
- (2). Click the “ 加熱” button again, heating area can be switched, the sequence is: waist heating, calf heating, waist + calf heating, off heating.

Button	Description	Heating Mode	Display
 加熱	Each time the button is pressed, the heating is adjusted. There are three modes to choose, which can be switched in a cycle.	Stop	 加熱
		Waist heating	 加熱
		Calf heating	 加熱
		Waist + Calf heating	 加熱

2. Foot Roller:

Click the foot roller button  to turn the foot roller massage on or off. It features three adjustable intensity levels—low, medium, and high—as well as an off mode, cycling through four states.

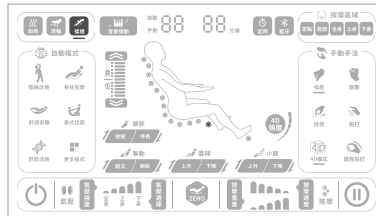
How to Use

Button	Description	Level	Display
	Press the key to turn the foot rollers on or off, with 3 adjustable speed levels, cycling through them.	Stop	
		1 Level	
		2 Level	
		3 Level	

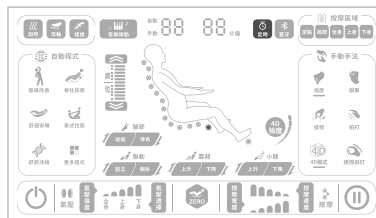
3. Calf kneading:

Click the calf kneading button  to turn on or off the calf kneading and function.

4. Time Adjustment

Click the calf kneading button  to adjust the timer in 5-minute increments, not exceeding 30 minutes. If it exceeds 30 minutes, it resets to the minimum time value.

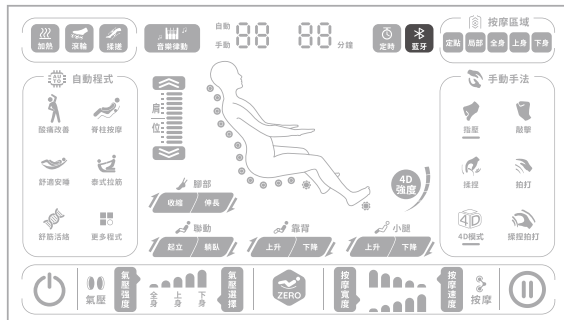
Button	Description	Display
	Press the key to adjust the timer in 5-minute increments, not exceeding 30 minutes. If it exceeds 30 minutes, it resets to the minimum time value.	自動 88 23 分鐘 手動 88 23 分鐘

How to Use

5. Bluetooth

- (1) Press the  key to turn on Bluetooth. The red box under the Bluetooth key will light up, indicating that Bluetooth is active.

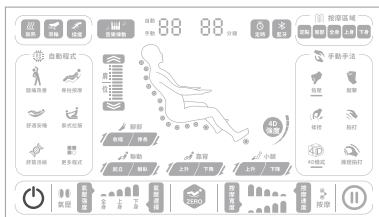
- (2) Enter the main menu of your smart device, go to "Settings" -- select the "Bluetooth" option -- turn on "Bluetooth" -- and the device will automatically "search for new Bluetooth devices."
- (3) Your smart device will find the corresponding Bluetooth model MAX-Z05A369-XXXX (XXXX represents the serial number). Click the pair button. When asked whether to pair with this device, select "Pair" to complete the pairing.
- (4) Open the music player on your smart device and select a song to play.
- (5) Press the  key again to turn off Bluetooth. The red box will disappear, indicating that Bluetooth is off.

How to Use

X. Ending Massage

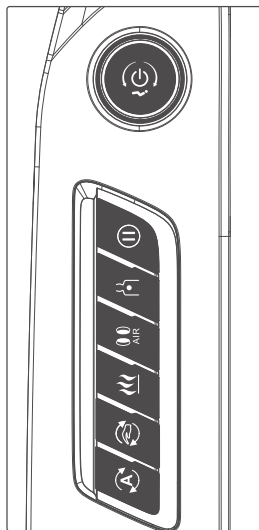
- (1) During the massage, press the power  on the touch screen control or the right armrest power button  to turn off. This will immediately disable all massage functions, reset the backrest and leg support, or if the massage timer runs out, all massage functions will automatically shut off without resetting the backrest and leg support.

(Touch screen)

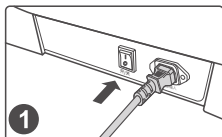
Power Switch Button
(Press this button to
power on/off)

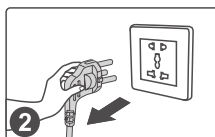
(Armrest short cut key)

- (2) Cut off the main power, as shown in the diagram (illustration for cutting off the main power).

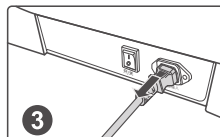
Turn off the power
switch ("O"position)

Unplug the power cord
(illustration for unplugging)

Disconnect the product's
power connection

How to Use

XI. APP download and Installation Instructions

1. For Apple (iOS system): Search for "Maxcare Massage Chair" in the App Store, download the iOS version and install it. Alternatively, scan the QR code below to download and install. Once installed, you can use the Bluetooth APP to operate the massage chair.
2. For Android system: Search for "Maxcare Massage Chair" in Google Play, download the Android version and install it. Alternatively, scan the QR code below to download and install. Once installed, you can use the Bluetooth APP to operate the massage chair.

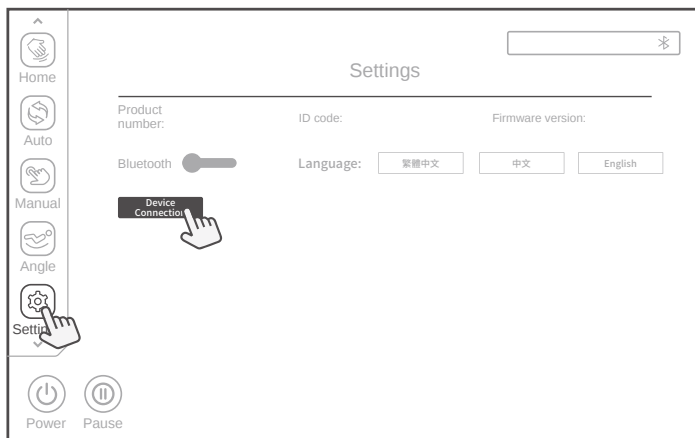
Maxcare massage chair

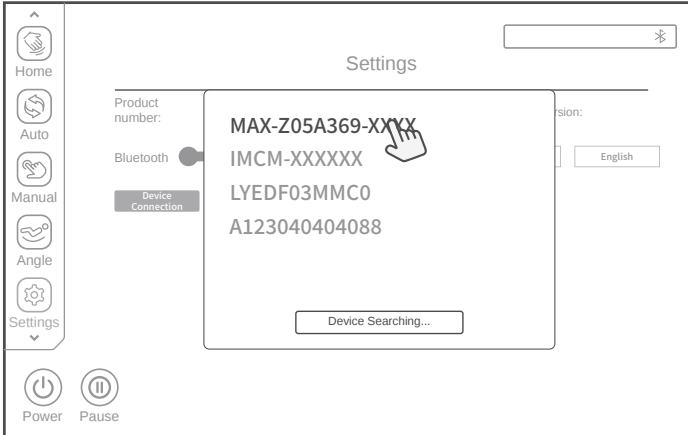
XII. APP Device Connection

1. Select the "⚙️" button in the left menu of the display to enter the settings interface, then select "Device Connection".

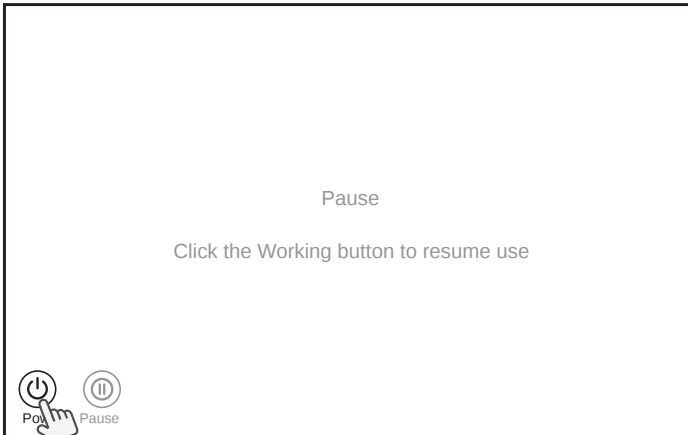
How to Use

2. APP will automatically search for the corresponding device number (e.g., MAX-Z05A369-XXXX). Tap the corresponding device number (it will show "Connected"), the connection will be completed.
3. You can operate the massage chair once the device is connected successfully.

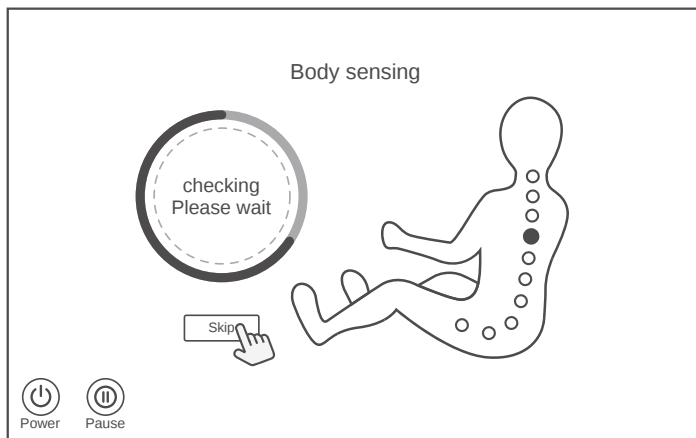
XIII. APP Instructions for controlling the massage chair

1. Press the "⏻" button on the interface to activate the massage function.

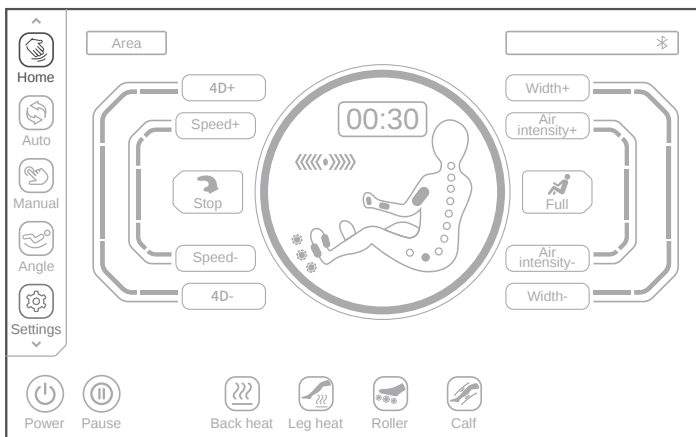
How to Use

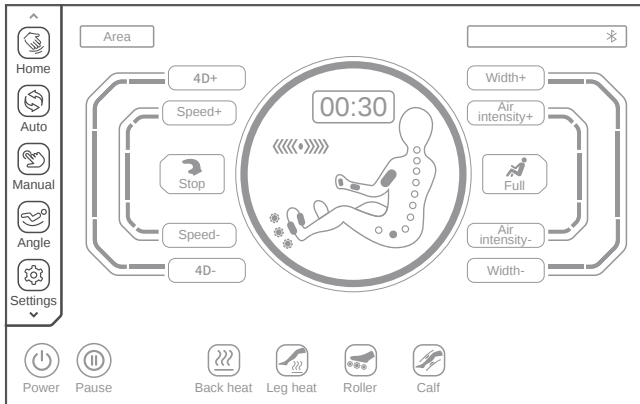
Note: Wait for the body shape detection to complete.

- After activation, enter the main interface and select the corresponding massage functions according to your needs.

How to Use

3. Title bar manual

Interface Diagram

Current Status: Displays the massage chair's operating status and allows operation of the chair.

Auto: Selects various built-in massage programs.

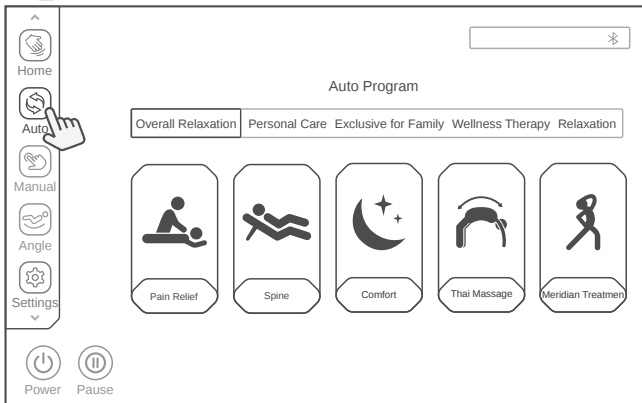
Manual: Manually adjusts various functions of the massage chair.

Seating Posture: Adjusts the seating posture of the massage chair.

Settings: View information such as the massage chair model and serial number, connect the chair to the App, set the software language, etc.

4. Select and Tap the Desired Auto Program to Start the Massage

Select the “” button in the left menu bar to enter the auto selection interface.

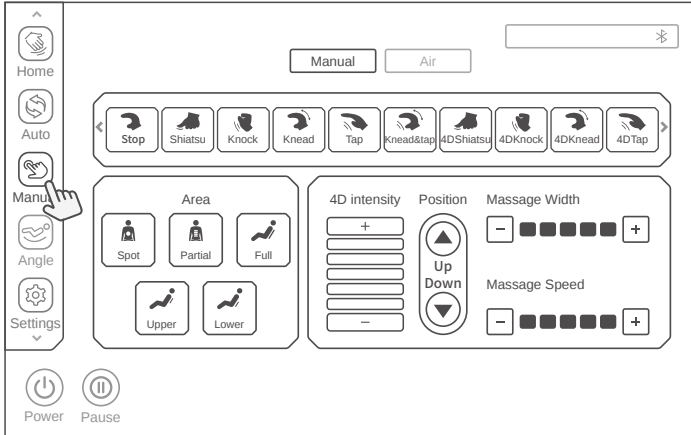
Auto Program Selection Diagram

How to Use

5. Manual Mode

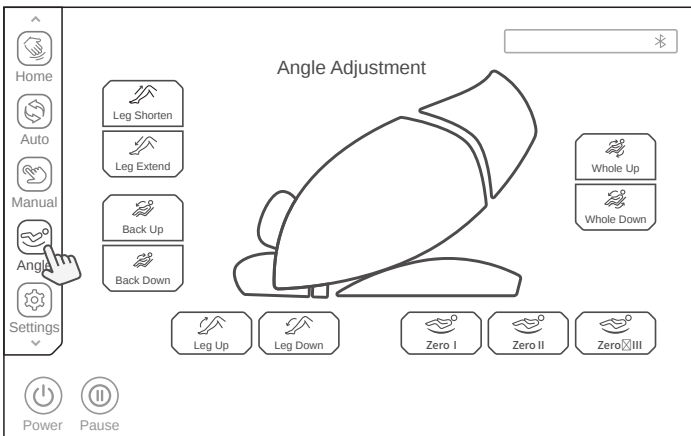
Select the “” in the left menu of the display to enter the manual program selection interface.

Manual Program Selection Diagram

6. Seating Posture Adjustment

Select the “” button in the left menu bar of the display to enter the seating posture adjustment interface. This interface includes functions such as Zero Gravity I/II/III, linked lift, linked lower, foot shorten, foot extend, backrest lift, backrest lower, calf lift, and calf lower.

Seating Posture Adjustment Diagram

Cleaning and Maintenance

Cleaning

- Before cleaning, please make sure the switch has been turned off and the power plug has been disconnected from the socket.
- Use soft cloths to clean the surface of product and control panel.

CAUTION

- Do not pull the plug from the socket with wet hands.
- Please do not rinse the product. Avoid water/ liquids from getting into the product.
- Do not use corrosive cleaning agents like benzene or thinner to clean the product, otherwise, the color of product will fade, and the internal parts may be corrupted or broken.
- Do not disassemble the product when it is humid. Simply let it dry.

Maintenance

- Make sure that the plug has been disconnected from the socket before maintenance.
- Keep the product clean and away from high-temperature/ humid areas, or any environment with volatile gases.
- If this product is to be left idle for long, properly cover it with a dust cover (self-purchase).

CAUTION

- Please keep the product away from direct sunlight, or high temperature, to avoid fading or product damage.

Kindly reminding

- If the product is used through the year, please check conditions as below:
 - * Burnt smell.
 - * Unstable power supply after contact with the power cord.
 - * Power cord heating up.
 - * Power cords damage.
 - * Other abnormalities.

CAUTION

- If above mentioned happens, please stop using the product.
- Switch off the power supply, pull out the plug, and contact us.
- Do not disassemble the product by yourself.

Maintenance and Cleaning of PU Leather

PU leather is made by the coating and drying of polyurethane. Compared with genuine leather, PU leather is more resistant to tear, heat and sunlight. Depending on its nature and chemical substances, however, it should be kept away from many dusty, high/low-temperature, humid environment, under direct sunlight or places with corrosive substance. Otherwise, the leather could tear, quickening aging and hydrolysis.

Cleaning

1. Please do not use organic reagent or grease solution, especially strong solvent, e.g. alcohol, turpentine and thinner, to clean the surface of the PU leather. Otherwise it will corrode the surface.
2. Please do not use water or washing powder to clean the surface of the PU leather. Otherwise there will be cracks.
3. When cleaning the surface of the PU leather, choose a cleaning agent that has a pH value of between 5~7, no dentifrice, and is harmless to the environment.
4. For dirty PU surface, apply egg whites to the surface using a piece of soft cloth. Avoid using a brush for cleaning.

Maintenance

1. Keep a good ventilation indoor, overly dry or humid environment can accelerate the aging of leather
2. Please keep the product away from direct sunlight or air vent, otherwise the leather will be hardened. Direct sunlight will also cause fading.
3. People are generally sweatier in summer. The holes on the surface of the leather will absorb sweat. While heat and high humidity prompt the chemical reaction between certain substances in sweat and the leather, and hence, causes smell. Please clean the surface more often with a piece of lightly wet cloth.
4. Please do not use alkaline cleaning agents on the leather. It will adversely affect the elasticity of the leather and lead to cracks as a result.

? FAQs and Solutions

The following symptoms do not indicate malfunction. Please check the below phenomenon to find out the solution before contact customer services department. If the problem still cannot be solved, please do not dismantle or repair the product. You should stop using it and contact us immediately.

Fault	When the calf mechanism or the backrest reaches a certain position, press the corresponding adjustment button and it will no longer change, and a continuous "beep, beep" sound will be emitted.
Solution	Due to product structure design and ergonomic requirements, this product has a limited protection design for the limit state of attitude adjustment. At the same time, a reminder sound will sound when the limit is reached.

Fault	The function does not work or some functions do not work after the machine is turned on.
Solution	After the machine is turned on by pressing the switch key, it is necessary to manually select an automatic massage program or other working state. If it is not operated, the machine will automatically shut down after 20 minutes. If the manual operation mode is selected, the functions that are not selected may be disabled, and the working mode needs to be manually set in the manual mode.

Fault	The machine makes a leather rubbing sound when adjusting its posture.
Solution	This product adopts high-grade anti-corrosion leather that meets environmental protection requirements. Due to the inherent characteristics of leather materials, there will be a normal slight noise when rubbing at the relative movement, which is a normal phenomenon. Note: If you hear fabric or leather tearing sound when the machine is running, please stop using it immediately and seek professional maintenance.

Fault	The massage wheel did not reach the shoulders or neck.
Solution	If the head does not touch the head cushion or the back does not touch the backrest, the shoulder position may be lower than the actual position during the body shape detection process. Sit in the deepest position of the seat, rest your head on the headrest, and start over from the beginning.

Fault	The unit is damaged. The power cord or power plug is abnormally hot.
Solution	To prevent accidents, be sure to contact an authorized service center.

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Trademark holder: Maxcare Industrial Limited