Phototherapy
Neckcare
頸康寶

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MAX-E03HXRE200

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User Manual

MAX-E03HXRE200

This product jointly developed by MTCURE (镁天健康) and The Chinese University of Hong Kong, and maxcare® provides engineering and manufacturing supports, effectively improves pain conditions while promoting metabolism and blood circulation.

MTCURE (镁天健康), an innovative biotechnology start-up incubated by The Chinese University of Hong Kong.

COMPANY PROFILE

After two decades of experience in the health industry, **maxcare**® founder utilized his knowledge and expertise in massage therapy from an Asian-Pacific perspective to establish our brand - **maxcare**®. Our massage therapy business was expanded to Hong Kong in 2012 and is conferred 2017 Hong Kong Top Brand Award. Based in Hong Kong, we are specializing in design & innovation, manufacturing, and retailing of massage and healthcare products.

Believing in the value of 'Technology improves our quality of life' and 'Caring for your family', **maxcare**® has been innovating highly qualified products with international patents. Besides, **maxcare**® products fulfill various international safety standards, namely EU CE Mark, US UL/ETL Mark, Germany GS Mark, TUV certification, and China CCC Mark.

Motivated by our value, **maxcare**® founder has been leading us to be one of the most competitive brands in Hong Kong health industry. **maxcare**® applies the contemporary fashion design – neat and modern – on the mechanical engineering of massage chairs. To deliver our value and offer incredible massage experience and healthcare solution to the public, **maxcare**® devotes to the innovation, design, manufacturing the top-notch family massage products. By integrating the company philosophy – Spirit Of Healthy Life - into our products, **maxcare**® brings health to everyone.

Phototherapy Neckcare

Thank you for purchasing Maxcare's products. Before using the product, please read this manual thoroughly, particularly the Safety Precautions section, for the correct use of the product.

This manual should be well kept for easy reference. It enables you to perform a safe and efficient operation.

CONTENT

Safety Precautions	P.01 - P.03
Product Introduction	P.04
Name and Function of Components	P.05
How to Use	P.06
Cleaning and Maintenance	P.07
Troubleshooting and FAQs	P.08
Technical Specifications	P.09

Safety Precautions

1. Important safety warning

- Please use a well-grounding power supply suited to this product.
- Please pull out the plug after use or before cleaning to avoid injury or damage to the product.
- Do not use unauthorized accessories.
- Please read this manual instruction carefully before operation.
- Any other usages that are not listed in this manual instruction are forbidden.
- Please check if any part of this product is damaged before use. If so, please do not use it and contact us for repair.
- Please do not use it if the product is broken or damaged.
- Please do not use it if the cover or surface is damaged.
- Please do not drop anything into this product.
- Please do not allow water or liquid into this product.
- Please do not fall asleep while using this product.
- Please do not use it if drunk or feeling unwell.
- Keep the plug clean from needles, rubbish, or water.
- If the power cord is damaged, do not attempt to change on your own. Please contact us for repair for safety reasons.

2. Environment for usage

- Please do not use it in a moist environment such as a bathroom.
- Please do not use it when there is a drastic change in temperature.
- Please do not use it under a heavy-dusty or caustic environment.
- Please do not use it where there is not enough space or not a good ventilation.
- Please switch off the power after use.
- Please do not use it when the surrounding temperature is up to 50°C.
- Turn off and pull out the plug during thunderstorms.
- Do not place it near fire or directly under the sun.
- When moved from a low-temperature storage to a warm environment, the moisture inside the machine may affect, or damage, its functionality. Please rest the product within room temperature for 1 hour before use.

Safety Precautions

3. People who have the below symptoms, please consult a doctor before use

- Children under the age of 13 and pregnant women should not use the light therapy product;
- Cancer patients;
- Tuberculosis patients;
- Acute suppurative inflammation;
- Hemorrhagic patients;
- Heart disease and pacemaker patients;
- Arteriosclerosis patients;
- Newborns and infants;
- Patients who are allergic to light or have skin that is insensitive to heat are considered in appropriate to use these products by doctors.

CAUTION

• Please note this is not a medical device.

4. Precautions before using the product and some safety instructions

- Do not use the Heat functions on area that have implanted metal or other materials, as well as areas with inflammation. It is recommended to use Photobiomodulation Treatment (PBMT) instead and keep the product at a distance of 3-5mm from the skin surface.
- Avoid using the Heat functions on area with wounds, pimples, or suppurative diseases. Once the wounds have stopped bleeding or the tissue sprain has healed, it is advisable to use only Photobiomodulation Treatment (PBMT) to promote faster healing. After 72 hours or when the wound has fully healed, it is safe to use the machine.
- Do not directly view the light. (Included wear sunglasses or protective spectacles)
- Please do not place heavy objects on this product.
- Do not use other therapy devices while using this product.
- Please note this is not a medical device, it is not to be used for medical treatment purposes. Only use for improvement and relief purposes.
- Recommended usage time of this product is 20 minutes/ session.
- Check if the voltage is suitable for the specification of this product.
- Do not use wet hands to pull the plug and use the machine.
- Do not allow water to get into the product, for the avoidance of electric shock or damage to this product.
- Do not pull or drag the power cord when unplugging. Please be gentle.
- Do not damage the wires or change the circuit and alter the structure or parts of this product.

Safety Precautions

- Do not use wet cloths to clean the electric parts such as the switch and plug.
- Stop using this product if it appears any irregularity, and please contact us immediately.
- Stop using this product if it causes discomfort, and please consult a doctor.
- This appliance is not intended for people (including children) with reduced physical, sensory or mental capabilities, or a lack of experience and knowledge, unless they are under supervision or instruction concerning the use of the appliance by a person responsible for their safety.
- Children should be supervised to ensure that they do not play with the appliance.
- If the power cord is damaged, it must be replaced by us in order to avoid a hazard.
- Keep the power cord away from hot or moving objects.
- Please cut off power and pull out the plug before moving.
- If the product is found damaged, broken, suffering from power leakage, or with exposed parts, please stop using it.
- This product contains a lithium-ion rechargeable battery.
- Do not place this product near fire or where it causes fire easily to avoid accidents; the product may also deform, fade, or break.

5. Maintenance

- This product should be maintained by us, users are prohibited to disassemble or maintain by themselves.
- If this product is to be left idle for long, please keep it in a dry and dustless environment.
- Do not keep this product at a high temperature or near a fire. Do not keep under direct sunlight.
- Please clean this product with dry cloths. Do not use thinner, benzene, or alcohol.
- Avoid contact with sharp edges or pointed objects which might cut or puncture the fabric surface. To avoid breakage.
- Do not fold during storage and use to avoid damage to the lamp beads.
- Please use it intermittently; do not keep it continuously running for long.

 **Product Introduction**

With the increasing pressure of life and work, people's bodies will have various problems. Phototherapy Neckcare is a new treatment that can be used in home or anywhere. It is a safe, efficient, reliable, painless, drug-free, cordless, and natural solution for body pain relief. It can transmit the Heat and Photobiomodulation Treatment (PBMT) to the cervical spine, and directly stimulate muscles, subcutaneous tissue, etc. It can be used separately or combined at the same time.

*This is not a medical device. This unit is intended to be used only for general physical conditioning.

Product Function

Heat	Hot compress can help to relax tired and sore muscles 42°C(+/- 3°C) constant temperature
MT-PBMT (Multi- Photobiomodulation Treatment)	Phototherapy combines red light in the range of visible light and invisible light, starting at superficial tissue to deeper tissue levels, it promotes wound healing and tissue repair, relieves pain and inflammation.

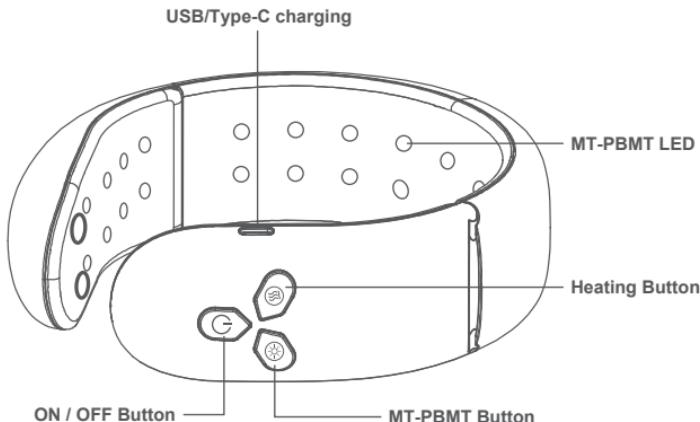
* MTCURE (MTcure Biotech Limited) biophysical combined therapy can effectively enhance the effects of various biological treatment signals, thereby avoiding side effects.

 **Name and Function of Components**

Phototherapy
Neckcare

USB/Type-C charging cable

 **How to Use**

	Power button: long press to turn on the device, heating and red lighting will turn on simultaneously
	Hot compress button: Short press to turn the heat off and back on again. 42°C (+/- 3°C) constant temperature
	Photobiotherapy button: short press to turn off the red light function, or turn the red light on again

Usage instructions

Long press the on/off button for 2 seconds to turn on the device. The power button will display a green light. The phototherapy neckcare will beep once and enter the red light mode, and the constant temperature hot compress function will be activated. If the power button shows a red light, it indicates low battery. Please charge it before use.

Device charging

If it is going to be used for the first time or has not been used for a long time, please charge it as per given instructions:

1. The USB cable in the accessory should be plugged into the charging socket on the bottom left side of the product.
2. Insert the USB cable into the power adapter and connect the power.
3. After charging begins, the product indicator light will remain constantly red. It takes about 2 hours to fully charge. Once fully charged, the charging indicator light will turn green and remain constantly on. After charging for a period of time, the indicator light will go out and the charging will automatically stop. Please unplug the power supply in time after it is fully charged.
4. If the red light on the power button of the phototherapy neckcare flashes during use, it indicates that the battery is low. Please charge the device in time.

*Note: Do not use this device during charging.

Cleaning and Maintenance

Cleaning

- Please ensure to turn off the power and unplug before cleaning.
- Use a soft dry cloth or slightly damp cloth to wipe the product surface.

CAUTION

- Do not touch the plug with wet hands.
- Do not wash or spill any liquid into the product.
- Do not use corrosive liquid such as benzene or thinner to clean the product, to prevent color fading and damage.
- When the product is damp, stop using it and do not dismantle. Leave it to dry naturally.

Maintenance

- Keep the whole set of the product clean and storage in a dust free and dry location under the condition. Storage temperature and humidity: -10oC to 60oC, 10% to 95% RH
- When storing the phototherapy neckcare, please avoid folding it or applying excessive pressure to prevent internal wire breakage and damage to the lamp beads.

CAUTION

- Do not place it near fire or direct sunlight, to prevent product color fading.

Kindly Reminder

- If the product has been used for a long time, please check the conditions as below:
 - * Burnt smell.
 - * Unstable power supply after connecting with the power cord
 - * Power cord overheat
 - * Power cords damage
 - * Other abnormalities

CAUTION

- If the above phenomenon occurs, please stop using the product to prevent malfunctions and accidents.
- Please turn off the power and unplug immediately. Users are prohibited to dismantle and repair by themselves. Please contact us for repair.

 **Troubleshooting and FAQs**

The following symptoms do not indicate malfunction. Please check the below phenomenon to find out the solution before contact customer services department. If the problem still cannot be solved, please do not dismantle or repair the product. You should stop using it and contact us immediately.

Issues	Possible Causes	Solutions
No Response The device is blinking	- Low battery	- Recharge the product
All buttons do not work well	- Low battery	- Recharge the product

 Q: After use, symptoms such as redness, inflammation, and itching may occur.
A: Individual responses to treatment may vary. If the symptoms do not improve the next day, please consult a healthcare professional.

 Q: Feeling discomfort and fatigue after use.
A: Individual responses to treatment may vary. Before the discomfort subsides, please temporarily reduce the frequency of use.

Technical Specifications

Product name	: Phototherapy Neckcare
Model no.	: MAX-E03HXRE200
Charging voltage	: 5V=-1.5A
Battery capacity	: 3.7V 1900mAh
Rated power	: 8W
Photon	: 30 pcs of 660nm red light (visible light) 90 pcs of 940nm near-infrared (invisible light)
Charge time and run time	: Charging for <2.5 hrs allows for continuous operation for 1-3 hrs (Depending on different modes of heating/red light)
Auto Timer	: 20 minutes
Product Size	: 156mm x 146mm x 46mm
Product weight	: 171g

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MTCURE BIOTECH LIMITED an innovative biotechnology start-up
incubated by The Chinese University of Hong Kong.