360

ShapeLift

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MAX-W02YC212W

maxcare®

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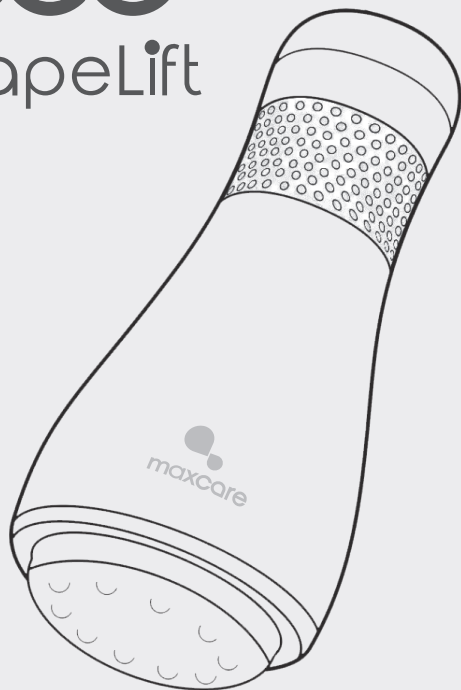
MAX-W02YC212W

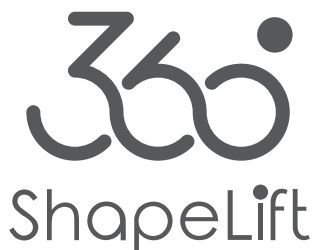
COMPANY PROFILE

After two decades of experience in the health industry, **maxcare** founder utilized his knowledge and expertise in massage therapy from an Asian-Pacific perspective to establish our brand - **maxcare**. Our massage therapy business was expanded to Hong Kong in 2012 and is conferred 2017 Hong Kong Top Brand Award. Based in Hong Kong, we are specializing in design & innovation, manufacturing, and retailing of massage and healthcare products.

Believing in the value of 'Technology improves our quality of life' and 'Caring for your family', **maxcare** has been innovating highly qualified products with international patents. Besides, **maxcare** products fulfill various international safety standards, namely EU CE Mark, US UL/ETL Mark, Germany GS Mark, TUV certification, and China CCC Mark.

Motivated by our value, **maxcare** founder has been leading us to be one of the most competitive brands in Hong Kong health industry. **maxcare** applies the contemporary fashion design – neat and modern – on the mechanical engineering of massage chairs. To deliver our value and offer incredible massage experience and healthcare solution to the public, **maxcare** devotes to the innovation, design, manufacturing the top-notch family massage products. By integrating the company philosophy –Spirit of Healthy Life - into our products, **maxcare** brings health to everyone.

Thank you for purchasing Maxcare's products. Before using the product, please read this manual thoroughly, particularly the Safety Precautions section, for the correct use of the product.

This manual should be well kept for easy reference. It enables you to perform a safe and efficient operation.

CONTENT

Safety Precautions	P.01 - P.03
Name and Function of Components	P.04
How to Use	P.04
Charging Method	P.05
Cleaning and Maintenance	P.06
Troubleshooting	P.07
Technical Specifications	P.08

Safety Precautions

Before use, please read this SAFETY PRECAUTIONS carefully to ensure proper operation.

The following SAFETY PRECAUTIONS help you to use the product correctly. Indicate the risk of danger and damage, the classes of accidents caused by incorrect use are classified into WARNING and CAUTION, hence every user must abide by it.

Indication for WARNING AND CAUTION	
WARNING	Indicating user may be seriously injured once the product is improperly used.
CAUTION	Indicating the user may be injured or the product may be damaged once the product is improperly used.

1. Important safety warning

- Please use power supply suited to the product.
- Don't use the product if the power source does not meet the rated voltage to avoid damage or fire hazards.
- Please pull out the plug after use or before cleaning to avoid injury or damage to the product.
- Don't use accessories which are not recommended.
- Any other usage that is not listed in this manual instruction is forbidden.
- 15 minutes of usage each time is recommended.
- Please don't use it if the case or parts are broken or damaged.
- Please don't use it if the radiator hole is covered or blocked.
- Please don't use it within one hour after having meal.
- Please don't adjust to the greatest intensity at the beginning of the massage to avoid the risk of being injured.
- If the product is damaged, don't attempt to change on your own. Please contact us for repair for safety reason.

2. Environment for usage

CAUTION

- Please don't use it under high temperature and moist environment such as bathrooms.
- Please don't use it immediately while the environmental temperature changes sharply.
- Please don't use it under heavy-dusty or caustic environment.
- Please don't use it where there is not enough space or not good ventilation.
- Please don't use it when surrounding temperature is up to 40°C.
- Please don't use it when surrounding temperature is below 5°C.
- Don't place it near fire or directly under the sun.
- When moved from a low-temperature storage to a warm environment, the moisture inside the machine may affect, or damage, its functionality. Please rest the product within room temperature for 1 hour before use.

Safety Precautions

3. People who are not suitable to use the product

CAUTION

- People who are receiving medical treatment, especially physically unfit, are not suitable to use the product.
- People who have cancers, acute illnesses, heart diseases, or severe high blood pressure, etc. are not suitable to use the product.
- Pregnant women or women that are in menses period are not suitable to use the product.
- People who are suffering from fever are not suitable to use the product.
- People who are suffering from heart diseases or carrying electronic medical equipment such as a heart pacemaker are not suitable to use the product.
- People who got injured or skin diseases are not suitable to use the product.
- People who have rheumatoid arthritis, rheumatism, gout are not suitable to use the product.
- People who have osteoporosis are not suitable to use the product.
- People who have blood circulation problems such as thrombosis, severe aneurysm, or acute venous tumor, or infectious skin diseases are not suitable to use the product.
- People who have perception disorders due to diabetes are not suitable to use the product.
- People who have open wounds are not suitable to use the product.
- People whose health condition is poor are not suitable to use the product.
- Patients with tenosynovitis.
- Children under the age of 14 or individuals with mental instability must not use this machine without supervision.
- People who aim for recovery are not suitable to use the product.

4. Safety

CAUTION

- Don't use other therapy devices while using the product.
- The product is not to be used for medical treatment purpose.
- Recommended usage time of the product 15 minutes/ session.
- Don't pull or drag the power cord when unplugging.
- Don't place heavy objects on the power cord.
- Before use, check if the machine and its parts are damaged. If so, don't use it and contact us for repair.
- Don't tie, bend, pull or twist the power cord to avoid damages.
- Stop using the product if there appears any irregularity, and please contact us immediately.
- Stop using the product if it causes discomfort, and please consult a doctor.
- This appliance is not intended for people (including children) with reduced physical, sensory or mental capabilities, or a lack of experience and knowledge, unless they are under supervision or instruction concerning the use of the appliance by a person responsible for their safety.
- Children should be supervised to ensure that they do not play with the appliance.
- Keep the power cord away from hot or moving objects.
- The product is not designed for long-time use. Allow a 30-minute break (power off) for every one hour use. Long-time use could cause overheating and malfunction.

Safety Precautions

FORBIDDEN

- It is forbidden to use if there appears dilapidation, crack, leakage or explosion of internal parts on the product.
- It is forbidden to place heavy objects on top of the product.
- It is forbidden to let children use or play with the product.
- It is forbidden to use the product while preparing to sleep or sleeping.
- It is forbidden to use while drunk or sick.
- It is forbidden to put fingers or other foreign objects into the gaps, to prevent harm or malfunction.
- In case of irregular situation (e.g. smoking or smelling), please immediately stop using and cut off power.
- It is forbidden to have water or any liquid on the main parts of the product, to avoid electric shock or fire.
- It is forbidden to operate the console or forcefully drag the power cord with wet hands.
- People whose bodies are wet are forbidden to use the product.

DISASSEMBLY FORBIDDEN

Anyone else, except maintenance personnel specified by our company, must not disassemble, maintain, or remodel the product. If any problem, please contact us immediately.

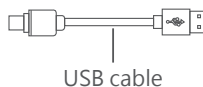
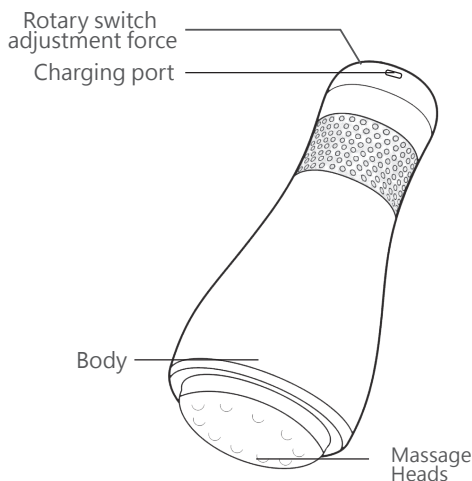
5. Maintenance

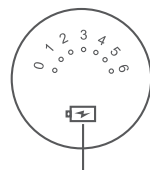
- The product should be maintained by us, users are prohibited to disassemble or maintain by themselves.
- If the product is to be left idle for long, put away the product and keep it in dry and dustless environment.
- Don't keep the product at a high temperature or near fire. Don't keep under direct sunlight.
- Please clean the product with dry cloths. Don't use thinner, benzene or alcohol.
- The mechanical components of the product are specially designed and made, it is not necessary to be specially maintained.
- Don't point sharp objects at the product.
- Please use it intermittently; don't keep it continuously running for long.

 **Name and Function of Components**

Product structure

Speed indicator

charging indicator

(Massage head configuration subject to actual packaging)

 **How to Use**

- 1 · When in use, press the top of the body to turn on the power.
- 2 · Rotate the top of the body to increase or decrease the massage intensity, with six levels of massage strength available for selection.
- 3 · For first use, start with lower massage intensity and gradually increase after adaptation.
- 4 · After use, press the power switch to turn off.
- 5 · It is recommended not to use the product for more than 15 minutes.
- 6 · Slight heating of the device during use is normal.

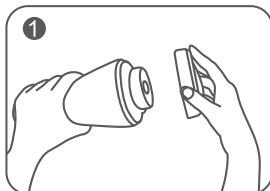
Proper use of the wear-resistant mesh cover protects hair.

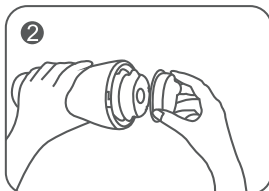
Improper use of this device may cause hair entanglement.

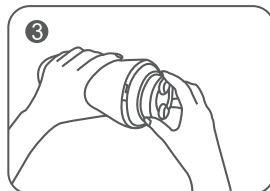
Massage head replacement steps

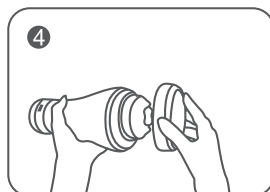
Hold the main unit with your left hand, grasp the ears on the mesh cover with your right hand, and pull the mesh cover off the device.

Find the two small gaps where the massage head is connected to the body and pull it out firmly.

Install your desired massage head and align it firmly into place.

Hold the body with your left hand, and with your right hand, press the cover on until it clicks securely into place on the device.

Charging Method

Remarks

1. Before first use or after prolonged non-use, please fully charge first.
2. With the product powered off, connect the charging cable to the massager's charging port.
3. Insert the other end of the charging cable into a 5V power source. Charge for approximately 3-5 hours. The massager displays a red battery icon while charging, and a green battery icon when fully charged.
4. After charging is complete, disconnect the charging cable from the power source.

Note:

The massager may become slightly warm during charging, which is normal. The device cannot operate while charging (automatically powers off).

Tips

It is recommended not to use the product for more than 15 minutes at a time;

 **Cleaning and Maintenance****Cleaning**

- Before cleaning, please make sure the switch has been turned off, and the power plug has been pulled out from socket.
- Use soft cloths to clean the surface of the product and the control interface.

**CAUTION**

- Do not pull the plug from socket with wet hands.
- Rinsing the product with water is not allowed.
- Do not use corrosive cleaning agents like benzene or thinner to clean the product, otherwise, the color of product will fade, and the internal parts may be corrupted or broken.
- Do not disassemble the product when it is wet.

Maintenance

- When not in use, turn off the main power, and remove the plug from socket.
- Keep the product clean and away from high-temperature/ humid areas, or any environment with volatile gases.
- If the product is to be left idle for long, properly cover it with a dust cover (self-purchase).

**CAUTION**

- Please keep the product away from direct sunlight, or high temperature, to avoid fading.

Kindly Reminder

If the product has been used for long, please check the conditions as below:

- Burnt smell.
- Unstable power supply after contact with the power cord.
- Power cord heating up.
- Power cords damage.
- Other abnormalities.

**CAUTION**

- If above mentioned happens, please stop using the product.
- Switch off the power supply, pull out the plug, and contact us.
- Don't disassemble the product by self.

Troubleshooting

If you encounter problems while using the product, please refer to the following tips and solutions. If the product remains faulty or unusable, please turn off the power. Users must not disassemble or repair the product themselves. Self-disassembly or repair causing damage will void the product warranty, and the company will not be responsible for any consequences arising from this situation.

Issue	Reason	Solution
Massager won't turn on	Power not started.	Short press the top body switch to power on, then rotate the switch to select massage mode.
	Low battery level.	Charge the product.
Cannot charge or no charging light displayed	Charging cable not properly inserted	Please check if the charging cable and device charging port are properly connected.
	Not using the designated charging cable.	Use the designated charging cable for this device.
Sudden stop during operation	15-minute massage program completed and auto power-off.	Restart
	Low battery	Charge the product

Technical Specifications

Product Name	: 360°ShapeLift
Product Model	: MAX-W02YC212W
Operating Voltage	: 5V
Rated Power	: 25W
Charging voltage	: 5V (A charging device that meets the product specifications must be used)
Charging Current	: 1A
Charge Time	: 4-5Hrs
Battery Capacity	: 2000mAh
Battery Life	: 1-2Hrs
Product Size	: 90x90x210mm
Product Net Weight	: 765g

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Trademark holder: Maxcare Industrial Limited