Phototherapy Hand Care 手康寶

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MAX-M01MC0653

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User Manual

MAX-M01MC0653

Phototherapy Hand Care, the product technology developed by MTCURE in collaboration with The Chinese University of Hong Kong, and maxcare provides engineering and manufacturing supports, is a product that effectively helps improve pain conditions while promoting metabolism and blood circulation.

MTCURE Biotech Ltd, an innovative biotechnology start-up incubated by The Chinese University of Hong Kong.

The clinical convenience sample testing of the related technology was jointly conducted by MTCURE(鎂天健康) and the Faculty of Medicine of The Chinese University of Hong Kong.

COMPANY PROFILE

After two decades of experience in the health industry, **maxcare** founder utilized his knowledge and expertise in massage therapy from an Asian-Pacific perspective to establish our brand - **maxcare**. Our massage therapy business was expanded to Hong Kong in 2012 and is conferred 2017 Hong Kong Top Brand Award. Based in Hong Kong, we are specializing in design & innovation, manufacturing, and retailing of massage and healthcare products.

Believing in the value of 'Technology improves our quality of life' and 'Caring for your family', **maxcare** has been innovating highly qualified products with international patents. Besides, **maxcare** products fulfill various international safety standards, namely EU CE Mark, US UL/ETL Mark, Germany GS Mark, TUV certification, and China CCC Mark.

Motivated by our value, **maxcare** founder has been leading us to be one of the most competitive brands in Hong Kong health industry. **maxcare** applies the contemporary fashion design – neat and modern – on the mechanical engineering of massage chairs. To deliver our value and offer incredible massage experience and healthcare solution to the public, **maxcare** devotes to the innovation, design, manufacturing the top-notch family massage products. By integrating the company philosophy – Spirit Of Healthy Life - into our products, **maxcare** brings health to everyone.

Phototherapy Hand Care

Thank you for purchasing Maxcare's products. Before using the product, please read this manual thoroughly, particularly the Safety Precautions section, for the correct use of the product.

This manual should be well kept for easy reference. It enables you to perform a safe and efficient operation.

CONTENT

Safety Precautions	P.01 - P.03
Product Introduction	P.04
Name and Function of Components	P.05
How to Use	P.06
Cleaning and Maintenance	P.07
Troubleshooting and FAQs	P.08
Technical Specifications	P.09

Safety Precautions

1. Important safety warning

- Please use a well-grounding power supply suited to this product.
- Please pull out the plug after use or before cleaning to avoid injury or damage to the product.
- Do not use unauthorized accessories.
- Please read this manual instruction carefully before operation.
- Any other usages that are not listed in this manual instruction are forbidden.
- Please check if any part of this product is damaged before use. If so, please do not use it and contact us for repair.
- Please do not use it if the product is broken or damaged.
- Please do not use it if the cover or surface is damaged.
- Please do not drop anything into this product.
- Please do not allow water or liquid into this product.
- Please do not fall asleep while using this product.
- Please do not use it if drunk or feeling unwell.
- Keep the plug clean from needles, rubbish, or water.
- If the power cord is damaged, do not attempt to change on your own. Please contact us for repair for safety reasons.

2. Environment for usage

- Please do not use it in a moist environment such as a bathroom.
- Please do not use it when there is a drastic change in temperature.
- Please do not use it under a heavy-dusty or caustic environment.
- Please do not use it where there is not enough space or not a good ventilation.
- Please switch off the power after use.
- Please do not use it when the surrounding temperature is up to 50°C.
- Please do not use it when the surrounding temperature is below -20°C. The lubricant can thicken quickly and increase the friction between parts.
- Turn off and pull out the plug during thunderstorms.
- Do not place it near fire or directly under the sun.
- When moved from a low-temperature storage to a warm environment, the moisture inside the machine may affect, or damage, its functionality. Please rest the product within room temperature for 1 hour before use.

Safety Precautions

3. People who have the below symptoms, please consult a doctor before use

- Children under the age of 13 and pregnant women should not use the light therapy product;
- Cancer patients;
- Tuberculosis patients;
- Acute suppurative inflammation;
- Hemorrhagic patients;
- Heart disease and pacemaker patients;
- Arteriosclerosis patients;
- Newborns and infants;
- Patients who are allergic to light or have skin that is insensitive to heat are considered in appropriate to use these products by doctors.

CAUTION

- Please note this is not a medical device.

4. Precautions before using the product and some safety instructions

- Do not use the NMES and Heat functions on area that have implanted metal or other materials, as well as areas with inflammation. It is recommended to use Red Light and Infra-Red Therapy instead and keep the product at a distance of 3-5mm from the skin surface.
- Avoid using the NMES and Heat functions on area with wounds, pimples, or suppurative diseases. Once the wounds have stopped bleeding or the tissue sprain has healed, it is advisable to use only Red Light and Infra-Red Therapy to promote faster healing. After 72 hours or when the wound has fully healed, it is safe to use the machine.
- Do not directly view the light. (Included wear sunglasses or protective spectacles)
- Please do not place heavy objects on this product.
- Do not use other therapy devices while using this product.
- Please note this is not a medical device, it is not to be used for medical treatment purposes. Only use for improvement and relief purposes.
- The product is set to automatically turn off after 10 minutes. Recommended usage time of this product is 10 minutes/ session.
- Check if the voltage is suitable for the specification of this product.
- Do not use wet hands to pull the plug and use the machine.
- Do not allow water to get into the product, for the avoidance of electric shock or damage to this product.
- Do not pull or drag the power cord when unplugging. Please be gentle.
- Do not damage the wires or change the circuit and alter the structure or parts of this product.

Safety Precautions

- Do not use wet cloths to clean the electric parts such as the switch and plug.
- Stop using this product if it appears any irregularity, and please contact us immediately.
- Stop using this product if it causes discomfort, and please consult a doctor.
- This appliance is not intended for people (including children) with reduced physical, sensory or mental capabilities, or a lack of experience and knowledge, unless they are under supervision or instruction concerning the use of the appliance by a person responsible for their safety.
- Children should be supervised to ensure that they do not play with the appliance.
- If the power cord is damaged, it must be replaced by us in order to avoid a hazard.
- Keep the power cord away from hot or moving objects.
- Please cut off power and pull out the plug before moving.
- If the product is found damaged, broken, suffering from power leakage, or with exposed parts, please stop using it.
- This product contains a lithium-ion rechargeable battery.
- Do not place this product near fire or where it causes fire easily to avoid accidents; the product may also deform, fade, or break.

5. Maintenance

- This product should be maintained by us, users are prohibited to disassemble or maintain by themselves.
- If this product is to be left idle for long, please keep it in a dry and dustless environment.
- Do not keep this product at a high temperature or near a fire. Do not keep under direct sunlight.
- Please clean this product with dry cloths. Do not use thinner, benzene, or alcohol.
- Avoid contact with sharp edges or pointed objects which might cut or puncture the fabric surface. To avoid breakage.
- Do not fold during storage and use to avoid damage to the lamp beads.
- Please use it intermittently; do not keep it continuously running for long.

6. Solution of ordinary malfunctions

- You will hear a "sizzling" sound when using it. This is the sound of current flowing. This is not a product fault.

Product Introduction

This is not a medical device. This unit is intended to be used only for general physical conditioning. No medical claims are implied or warranted by the use of this product.

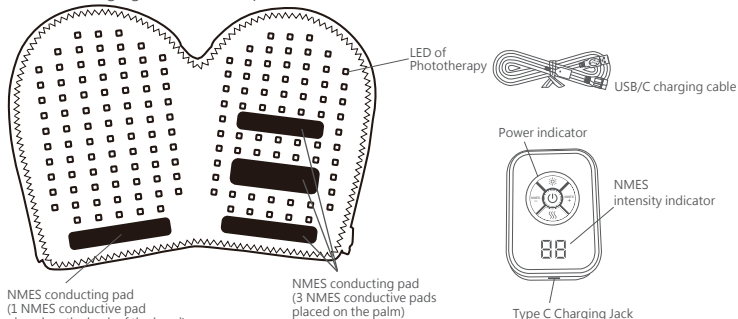
Product Function:

NMES	Electronic Muscle Stimulation helps to stimulates muscle blood flow, relief muscle fatigue, neural rehabilitation, enhance muscle strength.
Heat	Hot compress can help to relax tired and sore muscles, available with 2 levels from 41°C and 45°C(+/- 5°C)
MT-PBMT (MT-Photobiomodulation)	MT-PBMT combines light in the range of visible light and invisible light, starting at superficial tissue to deeper tissue levels, it promotes wound healing and tissue repair, relieves pain and inflammation.

*MTCURE (Mtcure Biotech Limited) biophysical combined therapy can effectively enhance the amplification effect of various biological therapy signals, thereby improving rehabilitation efficiency.

Name and Function of Components

1. Phototherapy Hand Care x 1 pc
2. USB/C charging cable x 1 pc

How to Use

Operation Instruction

	Switch button - press it for 3 seconds to switch on/off the product.
NMES +	NMES function increase intensity button, press it to turn on NMES function and increase the intensity from 0- 10 levels. If needs to stop NMES immediately, the Long press switch button for 3 seconds to turn off the device.
NMES -	NMES function decrease intensity button, press to decrease the intensity from 10-0 levels. Long press switch button for 3 seconds to turn off the NMES function.
	Heat button - when selecting heat function, this icon will light up in blue for 1st Level 41°C> red 2nd Level 45°C (+/- 5°C)
	Phototherapy when select photo treatment, this icon will light up in red. When selected this function, the keypad will light up in red, press once to turn on, press again to off photon.

How to Use

How to use

Removal: Remove the Phototherapy Hand Care from the packaging and unfold it to confirm the position of the black microcurrent conduction pad. It can be used with either hand.

Positioning: Align your palm and the back of your hand with the corresponding position and place it inside the Phototherapy Hand Care.

Secure: Then zip up the zipper to secure the product and ensure optimal performance.

1. Press and hold the power button on the controller for 3 seconds to turn on the device. The power indicator light will turn blue, indicating that the product is in standby mode.
2. Select the function that you want. If you have selected the NMES function, suggest you apply some water base gel or water on your hand to enhance the conductivity of your skin. And make sure the Phototherapy Hand Care is contact the skin firmly, Otherwise the electric impulse cannot be transmitted to your muscle.
3. Please don't use 1 position for over 20 mins.
4. Low voltage warning, when the battery voltage is lower than 3.35V, the power indicating light will flash, flashing for 15 seconds to shut down. Then you should charge the product.

Remark: It is recommended to turn off the device to allow it to cool down after using it continuously for 2 sessions or being powered on for more than 20 minutes, in order to prolong its service life.

Charging instruction

1. Before the first use, charge the product for about 4 hours and fully charge it.
2. Plug in the USB/C cable to the type C charging jack of the device and the other end of the USB plug into the adaptor with 5V/2A output. The power indicator blinking in blue.
3. When all the power indicator stop flashing and keep on, it means the battery is fully charged.
4. Charge or discharge batteries to approximately 50% capacity before storage.
5. Recharge the battery to approximately 50% capacity at least every 3-6 months.
6. Store product at a temperature between 5 °C and 25 °C.
7. Do not use while charging.

Cleaning and Maintenance

Cleaning

- Please ensure to turn off the power and unplug before cleaning.
- Use a soft dry cloth or slightly damp cloth to wipe the product surface.

CAUTION

- Do not touch the plug with wet hands.
- Do not wash or spill any liquid into the product.
- Do not use corrosive liquid such as benzene or thinner to clean the product, to prevent color fading and damage.
- When the product is damp, stop using it and do not dismantle. Leave it to dry naturally.

Maintenance

- Keep the whole set of the product clean and storage in a dust free and dry location under the condition. Storage temperature and humidity: -10oC to 60oC, 10% to 95% RH
- When storing the product, please avoid folding it or applying excessive pressure to prevent internal wire breakage and damage to the lamp beads.

CAUTION

- Do not place it near fire or direct sunlight, to prevent product color fading.

Kindly Reminder

- If the product has been used for a long time, please check the conditions as below:
 - * Burnt smell.
 - * Unstable power supply after connecting with the power cord
 - * Power cord overheat
 - * Power cords damage
 - * Other abnormalities

CAUTION

- If the above phenomenon occurs, please stop using the product to prevent malfunctions and accidents.
- Please turn off the power and unplug immediately. Users are prohibited to dismantle and repair by themselves. Please contact us for repair.

Troubleshooting and FAQs

The following symptoms do not indicate malfunction. Please check the below phenomenon to find out the solution before contact customer services department. If the problem still cannot be solved, please do not dismantle or repair the product. You should stop using it and contact us immediately.

Problem	Possible Causes	Solutions
No Response The device blinking	- Low in Battery	- Charge the product
Can't feel NMES upon select NMES function	- Poor skin conductivity - Too low intensity for you	- Use water-based gel on your skin and make sure the black color NMES conductive pad is contact on your skin firmly. - Increase the intensity level gradually.
NMES intensity cannot goes up, always drop to 00	- No or Poor skin contact	- Use water-based gel on conductive area and place pad directly onto skin. - Wrap firmly the product. Make sure the black color NMES conductive pads are contact on your skin firmly.
All buttons do not work well.	- Microprocessor hang up - Battery is low	- Restart the unit 1 minute after turning it off. - Recharge the product

FAQs

Q : I feel a bit warm during photon therapy, why ?

A : the warm is from the light during activating, which is normal.

Q : When I wipe the electrode area, it becomes black and dirty.

A : If I rub the electrode area forcefully, the color may transfer. However, it will not dirty the skin during the care process, so it is safe.

Q : After use, symptoms such as redness, inflammation, and itching may occur.

A : Individual responses to NMES may vary. If the symptoms do not improve the next day, please consult a healthcare professional.

Q : Feeling discomfort and fatigue after use.

A : Individual responses to NMES may vary. Before the discomfort subsides, please temporarily reduce the frequency of use.

**Technical Specifications**

Product name	: Phototherapy Hand Care
Model no.	: MAX-M01MC0653
Charging voltage	: 5V USB 2A (A charging device that meets the product specifications must be used)
Battery capacity	: 3.7V 5000mAh
Max. Rated power	: 13W
Max. output Current	: 3500 mA
NMES – Max. Intensity	: 1 mode with 10 intensity
NMES – Frequency Generation	: 1.66 – 16.67 Hz
NMES - Max. output voltage	: 16.0V – 42.1V (under 1KΩ)
Heat	: 2 levels : Blue for 1st Level 41°C> Red 2nd Level 45°C (+/- 5°C)
Photon	: 151 pcs of 660nm red light (visible light), 151 pcs of 810nm near-infrared (invisible light)
Charge time	: About 180 mins
Work for when full charge	: About 100 mins
Auto Timer	: 10 minutes

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客戶服務熱線 Customer Service Hotline : (852) 2950-9383

www.maxcareonlineshop.com**鎂天健康生物科技有限公司****MTCURE BIOTECH LIMITED** an innovative biotechnology start-up
incubated by The Chinese University of Hong Kong.